

Lab_Notebook_2024

Callum Malcolm

Contents

January 2024	7
Tuesday 16-01-2024	7
WILDseq - Mouse Experiment 1: PCR 1 Samples 1-8 / PCR 2 Sample 5	7
Overview:	7
PCR 1: Samples 1-8	8
Bead Clean Up	8
PCR 2: Samples 5	9
Wednesday 17-01-2024	12
WILDseq - Mouse Experiment 1: PCR 1 Samples 9-30	12
Overview:	12
PCR 1: Samples 9-30	12
Bead Clean Up	12
Thursday 18-01-2024	15
WILDseq - Mouse Experiment 1: PCR 1 Samples 17/18/29	15
Overview:	15
PCR 1: Samples 17/18/29	15
Bead Clean Up	15
Friday 19-01-2024	16
WILDseq - Mouse Experiment 1: RNA Extraction/RT Samples 31-33	16
RNA extraction Samples 31-33	16
RT Protocol Samples 31-33	16
Monday 22-01-2024	17
PCR 1: Samples 31-33	17
Bead Clean Up	17
PCR 1: Negative Control	19
Tuesday 23-01-2024	22
Made New Media: AR-5	22
Cell Culture	23
A20	23
HEK	23
WS-ME1 Library Prep: Control Test 3	23
PCR 1: Negative Control	23
Wednesday 24-01-2024	25
Cell Culture	25
A20	25
HEK	26

WILDseq - Mouse Experiment 1: Library Prep Attempt 5 - RT	26
RT Protocol	26
Thursday 25-01-2024	27
Cell Culture	27
A20	27
HEK	27
WS-ME1 Library Prep - Attempt 5 PCR 1	27
PCR 1: Negative Control	27
Bead Clean Up	28
Tapestation PCR 1 Samples 1-33	28
Friday 26-01-2024	32
Cell Culture	32
A20	32
HEK	32
WS-ME1 Library Prep - Attempt 5 PCR 1 repeat and PCR 2 prep	32
PCR2 Prep	32
PCR1 Samples 18-20	32
Bead Clean Up	33
Tapestation PCR 1 Samples 18-20	33
Monday 29-01-2024	33
Cell Culture	33
A20	33
HEK	35
PCR Step 2 - Samples 1-24	35
Overview	35
PCR2 Protocol	35
Tapestation	36
Tuesday 30-01-2024	37
Cell Culture	37
HEK	37
A20	37
February 2024	38
Thursday 01-02-2024	38
PCR Step 2 - Samples 1-8, 25-39	38
Overview	38
PCR2 Protocol	38
Tapestation	39
Monday 05-02-2024	40
Cell Culture	40
Made DMEM	40
HEK - Seeded	40
A20 - Split	40
Tuesday 06-02-2024	40
WS-ME1 Pooling for Submission	40
Wednesday 07-02-2024	41
Cell Culture	41

HEK	41
A20	41
Friday 09-02-2024	42
Cell Culture	42
A20	42
Monday 12-02-2024	42
Cell Culture	42
Made DMEM-V1	42
Thawed HEK	42
A20	42
Processing scRNAseq	42
Tuesday 13-02-2024	42
Cell Culture	42
HEK	42
Processing scRNAseq	42
Wednesday 14-02-2024	42
Cell Culture	42
A20	42
Processing scRNAseq	43
Friday 16-02-2024	43
A20 Xenograft EXP 1	43
Cell Prep	43
Injection	43
Monday 19-02-2024	43
Cell Culture	43
A20	43
Wednesday 21-02-2024	44
Cell Culture	44
A20	44
HEK	44
Friday 23-02-2024	44
Cell Culture	44
A20	44
HEK	44
Monday 26-02-2024	44
Cell Culture	44
A20	44
HEK	44
WILDseq Virus Production Day 1	45
March 2024	45
Wednesday 06-03-2024	45
Wednesday 13-03-2024	45
Thursday 14-03-2024	45

Friday 15-03-2024	45
Cell Culture	45
N2-BC4 RTX CDC Assay - version 1	45
Monday 18-03-2024	45
Cell Culture - N2 -C4	45
N2-BC4 RTX CDC Assay - version 2	45
Friday 22-03-2024	45
Cell Culture	45
N2-BC4	45
EC50 RTX N2-BC4 22324 - Seeding	46
Monday 25-03-2024	47
EC50 RTX N2-BC4 22324 - Collection	47
EC50 RTX N2-BC4 25324 - Seeding (DNW)	47
Cell Culture	48
N2-BC4	48
Wednesday 27-03-2024	49
Thursday 28-03-2024	49
Friday 29-03-2024	49
EC50 RTX N2-BC4/RAMOS-BC5 25324 - Seeding (DNW)	49
Sunday 31-03-2024	50
Cell Culture	50
N2-BC4	50
Ramos-BC5	50
Ramos CDC Testing	50
EC50 RTX N2-BC4/RAMOS-BC5 24331 - Seeding (DNW)	51
April	52
Monday 01-04-2024	52
Cell Culture	52
N2	52
Ramos	53
EC50 Collection: 24329/24331	53
EC50 RTX N2-BC4/RAMOS-BC5 24401 - Seeding (DNW)	53
Tuesday 02-04-2024	54
N2/Ramos CD20 Flow Cytometry	54
Thursday 04-04-2024	55
Cell Culture	55
N2	55
Ramos	55
EC50 Collection: 240401	55
Monday 08-04-2024	55
Cell Culture	55
N2	55
Ramos	55

Tuesday 09-04-2024	56
Perla Drug Randomizer	56
Saturday 14-04-2024	56
Cell Culture	56
N2-BC4	56
Ramos-BC5	56
EC50 RTX N2-BC4/RAMOS-BC5 24414 - Seeding (DNW)	56
Wednesday 17-04-2024	57
EC50 Collection: 240414	57
Monday 22-04-2024	58
Cell Culture	58
N2-BC4	58
Ramos-BC5	58
Thawed NK-92	58
Wednesday 24-04-2024	58
Cell Culture	58
N2-BC4	58
Ramos-BC5	58
NK-92	58
Friday 26-04-2024	58
Cell Culture	58
N2-BC4	58
Ramos-BC5	58
NK-92	58
Monday 29-04-2024	59
Cell Culture	59
N2-BC4	59
Ramos-BC5	59
Tuesday 30-04-2024	59
Cell Culture	59
N2-BC4	59
Ramos-BC5	59
Human Serum Reciept	59
Ramos/N2 CDC Testing	59
EC50 RTX RAMOS-BC5 24430 - Seeding	60
May	61
Wednesday 01-05-2024	61
Cell Culture	61
N2-BC4	61
Ramos-BC5	61
June	61
Monday 10-06-2024	61
Cell Culture	61
EC50 RTX RAMOS-BC5 240610 - Seeding	62

Wednesday 13-06-2024	63
Cell Culture	63
EC50 Collection: EC50_RTX_N2-BC4_10064	63
Thursday 13-06-2024	64
Cryopreservation - RAMOS-BC 1/3/5	64
Protocol	64
Monday 24-06-2024	64
Cell Culture	64
Making RPMI	64
Splitting Ramos BC 3	64
Monday 24-06-2024	65
Cell Culture	65
CDC Tests	65
July	65
Wednesday 10-07-2024	65
Cell Culture	65
RAMOS - RTX In Vitro CDC Pressure	66
Freezing Down Cells Protocol	66
RAMOS RTX DP - Dose 1	66

January 2024

Tuesday 16-01-2024

WILDseq - Mouse Experiment 1: PCR 1 Samples 1-8 / PCR 2 Sample 5

Overview:

- PCR tests to see if PCR would work for WS-ME1 cDNA previously made
- Started with samples 1-8 so as not to waste sample/reagents if PCR 1 did not work

Sample List

Sample ID	Treatment/Sample	ng/uL	i7 index	i5 index
1	Cyclophosphamide	773.7	N701	S502
2	Cyclophosphamide	545.3	N702	S502
3	Cyclophosphamide	903.6	N703	S502
4	Cyclophosphamide	1056.1	N704	S502
5	Cyclophosphamide	959.1	N705	S502
6	Combination	730.1	N706	S502
7	Combination	602.3	N707	S502
8	Combination	449.6	N710	S502
9	Combination	1002.1	N701	S503
10	Combination	1929.2	N702	S503
11	Methotrexate	861.5	N703	S503
12	Methotrexate	1110.4	N704	S503
13	Methotrexate	1171.0	N705	S503
14	Methotrexate	1347.8	N706	S503
15	Methotrexate	891.4	N707	S503
16	Vehicle	374.5	N710	S503
17	Vehicle	911.2	N701	S505
18	Vehicle	829.8	N702	S505
19	Vehicle	600.1	N703	S505
20	Vehicle	750.3	N704	S505
21	Baseline	401.5	N705	S505
22	Baseline	443.3	N706	S505
23	Baseline	373.4	N707	S505
24	Baseline	444.5	N710	S505
25	Baseline	267.6	N701	S506
26	BLLW 14K	378.7	N702	S506
27	BLLW 2K	362.7	N703	S506
28	BLLW 1K	563.8	N704	S506
29	Methotrexate (outlier)	348.1	N705	S506
30	Combo	380.8	N706	S506
31	RBL2P 2K	173.3	N707	S506
32	RBL2P 7K	2708.0	N7010	S506
33	RBL2P 250K	1418.5	N701	S507
34	Mock_direct_1	-	N702	S507
35	Mock_direct_2	-	N703	S507
36	Mock_direct_3	-	N704	S507
37	Mock_culture_1	-	N705	S507
38	Mock_culture_2	-	N706	S507
39	Mock_culture_3	-	N707	S507

PCR 1: Samples 1-8

- Checked [cDNA] on NanoDrop of sample 8 to make sure samples were still intact
 - Sample 8 $\text{ng}/\mu\text{L} = 1042.3$
 - It is assumed all other cDNA is of similar quality

1. Made a master mix of PCR1 reagents

- Made enough for 31 samples
- Primer mix was made earlier

Component	Volume	Master Mix
10uM WS PCR1 Primer Mix	3 μL	93 μL
DNase/RNase H20	12 μL	372 μL
Kapa Hifi HotStart Ready Mix (2X)	25 μL	775 μL

2. Add following components to tubes

Component	Volume
PCR1 MasterMix	40 μL
cDNA	10 μL

3. Performed PCR using the following parameters:

Step Name	Steps	Time
Initial Denaturation	Step 1: 95C	3mins
Denaturation	Step 2: 98C	20s
Annealing	Step 3: 60C	15s
Extension	Step 4: 72C	15s
Final Extension	Step 5: 72C	1 min
Hold Step	Step 6: 12C	Hold

10-25 cycles of steps 2-4

Bead Clean Up

Overview

- Clean PCR reaction by removing leftover primers
- This prevents false reads, cleaner indexing in later steps, and higher quality samples

Materials

- For 40 samples


Component	Expected Volume/experiment	Material ID
Ampure Bead Kit	200 μL	

Protocol

1. Equilibrate Ampure beads to room temp (available from supplies cold room)
2. Mix well and vortex for 30sec to ensure uniform distribution.
3. Add 90ul (1.8X volume) to PCR reaction and mix thoroughly by pipetting.

4. Incubate 5min, room temp
 5. Place on magnet for 2 min and keep on magnet until final elution
 6. Remove supernatant and add 200ul 70% ethanol (make fresh)
 7. Incubate 30sec and remove
 8. Repeat 70% ethanol wash
 9. Air dry beads 5-10min room temp. Once dry the beads will go from shiny to matt in appearance. Avoid over drying (cracked appearance).
 10. Add 15 ul EB buffer and resuspend beads by pipetting
 11. Incubate 2min room temp
 12. Place on magnet and remove eluate
 13. Assess PCR product size, contamination and concentration on tapestation.
- Expected size = 165-172bp**

Tapestation PCR 1 - Samples 1-8

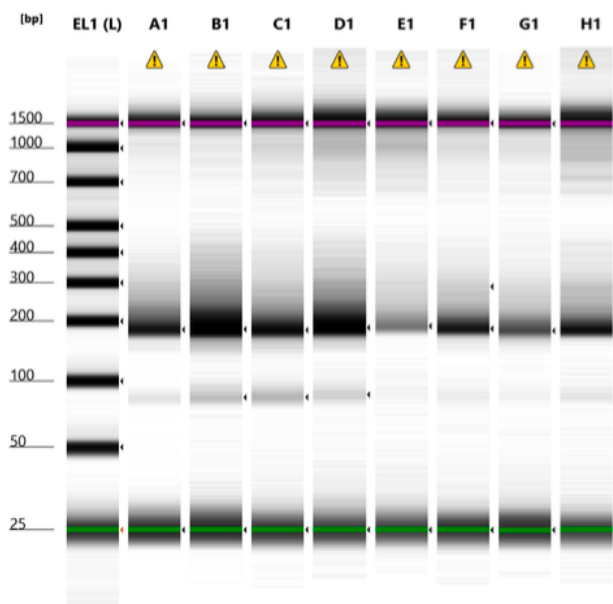

Default image (Contrast 100%)

Sample Info

Well	Conc. [pg/ul]	Sample Description	Alert	Observations
EL1	2350	Electronic Ladder		Ladder
A1	310	1 PCR1	⚠	Caution! Expired ScreenTape device
B1	631	2 PCR1	⚠	Caution! Expired ScreenTape device
C1	333	3 PCR1	⚠	Caution! Expired ScreenTape device
D1	369	4 PCR1	⚠	Caution! Expired ScreenTape device
E1	90.0	5 PCR1	⚠	Caution! Expired ScreenTape device
F1	266	6 PCR1	⚠	Caution! Expired ScreenTape device
G1	181	7 PCR1	⚠	Caution! Expired ScreenTape device
H1	206	8 PCR1	⚠	Caution! Expired ScreenTape device

Initial Test PCR 1

PCR 2: Samples 5

Overview

- Yield from PCR1 is low (.09 - .631 ng/ μ L)
- Sample with lowest concentration was taken ahead to PCR 2 (E1 - Sample 5)
- PCR to index samples for sequencing by attaching UMI

Materials

Component	Expected Volume/experiment	Material ID
10uM Nxxx Nextera i7 adapter	-	-
10uM Sxxx Nextera i5 adapter	-	-
DNase/RNase H20	1000 μ L	
Kapa Hifi HotStart Ready Mix (2X)	1200 μ L	KK2601

Protocol

1. Prepared samples according to the following table:

- Sample - Adaptor pairs in table above

Component	Volume
10uM N705 Nextera i7 adapter	1.5 μ L
10uM S502 Nextera i5 adapter	1.5 μ L
DNase/RNase H20	21 μ L
Kapa Hifi HotStart Ready Mix (2X)	25 μ L
10ng/ μ L PCR1	1 μ L

2. Performed PCR using the following parameters:

Step Name	Steps	Time
Initial Denaturation	Step 1: 95C	3mins
Denaturation	Step 2: 98C	20s
Annealing	Step 3: 55C	15s
Extension	Step 4: 72C	15s
Final Extension	Step 5: 72C	2 min
Hold Step	Step 6: 12C	Hold

8 cycles of steps 2-4

3. Performed Ampure bead clean up as above with 90ul beads per 50ul PCR reaction

4. Checked size and concentration on tapestation


Tapestation PCR 2

Review

- Contamination/PCR bubble?
- Yield is adequate (assuming it is representative of the correct peak)
- Chris wants to repeat PCR1 for remaining samples, take PCR2 ahead, and re-run same with D1000 tape (non High-Sensitivity)

Plan for tomorrow

- PCR 1 samples 9-30
- Bead Clean Up
- Tape Station

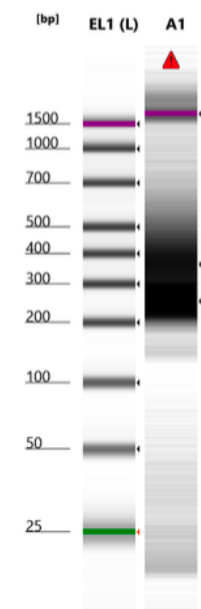

Default image (Contrast 100%)

Sample Info

Well	Conc. [pg/ul]	Sample Description	Alert	Observations
EL1	2350	Electronic Ladder		Ladder
A1	3180	5 PCR2	▲	Marker(s) not detected! Caution! Expired ScreenTape device

Figure 1: 2024-01-16 Tapestation PCR 2 Test Samples 5

Wednesday 17-01-2024

WILDseq - Mouse Experiment 1: PCR 1 Samples 9-30

Overview:

- Completed PCR 1 for Samples 9-30

PCR 1: Samples 9-30

- Used master mix of PCR1 reagents made 2024-01-16
- Add following components to tubes |Component |Volume| |:————|—| | PCR1 MasterMix | 40 μ L |
| cDNA | 10 μ L |
- Performed PCR using the following parameters:

Step Name	Steps	Time
Initial Denaturation	Step 1: 95C	3mins
Denaturation	Step 2: 98C	20s
Annealing	Step 3: 60C	15s
Extension	Step 4: 72C	15s
Final Extension	Step 5: 72C	1 min
Hold Step	Step 6: 12C	Hold

10-25 cycles of steps 2-4

Bead Clean Up

Overview

- Clean PCR reaction by removing leftover primers
- This prevents false reads, cleaner indexing in later steps, and higher quality samples

Materials


- For 40 samples

Component	Expected Volume/experiment	Material ID
Ampure Bead Kit	200 μ L	

Protocol

- Equilibrate Ampure beads to room temp (available from supplies cold room)
- Mix well and vortex for 30sec to ensure uniform distribution.
- Add 90ul (1.8X volume) to PCR reaction and mix thoroughly by pipetting.
- Incubate 5min, room temp
- Place on magnet for 2 min and keep on magnet until final elution
- Remove supernatant and add 200ul 70% ethanol (make fresh)
- Incubate 30sec and remove
- Repeat 70% ethanol wash
- Air dry beads 5-10min room temp. Once dry the beads will go from shiny to matt in appearance.
Avoid over drying (cracked appearance).
- Add 15 ul EB buffer and resuspend beads by pipetting
- Incubate 2min room temp

12. Place on magnet and remove eluate
13. Assess PCR product size, contamination and concentration on tapestation.
Expected size = 165-172bp

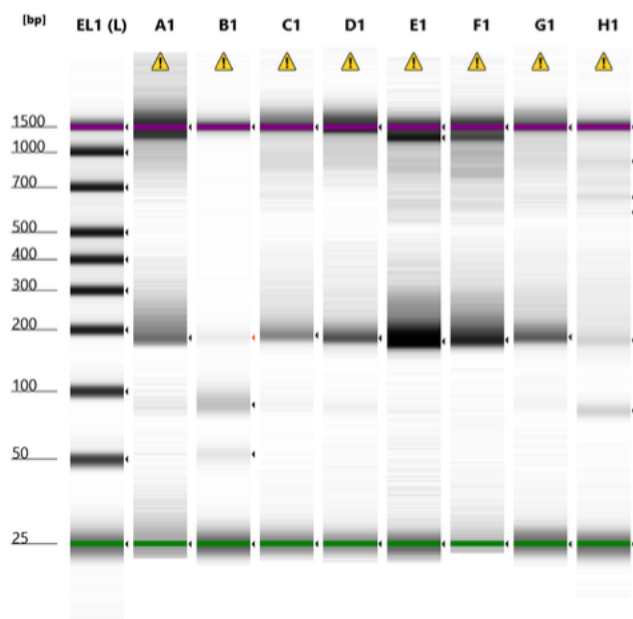

Default image (Contrast 100%)


Sample Info

Well	Conc. [pg/μl]	Sample Description	Alert	Observations
EL1	2350	Electronic Ladder		Ladder
A1	88.4	9 PCR 1	⚠	Caution! Expired ScreenTape device
B1	164	10 PCR 1	⚠	Caution! Expired ScreenTape device
C1	142	11 PCR 1	⚠	Caution! Expired ScreenTape device
D1	140	12 PCR 1	⚠	Caution! Expired ScreenTape device
E1	1070	13 PCR 1	⚠	Caution! Expired ScreenTape device
F1	354	14 PCR 1	⚠	Caution! Expired ScreenTape device
G1	203	15 PCR 1	⚠	Caution! Expired ScreenTape device
H1	203	16 PCR 1	⚠	Caution! Expired ScreenTape device

Figure 2: 2024-01-17 Tapestation PCR 1 Samples 9-16

Tapestation PCR 1

Tapestation PCR 1

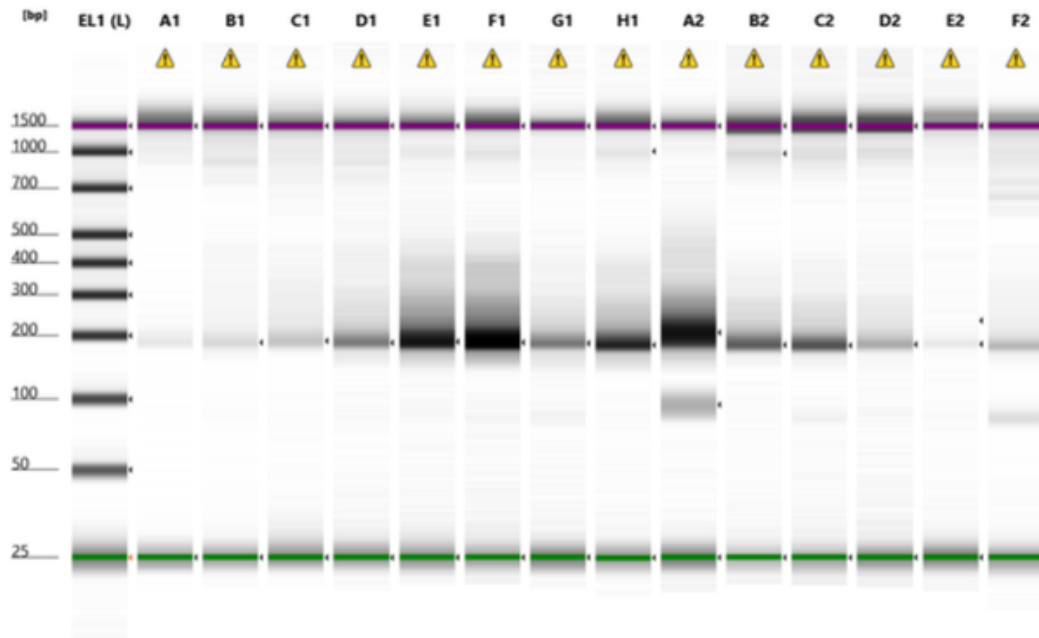

Default image (Contrast 100%)

Sample Info

Well	Conc. [pg/ul]	Sample Description	Alert	Observations
EL1	2350	Electronic Ladder		Ladder
A1		17 PCR 1	⚠	Caution! Expired ScreenTape device
B1	40.7	18 PCR 1	⚠	Caution! Expired ScreenTape device
C1	84.0	19 PCR 1	⚠	Caution! Expired ScreenTape device
D1	226	20 PCR 1	⚠	Caution! Expired ScreenTape device
E1	665	21 PCR 1	⚠	Caution! Expired ScreenTape device
F1	820	22 PCR 1	⚠	Caution! Expired ScreenTape device
G1	270	23 PCR 1	⚠	Caution! Expired ScreenTape device
H1	423	24 PCR 1	⚠	Caution! Expired ScreenTape device
A2	1290	25 PCR 1	⚠	Caution! Expired ScreenTape device
B2	223	26 PCR 1	⚠	Caution! Expired ScreenTape device
C2	170	27 PCR 1	⚠	Caution! Expired ScreenTape device
D2	71.5	28 PCR 1	⚠	Caution! Expired ScreenTape device
E2	30.5	29 PCR 1	⚠	Caution! Expired ScreenTape device
F2	157	30 PCR 1	⚠	Caution! Expired ScreenTape device

Review

- Samples 10/17/18/29 don't have adequate yield
 - Will need to redo PCR 1
 - Sample 10 has acceptable PCR1 sample from previous run
- Yield is adequate for the rest to move ahead

Plan for tomorrow

- PCR 1 samples 17/18/29
 - Bead Clean Up
 - Tape Station
- PCR 2 for all samples
 - Double check reagents amounts
 - Running low on KAPA (re-ordered) and Nextera Primers

- Chris wants to add 2 additional samples
 - To Do - RNA Extraction / RT / PCR 1 / PCR 2
 - Wait until majority of samples have acceptable PCR 2 product

Thursday 18-01-2024

WILDseq - Mouse Experiment 1: PCR 1 Samples 17/18/29

Overview:

- Repeated PCR 1 for Samples 17/18/29

PCR 1: Samples 17/18/29

1. Used master mix of PCR1 reagents made 2024-01-16
2. Add following components to tubes |Component |Volume| |—————:|:—:| | PCR1 MasterMix | 40 μ L | | cDNA | 10 μ L |
3. Performed PCR using the following parameters:

Step Name	Steps	Time
Initial Denaturation	Step 1: 95C	3mins
Denaturation	Step 2: 98C	20s
Annealing	Step 3: 60C	15s
Extension	Step 4: 72C	15s
Final Extension	Step 5: 72C	1 min
Hold Step	Step 6: 12C	Hold

10-25 cycles of steps 2-4

Bead Clean Up

Overview

- Clean PCR reaction by removing leftover primers
- This prevents false reads, cleaner indexing in later steps, and higher quality samples

Materials

- For 40 samples

Component	Expected Volume/experiment	Material ID
Ampure Bead Kit	200 μ L	

Protocol

1. Equilibrate Ampure beads to room temp (available from supplies cold room)
2. Mix well and vortex for 30sec to ensure uniform distribution.
3. Add 90ul (1.8X volume) to PCR reaction and mix thoroughly by pipetting.
4. Incubate 5min, room temp
5. Place on magnet for 2 min and keep on magnet until final elution
6. Remove supernatant and add 200ul 70% ethanol (make fresh)
7. Incubate 30sec and remove
8. Repeat 70% ethanol wash

9. Air dry beads 5-10min room temp. Once dry the beads will go from shiny to matt in appearance.
Avoid over drying (cracked appearance).
10. Add 15 μ l EB buffer and resuspend beads by pipetting
11. Incubate 2min room temp
12. Place on magnet and remove eluate
13. Assess PCR product size, contamination and concentration on tapestation.
Expected size = 165-172bp

Tapestation PCR 1

Review

- Samples 10/17/18/29 don't have adequate yield
 - Will need to redo PCR 1
 - Sample 10 has acceptable PCR1 sample from previous run
- Yield is adequate for the rest to move ahead

Plan for tomorrow

- PCR 1 samples 17/18/29
 - Bead Clean Up
 - Tape Station
- PCR 2 for all samples
 - Double check reagents amounts
 - Running low on KAPA (re-ordered) and Nextera Primers
- Chris wants to add 2 additional samples
 - To Do - RNA Extraction / RT / PCR 1 / PCR 2
 - Wait until majority of samples have acceptable PCR 2 product

Friday 19-01-2024

WILDseq - Mouse Experiment 1: RNA Extraction/RT Samples 31-33

Overview: Extra samples given by Chris for whitelist

RNA extraction Samples 31-33

- Performed RNA extraction of tissue/cell samples according to the following kit:

RT Protocol Samples 31-33

1. In PCR strip tubes, prepared 5 μ g of RNA in a total volume of 10 μ l of RNase/DNAase-free water.
2. Add 1 μ l of WS_RT_UMI_NexteraR2 primer (2 μ M stock)
3. Add 1 μ l dNTPs (10 mM each).
4. Denature at 65 C for 5 mins in the PCR machine, then straight onto ice for at least 2 mins, spin briefly to get liquid to the bottom of the tube.
5. Used Master Mix previously made
6. Add 7 μ l of RT MM prepared above to each sample and mix, spin briefly to get liquid to the bottom of the tube.
 - Spin briefly to get liquid to bottom of the tube
7. In the PCR machine incubate at 53 C for 10 mins followed by 80 C for 10 mins.
8. Add 1 μ l Themolabile Exonuclease I (NEB M0568) to remove excess RT primer.

Component	Volume	MM volume
Thermolabile Exonuclease I	1 μ L	
NEBuffer r3.1*	2 μ L	

Most PCR buffers are compatible

9. In PCR Machine: heat at 37C for 4 mins followed by 80C for 1 min
 - Spin briefly to get liquid to bottom of the tube
10. Add 1 μ l of RNase H and incubate at 37 C for 20 mins.
11. Checked [cDNA] on NanoDrop
 - Added to Sample List

Monday 22-01-2024

PCR 1: Samples 31-33

- Included negative control (MM + H2O)
1. Used master mix of PCR1 reagents made 2024-01-16
 2. Add following components to tubes

Component	Volume
PCR1 MasterMix	40 μ L
cDNA	10 μ L

3. Performed PCR using the following parameters:

Step Name	Steps	Time
Initial Denaturation	Step 1: 95C	3mins
Denaturation	Step 2: 98C	20s
Annealing	Step 3: 60C	15s
Extension	Step 4: 72C	15s
Final Extension	Step 5: 72C	1 min
Hold Step	Step 6: 12C	Hold

10-25 cycles of steps 2-4

Bead Clean Up

Overview

- Clean PCR reaction by removing leftover primers
- This prevents false reads, cleaner indexing in later steps, and higher quality samples

Materials

- For 40 samples


Component	Expected Volume/experiment	Material ID
Ampure Bead Kit	200 μ L	

Protocol

1. Equilibrate Ampure beads to room temp (available from supplies cold room)
2. Mix well and vortex for 30sec to ensure uniform distribution.
3. Add 90ul (1.8X volume) to PCR reaction and mix thoroughly by pipetting.
4. Incubate 5min, room temp
5. Place on magnet for 2 min and keep on magnet until final elution
6. Remove supernatant and add 200ul 70% ethanol (make fresh)
7. Incubate 30sec and remove
8. Repeat 70% ethanol wash
9. Air dry beads 5-10min room temp. Once dry the beads will go from shiny to matt in appearance.
Avoid over drying (cracked appearance).
10. Add 15 ul EB buffer and resuspend beads by pipetting
11. Incubate 2min room temp
12. Place on magnet and remove eluate
13. Assess PCR product size, contamination and concentration on tapestation.
Expected size = 165-172bp

Tapestation PCR 1 31-33

Tapestation PCR 1

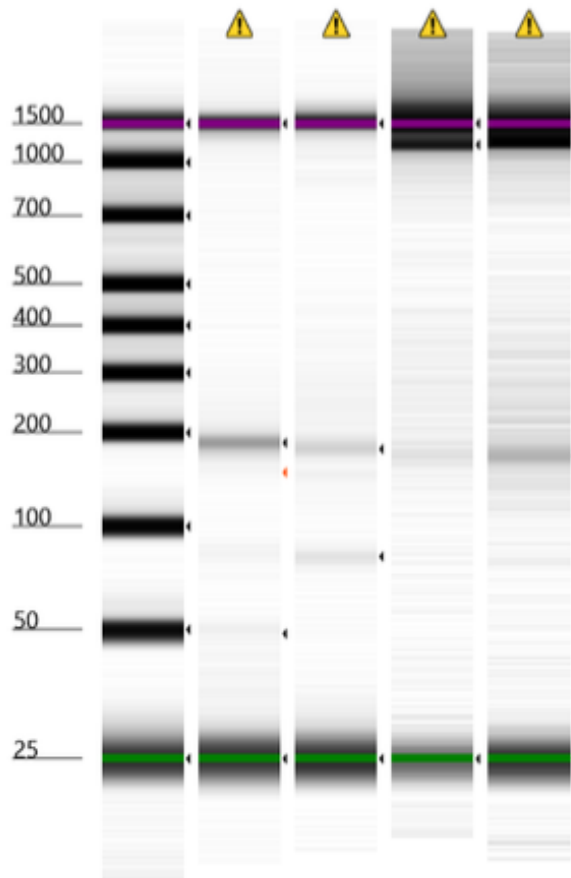

Default image (Contrast 100%)

Sample Info

Well	Conc. [pg/μl]	Sample Description
EL1	2350	Electronic Ladder
A1	89.4	-RT Control
B1	46.1	31 PCR 1
C1	87.3	32 PCR 1
D1	195	33 PCR 1

!!!! Band present in - Control lane !!!!

PCR 1: Negative Control

- Remade negative control (MM + H2O)
1. Used master mix of PCR1 reagents made 2024-01-16
 2. Add following components to tubes

Component	Volume
PCR1 MasterMix	40 μ L
cDNA	10 μ L


3. Performed PCR using the following parameters:

Step Name	Steps	Time
Initial Denaturation	Step 1: 95C	3mins
Denaturation	Step 2: 98C	20s
Annealing	Step 3: 60C	15s
Extension	Step 4: 72C	15s
Final Extension	Step 5: 72C	1 min
Hold Step	Step 6: 12C	Hold

10-25 cycles of steps 2-4

Tapestation Control Test 1

- Compared New -Control to old
 - Wanted to determine if this was a pipetting error

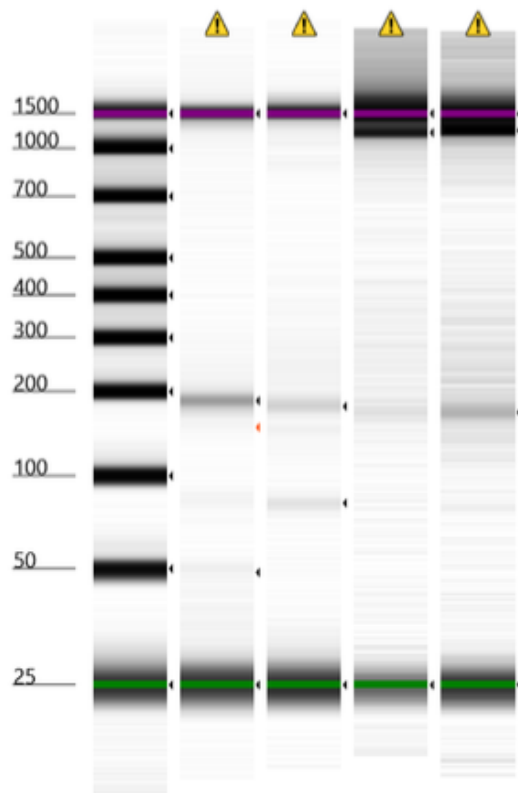


Default image (Contrast 100%)

Sample Info

Well	Conc. [pg/μl]	Sample Description
EL1	2350	Electronic Ladder
A1	89.4	-RT Control
B1	46.1	31 PCR 1
C1	87.3	32 PCR 1
D1	195	33 PCR 1

Tapestation Control Test 2

- Compared New -Control to H2O
 - Wanted to determine if water was contaminated

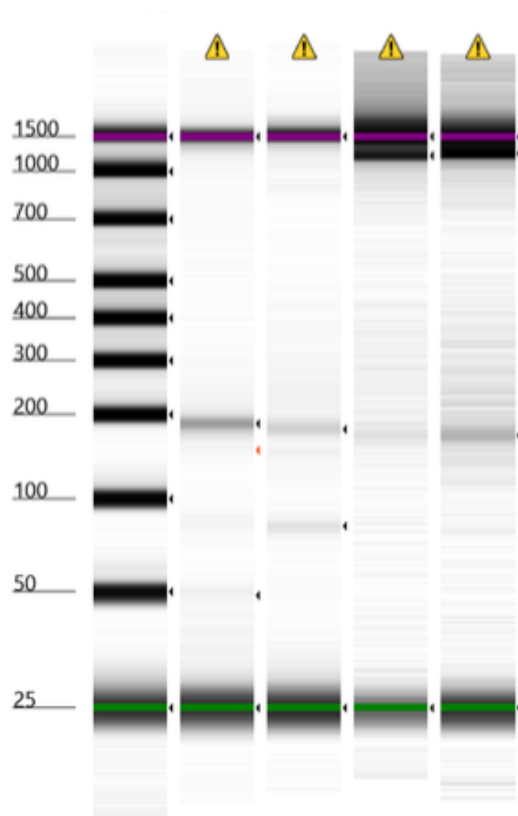

Default image (Contrast 100%)

Sample Info

Well	Conc. [pg/μl]	Sample Description
EL1	2350	Electronic Ladder
A1	89.4	-RT Control
B1	46.1	31 PCR 1
C1	87.3	32 PCR 1
D1	195	33 PCR 1

Tuesday 23-01-2024

Made New Media: AR-5

1. Made new media: AR-5

Solution	ID code	Volume	% Total volume
Advanced RMPI	—	500 mL	78%
Glutamax	—	6.5 mL	20%
Pen-Strep	—	6.5 mL	1%
FBS	—	128 mL	1%

Cell Culture

A20

- Thawed A20 and seeded in T25
- Spun down CS to remove DMSO, cell pellet small
- Cells didn't look super healthy
- Monitor tomorrow

HEK

- Got T75 flask from Chris
- 70% confluent
- Added 10mL of trypsin and tapped flask until cells were dislodged
- Added 10mL of DMEM (20% FBS) and added to 50mL tube
- Added 25mL of DMEM to 5 T75 flasks
- Mixed CS to get single cell suspension and remove clumps
- Transferred 4mL of CS to each flask
- Will grow up for 1-2 days and freeze down 4 T75 flasks
 - Other flask will be used to grow WILDseq library virus for A20's

WS-ME1 Library Prep: Control Test 3

Overview: question around contaminated primers - Ran PCR1 with new primer/PCR 1 mix

PCR 1: Negative Control

- Remade negative control (MM + H2O)
- Remade Primer mix and PCR1 MM

PCR1 Primer Mix

- WILDseq protocol from Kirsty dictates 8 fwd primers mixed in equal amounts with equivalent amount of rev primer
- Primers were reconstituted at 100 μ M so needed to be diluted to 10 μ

Component	Volume
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	16 μ L
DNase/RNase H2O	320 μ L

PCR1 MM

Component	Volume	Master Mix
10uM WS PCR1 Primer Mix	3 μ L	
DNase/RNase H2O	12 μ L	
Kapa Hifi HotStart Ready Mix (2X)	25 μ L	

1. Add following components to tubes

Component	Volume
PCR1 MasterMix	40 μ L
H2O	10 μ L


2. Performed PCR using the following parameters:

Step Name	Steps	Time
Initial Denaturation	Step 1: 95C	3mins
Denaturation	Step 2: 98C	20s
Annealing	Step 3: 60C	15s
Extension	Step 4: 72C	15s
Final Extension	Step 5: 72C	1 min
Hold Step	Step 6: 12C	Hold

10-25 cycles of steps 2-4

Tapestation Control Test 3

- Compared New -Control to old -control
 - Wanted to determine if water was contaminated

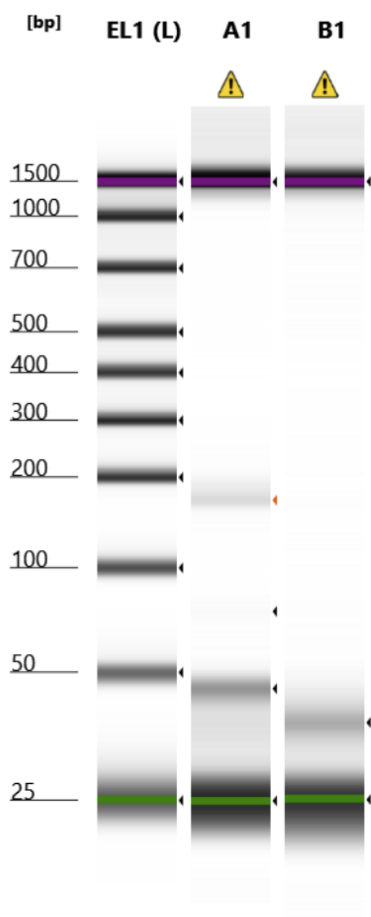

Default image (Contrast 100%)

Sample Info

Well	Conc. [ng/μl]	Sample Description
EL1	20.3	Electronic Ladder
A1	2.13	Control 1
B1	2.01	Control 2

Wednesday 24-01-2024

Cell Culture

A20

- Checked cells ~10% confluent
- Cells looked healthy but sparse
- Grow up for injection
- Ask Chris about cell numbers needed for injecting into 12 mice

HEK

- 20% confluent
- Check tomorrow and split on Friday

WILDseq - Mouse Experiment 1: Library Prep Attempt 5 - RT

- Decided to repeat the entire library prep process using RNA extracted previously
- Everything done in PCR room RNA hood

RT Protocol

1. In PCR strip tubes, prepared 5 μg of RNA in a total volume of 10 μl of RNase/DNAase-free water.
2. Add 1 μl of WS_RT_UMI_NexteraR2 primer (2 μM)
 - Diluted primer from stock (100 μM)
 - Added 1 μL RT-Primer stock into 49 μL water
3. Add 1 μl dNTPs (10 mM each).
4. Denature at 65 C for 5 mins in the PCR machine, then straight onto ice for at least 2 mins, spin briefly to get liquid to the bottom of the tube.
5. Prepare mastermix of RT enzyme and buffers. 1.1x the number of samples.

Component	Volume	MM volume
5x SSIV Buffer	4 μL	134 μL
SSIV RT	1 μL	33.5 μL
100 mM DTT	1 μL	33.5 μL
RNase Out	1 μL	33.5 μL

6. Add 7 μl of RT MM prepared above to each sample and mix, spin briefly to get liquid to the bottom of the tube.
 - Spin briefly to get liquid to bottom of the tube
7. In the PCR machine incubate at 53 C for 10 mins followed by 80 C for 10 mins.
8. Add 1 μl Themolabile Exonuclease I (NEB M0568) to remove excess RT primer.

Component	Volume	MM volume
Themolabile Exonuclease I	1 μL	33.5 μL
NEBuffer r3.1*	2 μL	67 μL

Most PCR buffers are compatible

9. In PCR Machine: heat at 37C for 4 mins followed by 80C for 1 min
 - Spin briefly to get liquid to bottom of the tube
10. Add 1 μl of RNase H and incubate at 37 C for 20 mins.
 - Stored at -4C

Plan for tomorrow

- PCR 1 (check for KAPA HIFI amount)
 - Bead Clean Up
 - Tape Station
- Check on cells
- Plan cell numbers needed for A20 mouse injections

Thursday 25-01-2024

Cell Culture

A20

- Checked cells ~50% confluent
- Cells looked healthy but sparse
- Added 10ml of AR-5
- Grow up for injection
- Ask Chris about cell numbers needed for injecting into 12 mice

HEK

- Collected all T175 and froze down
- Left one flask going at 20% confluency

WS-ME1 Library Prep - Attempt 5 PCR 1

Ran PCR1 for attempt 5

PCR 1: Negative Control

- Remade negative control (MM + H2O)
- Remade Primer mix and PCR1 MM

PCR1 Primer Mix

- WILDseq protocol from Kirsty dictates 8 fwd primers mixed in equal amounts with equivalent amount of rev primer
- Primers were reconstituted at 100 μ M so needed to be diluted to 10 μ M

Component	Volume
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	2 μ L
100uM WS PCR1 Primer Mix	16 μ L
DNase/RNase H2O	144 μ L

PCR1 MM

Component	Volume	Master Mix
10uM WS PCR1 Primer Mix	3 μ L	
DNase/RNase H2O	12 μ L	
Q5 Master Mix	25 μ L	

- Used Q5 instead of KAPA
1. Add following components to tubes

Component	Volume
PCR1 MasterMix	40 μ L
H2O	10 μ L

2. Performed PCR using the following parameters:

Step Name	Steps	Time
Initial Denaturation	Step 1: 95C	3mins
Denaturation	Step 2: 98C	20s
Annealing	Step 3: 60C	15s
Extension	Step 4: 72C	15s
Final Extension	Step 5: 72C	1 min
Hold Step	Step 6: 12C	Hold

10-25 cycles of steps 2-4

Bead Clean Up

Overview

- Clean PCR reaction by removing leftover primers
- This prevents false reads, cleaner indexing in later steps, and higher quality samples

Materials

- For 40 samples

Component	Expected Volume/experiment	Material ID
Ampure Bead Kit	200 μ L	


Protocol

1. Equilibrate Ampure beads to room temp (available from supplies cold room)
2. Mix well and vortex for 30sec to ensure uniform distribution.
3. Add 90ul (1.8X volume) to PCR reaction and mix thoroughly by pipetting.
4. Incubate 5min, room temp
5. Place on magnet for 2 min and keep on magnet until final elution
6. Remove supernatant and add 200ul 70% ethanol (make fresh)
7. Incubate 30sec and remove
8. Repeat 70% ethanol wash
9. Air dry beads 5-10min room temp. Once dry the beads will go from shiny to matt in appearance. Avoid over drying (cracked appearance).
10. Add 15 ul EB buffer and resuspend beads by pipetting
11. Incubate 2min room temp
12. Place on magnet and remove eluate
13. Assess PCR product size, contamination and concentration on tapestation.
Expected size = 165-172bp

Tapestation PCR 1 Samples 1-33

Tapestation PCR 1 1-16

Tapestation PCR 1

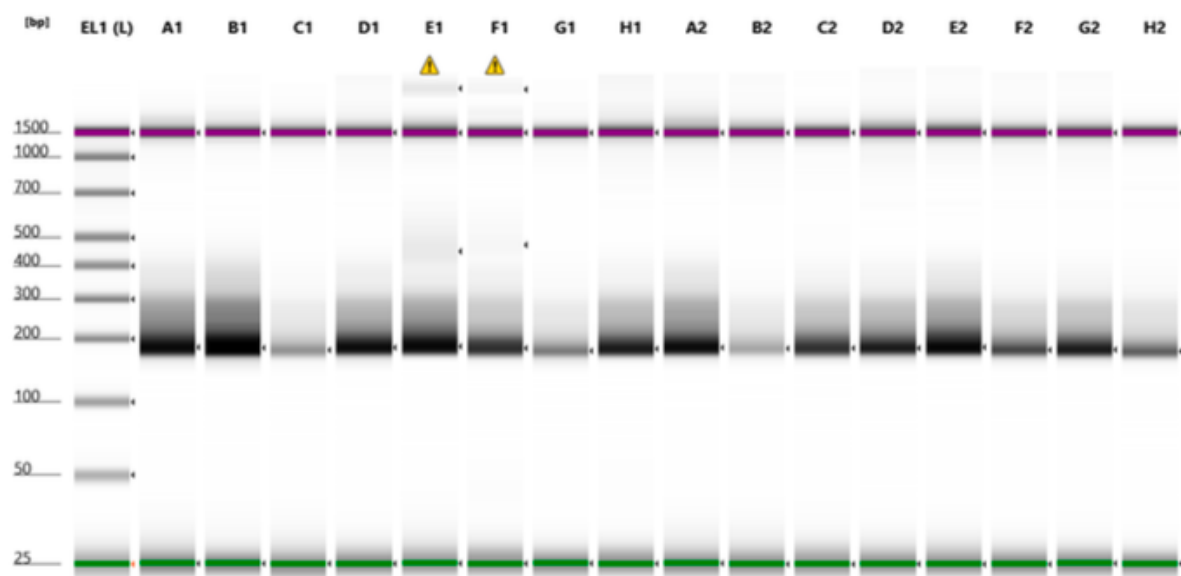


Default image (Contrast 100%)

Sample Info

Well	Conc. [ng/ul]	Sample Description	Alert	Observations
EL1	20.3	Electronic Ladder		Ladder
A1	18.8	1 PCR 1		
B1	28.9	2 PCR 1		
C1	4.49	3 PCR 1		
D1	17.5	4 PCR 1		
E1	19.4	5 PCR 1	⚠	Peak out of Sizing Range
F1	12.7	6 PCR 1	⚠	Peak out of Sizing Range
G1	5.42	7 PCR 1		
H1	13.6	8 PCR 1		
A2	15.8	9 PCR 1		
B2	3.52	10 PCR 1		
C2	12.2	11 PCR 1		
D2	13.0	12 PCR 1		
E2	18.8	13 PCR 1		
F2	8.74	14 PCR 1		
G2	13.2	15 PCR 1		
H2	6.68	16 PCR 1		

Figure 3: 2024-01-25 Tapestation PCR 1 Samples 1-16

Tapestation PCR 1 1-32

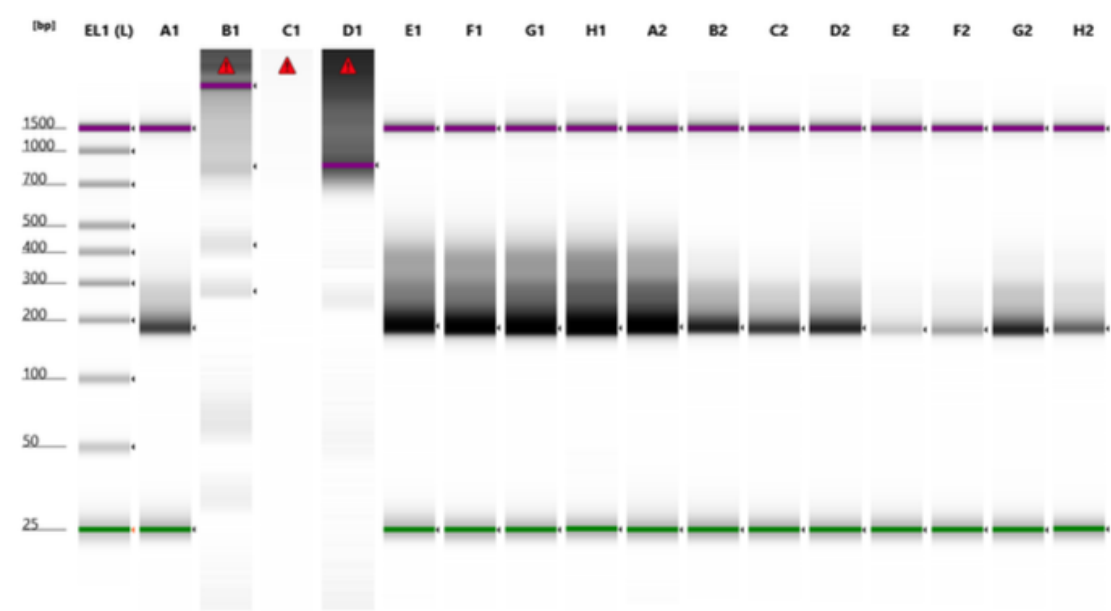

Default image (Contrast 100%)

Sample Info

Well	Conc. [ng/ul]	Sample Description	Alert	Observations
EL1	20.3	Electronic Ladder		Ladder
A1	15.5	17 PCR 1		
B1	2.76	18 PCR 1	▲	Marker(s) not detected
C1		19 PCR 1	▲	Marker(s) not detected
D1		20 PCR 1	▲	Marker(s) not detected
E1	39.1	21 PCR 1		
F1	43.3	22 PCR 1		
G1	46.2	23 PCR 1		
H1	49.7	24 PCR 1		
A2	47.3	25 PCR 1		
B2	20.2	26 PCR 1		
C2	15.7	27 PCR 1		
D2	18.5	28 PCR 1		
E2	2.85	29 PCR 1		
F2	5.29	30 PCR 1		
G2	20.4	31 PCR 1		
H2	10.8	32 PCR 1		

Figure 4: 2024-01-25 TapeStation PCR 1 Samples 17-32


Tapestation PCR 1

- Samples 18, 19, and 20 DNW

Tapestation PCR 1 18-20, 33

Tapestation PCR 1

- Samples 18, 19, and 20 DNW again

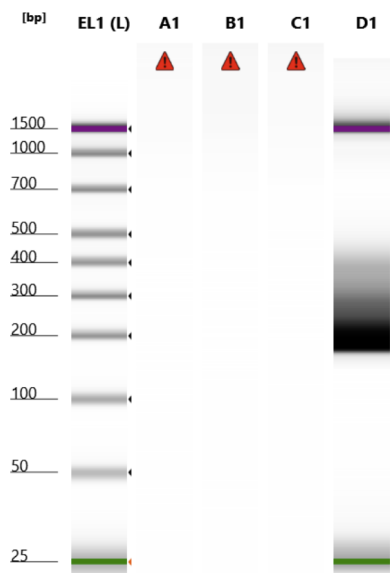

Default image (Contrast 100%)

Sample Info

Well	Conc. [ng/μl]	Sample Description	Alert	Observations
EL1	20.3	Electronic Ladder		Ladder
A1		17 PCR 1 v2	▲	Marker(s) not detected
B1		18 PCR 1 v2	▲	Marker(s) not detected
C1		20 PCR 1 v2	▲	Marker(s) not detected
D1	36.1	33 PCR 1		

Figure 5: 2024-01-25 Tapestation PCR 1 Samples 33 and redo 18-20

Friday 26-01-2024

Cell Culture

A20

- Expanded to T175
- Grow up for injection
- Ask Chris about cell numbers needed for injecting into 12 mice

HEK

- Split 1/2
- DMEM

WS-ME1 Library Prep - Attempt 5 PCR 1 repeat and PCR 2 prep

- Reran PCR1 for samples 18,19,20
- Prepped PCR2 primers

PCR2 Prep

- Reconstituted PCR2 primers bought by Chris
 - Used DNase/RNase Free H2O and diluted to 100 μ M
- Diluted into strip tubes working stock (10 μ M)
 - 90 μ L H2O + 10 μ L primer stock

PCR1 Samples 18-20

PCR1 Protocol

- Used Primer Mix previously made on 25-01-2024

Component	Volume	Master Mix
10uM WS PCR1 Primer Mix	3 μ L	
DNase/RNase H2O	12 μ L	
Q5 Master Mix	25 μ L	

- Used Q5 instead of KAPA

1. Add following components to tubes

Component	Volume
PCR1 MasterMix	40 μ L
H2O	10 μ L

2. Performed PCR using the following parameters:

Step Name	Steps	Time
Initial Denaturation	Step 1: 95C	3mins
Denaturation	Step 2: 98C	20s
Annealing	Step 3: 60C	15s
Extension	Step 4: 72C	15s
Final Extension	Step 5: 72C	1 min

Step Name	Steps	Time
Hold Step	Step 6: 12C	Hold

10-25 cycles of steps 2-4

Bead Clean Up

Overview

- Clean PCR reaction by removing leftover primers
- This prevents false reads, cleaner indexing in later steps, and higher quality samples

Materials

- For 40 samples

Component	Expected Volume/experiment	Material ID
Ampure Bead Kit	200 μ L	

Protocol

1. Equilibrate Ampure beads to room temp (available from supplies cold room)
2. Mix well and vortex for 30sec to ensure uniform distribution.
3. Add 90ul (1.8X volume) to PCR reaction and mix thoroughly by pipetting.
4. Incubate 5min, room temp
5. Place on magnet for 2 min and keep on magnet until final elution
6. Remove supernatant and add 200ul 70% ethanol (make fresh)
7. Incubate 30sec and remove
8. Repeat 70% ethanol wash
9. Air dry beads 5-10min room temp. Once dry the beads will go from shiny to matt in appearance. Avoid over drying (cracked appearance).
10. Add 15 ul EB buffer and resuspend beads by pipetting
11. Incubate 2min room temp
12. Place on magnet and remove eluate
13. Assess PCR product size, contamination and concentration on tapestation.
Expected size = 165-172bp

Tapestation PCR 1 Samples 18-20

Tapestation PCR 1 18-20


Tapestation PCR 1

Monday 29-01-2024

Cell Culture

A20

- Split 1:6
- AR-5
- Ask Chris about cell numbers needed for injecting into 12 mice

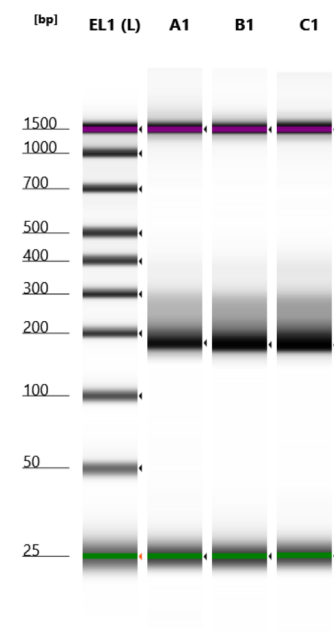

Default image (Contrast 100%)

Sample Info

Well	Conc. [ng/ul]	Sample Description	Alert	Observations
EL1	20.3	Electronic Ladder		Ladder
A1	8.30	18 PCR 1		
B1	9.76	19 PCR 1		
C1	10.3	20 PCR 1		

Figure 6: 2024-01-26 Tapestation PCR 1 Samples 18-20

HEK

- Split 1/6
- DMEM

PCR Step 2 - Samples 1-24

Overview

- PCR to index samples for sequencing by attaching UMI
- Didn't have enough Q5 to complete all the samples

PCR2 Protocol

1. Prepare samples according to the following table:
 - Keep careful note of which adaptors are added to which sample

Component	Volume
10uM Nxxx Nextera i7 adapter	1.5 μ L
10uM Sxxx Nextera i5 adapter	1.5 μ L
DNAse/RNAse H20	21 μ L
Q5	25 μ L
10ng/ μ L PCR1	1 μ L


2. Perform PCR using the following parameters:

Step Name	Steps	Time
Initial Denaturation	Step 1: 95C	3mins
Denaturation	Step 2: 98C	20s
Annealing	Step 3: 55C	15s
Extension	Step 4: 72C	15s
Final Extension	Step 5: 72C	2 min
Hold Step	Step 6: 12C	Hold

8 cycles of steps 2-4

3. Perform Ampure bead clean up as above with 90ul beads per 50ul PCR reaction.
 - If you don't have a crazy number of samples you might also want to consider doing up to 4 PCR2 reactions per sample and pooling them afterwards just to make sure all your sequences are captured.
4. Check size and concentration on tapestation. Dilute and pool samples for sequencing.

Tapestation

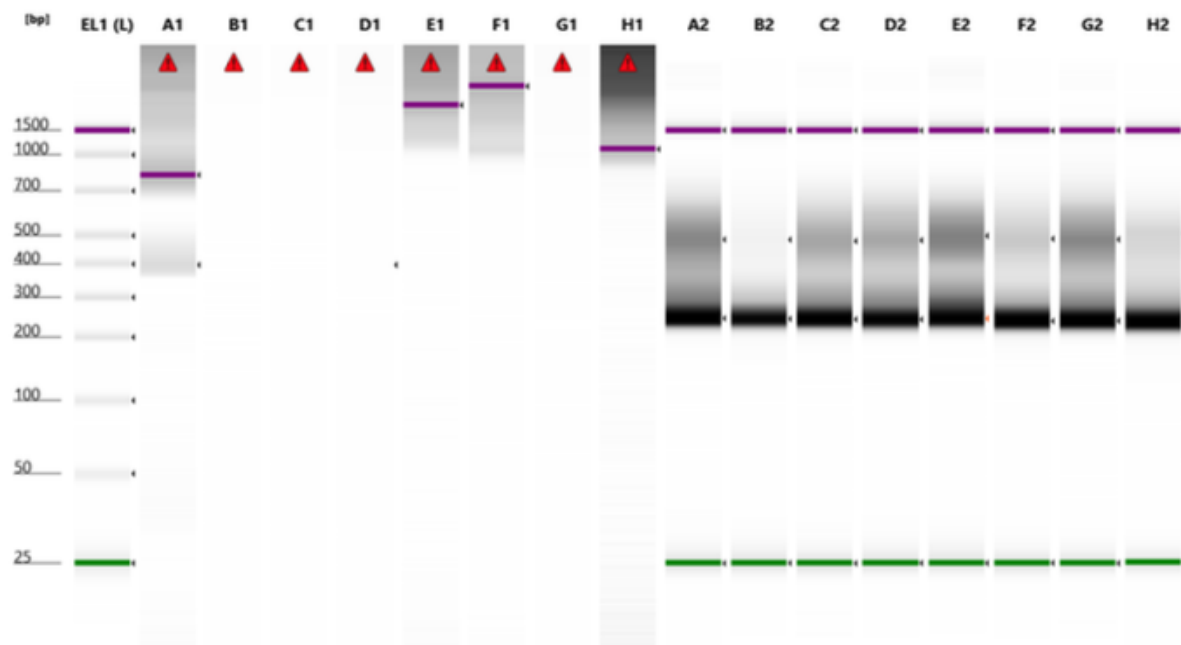


Default image (Contrast 100%)

Sample Info

Well	Conc. [ng/ul]	Sample Description	Alert	Observations
EL1	20.3	Electronic Ladder		Ladder
A1	2.72	1 PCR 2	▲	Marker(s) not detected
B1		2 PCR 2	▲	Marker(s) not detected
C1		3 PCR 2	▲	Marker(s) not detected
D1		4 PCR 2	▲	Marker(s) not detected
E1		5 PCR 2	▲	Marker(s) not detected
F1		6 PCR 2	▲	Marker(s) not detected
G1		7 PCR 2	▲	Marker(s) not detected
H1		8 PCR 2	▲	Marker(s) not detected
A2	206	9 PCR 2		
B2	88.9	10 PCR 2		
C2	173	11 PCR 2		
D2	159	12 PCR 2		
E2	160	13 PCR 2		
F2	164	14 PCR 2		
G2	214	15 PCR 2		
H2	149	16 PCR 2		

- 1-8 PCR DNW (suspected issue with bead clean-up)

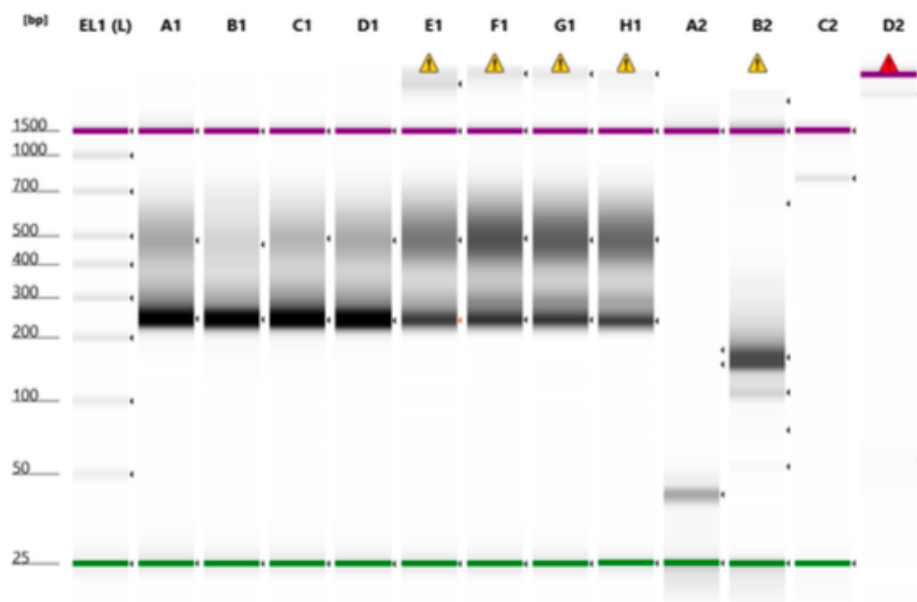

Default image (Contrast 100%)

Sample Info

Well	Conc. [ng/ul]	Sample Description	Alert	Observations
EL1	20.3	Electronic Ladder		Ladder
A1	167	17 PCR 2		
B1	155	18 PCR 2		
C1	186	19 PCR 2		
D1	200	20 PCR 2		
E1	177	21 PCR 2	⚠	Peak out of Sizing Range
F1	224	22 PCR 2	⚠	Peak out of Sizing Range
G1	210	23 PCR 2	⚠	Peak out of Sizing Range
H1	184	24 PCR 2	⚠	Peak out of Sizing Range
A2	13.5			
B2	37.3		⚠	Peak out of Sizing Range
C2	2.29			
D2	3.24		⚠	Marker(s) not detected

- Last 4 samples are from Chris (seperate experiment)

Tuesday 30-01-2024

Cell Culture

HEK

- Unhealthy
- Discarded

A20

- Split 1/6

February 2024

Thursday 01-02-2024

PCR Step 2 - Samples 1-8, 25-39

- Repeated PCR2 for samples 1-8

Overview

- PCR to index samples for sequencing by attaching UMI
- Didn't have enough Q5 to complete all the samples

PCR2 Protocol

1. Prepare samples according to the following table:
 - Keep careful note of which adaptors are added to which sample

Component	Volume
10uM Nxxx Nextera i7 adapter	1.5 μ L
10uM Sxxx Nextera i5 adapter	1.5 μ L
DNase/RNase H20	21 μ L
Q5	25 μ L
10ng/ μ L PCR1	1 μ L


2. Perform PCR using the following parameters:

Step Name	Steps	Time
Initial Denaturation	Step 1: 95C	3mins
Denaturation	Step 2: 98C	20s
Annealing	Step 3: 55C	15s
Extension	Step 4: 72C	15s
Final Extension	Step 5: 72C	2 min
Hold Step	Step 6: 12C	Hold

8 cycles of steps 2-4

3. Perform Ampure bead clean up as above with 90ul beads per 50ul PCR reaction.
 - If you don't have a crazy number of samples you might also want to consider doing up to 4 PCR2 reactions per sample and pooling them afterwards just to make sure all your sequences are captured.
4. Check size and concentration on tapestation. Dilute and pool samples for sequencing.


Tapestation

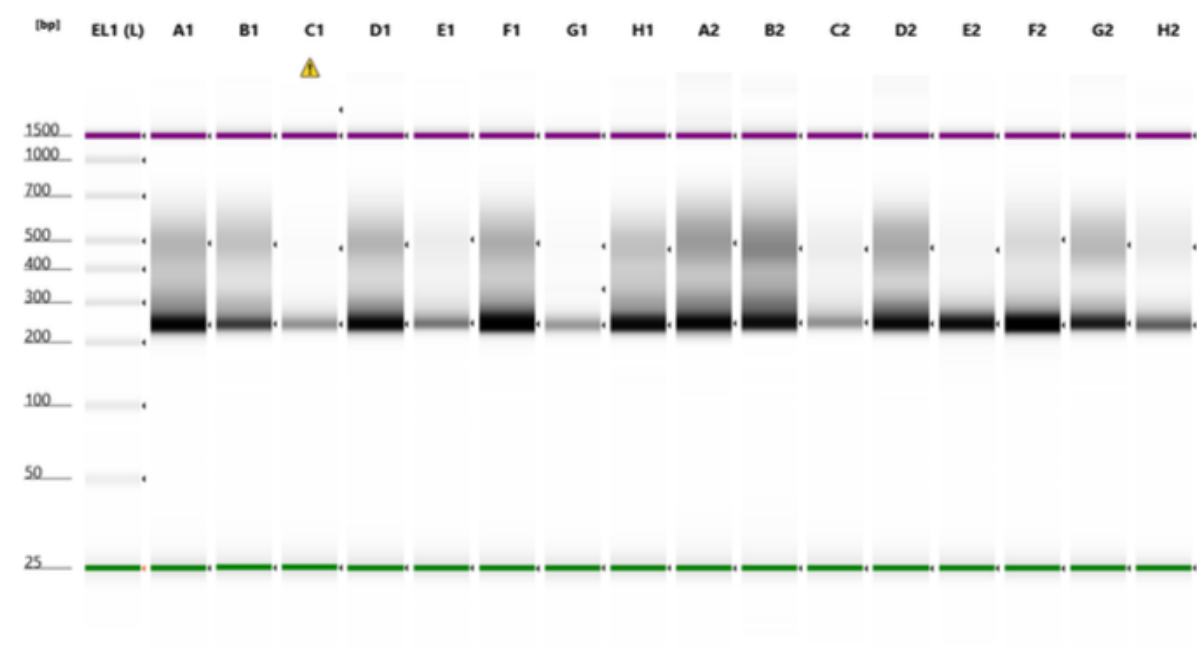

Default image (Contrast 100%)

Sample Info

Well	Conc. [ng/ul]	Sample Description	Alert	Observations
EL1	20.3	Electronic Ladder		Ladder
A1	133	1 PCR 2		
B1	83.8	2 PCR 2		
C1	20.0	3 PCR 2	⚠	Peak out of Sizing Range
D1	130	4 PCR 2		
E1	34.0	5 PCR 2		
F1	166	6 PCR 2		
G1	21.2	7 PCR 2		
H1	94.9	8 PCR 2		
A2	153	25 PCR 2		
B2	147	26 PCR 2		
C2	29.1	27 PCR 2		
D2	149	28 PCR 2		
E2	81.6	29 PCR 2		
F2	141	30 PCR 2		
G2	104	31 PCR 2		
H2	45.9	32 PCR 2		

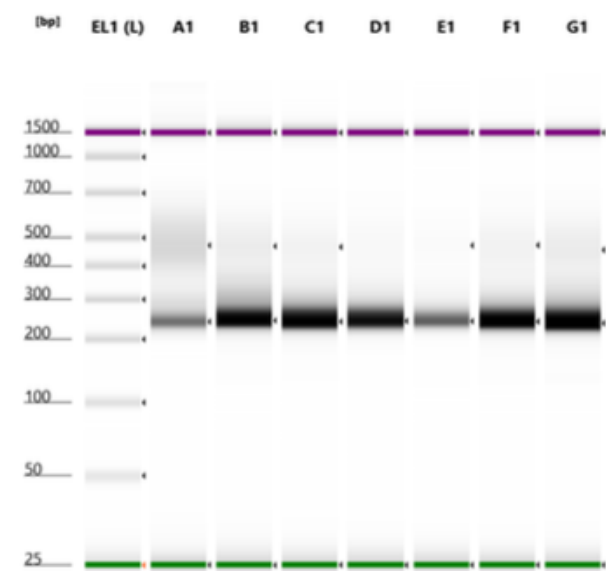

Monday 05-02-2024

Cell Culture

Made DMEM

Solution	ID code	Volume	% Total volume
DMEM	—	500 mL	78%
Pen-Strep	—	5.5 mL	1%
FBS	—	55 mL	10%

HEK - Seeded

- Seeded New HEK
 - HEK stock are in CM Box 1 -80
1. Thawed vial
 2. Added 1mL DMEM to vial
 3. Transferred to 15mL falcon tube
 4. Slowly added 5mL DMEM
 5. Spun down at 650rpm for 4min
 6. Removed supernant and resuspended vial flicking
 7. Added 1mL DMEM and pipetted multiple times to create single cell suspension
 8. Added to T25 flask
 9. Added 14mL of DMEM
 10. Placed in incubator

A20 - Split

- Split 1/6
- Split T25 vial for Emily James (EJ)

Tuesday 06-02-2024

WS-ME1 Pooling for Submission

- Genomics core indicated that library should be submitted at 5nM-10nM in a volume of 20-30uL
- Cubit PCR2 Samples to determine concentration

Sample ID	Submission ID	Nextera Adaptor Indices	Desired ng/uL	Base Pair	[Final] nM	Sample Volume (uL)	H2O (uL)
1	1_Cyclophosphamidei701-i502		1.5	241	9.43	1	29.60
2	2_Cyclophosphamidei702-i502		1.5	241	9.43	1	19.20
3	3_Cyclophosphamidei703-i502		1.5	241	9.43	1	4.52
4	4_Cyclophosphamidei704-i502		1.5	241	9.43	1	25.90
5	5_Cyclophosphamidei705-i502		1.5	241	9.43	1	6.40
6	6_Combination	i706-i502	1.5	241	9.43	1	31.90
7	7_Combination	i707-i502	1.5	241	9.43	1	4.50
8	8_Combination	i710-i502	1.5	241	9.43	1	24.70
9	9_Combination	i701-i503	1.5	241	9.43	1	45.70
10	10_Combination	i702-i503	1.5	241	9.43	1	22.30
11	11_Methotrexate	i703-i503	1.5	241	9.43	1	43.90
12	12_Methotrexate	i704-i503	1.5	241	9.43	1	32.70

Sample ID	Submission ID	Nextera Adaptor Indices	Desired ng/uL	Base Pair	[Final] nM	Sample Volume (uL)	H2O (uL)
13	13_Methotrexate	i705-i503	1.5	241	9.43	1	41.20
14	14_Methotrexate	i706-i503	1.5	241	9.43	1	34.80
15	15_Methotrexate	i707-i503	1.5	241	9.43	1	46.10
16	16_Vehicle	i710-i503	1.5	241	9.43	1	28.30
17	17_Vehicle	i701-i505	1.5	241	9.43	1	39.70
18	18_Vehicle	i702-i505	1.5	241	9.43	1	32.00
19	19_Vehicle	i703-i505	1.5	241	9.43	1	34.90
20	20_Vehicle	i704-i505	1.5	241	9.43	1	34.00
21	21_Baseline	i705-i505	1.5	241	9.43	1	34.10
22	22_Baseline	i706-i505	1.5	241	9.43	1	34.30
23	23_Baseline	i707-i505	1.5	241	9.43	1	28.90
24	24_Baseline	i710-i505	1.5	241	9.43	1	31.60
25	25_Baseline	i701-i506	1.5	241	9.43	1	27.60
26	26_BLLW 14K	i702-i506	1.5	241	9.43	1	24.10
27	27_BLLW 2K	i703-i506	1.5	241	9.43	1	5.20
28	28_BLLW 1K	i704-i506	1.5	241	9.43	1	30.50
29	29_Methotrexate (outlier)	i705-i506	1.5	241	9.43	1	14.00
30	30_Combo	i706-i506	1.5	241	9.43	1	28.10
31	31_RBL2P_2K	i707-i506	1.5	241	9.43	1	15.10
32	32_RBL2P_7K	i710-i506	1.5	241	9.43	1	49.40
33	33_RBL2P_250K	i701-i507	1.5	241	9.43	1	41.50
34	34_Mock_direct_1	i702-i507	1.5	241	9.43	1	15.30
35	35_Mock_direct_2	i703-i507	1.5	241	9.43	1	15.90
36	36_Mock_direct_3	i704-i507	1.5	241	9.43	1	10.90
37	37_Mock_culture_1	i705-i507	1.5	241	9.43	1	4.48
38	38_Mock_culture_2	i706-i507	1.5	241	9.43	1	18.00
39	39_Mock_culture_3	i707-i507	1.5	241	9.43	1	17.70

- Used IDT Library Concentration Conversion Calculator
 - Link: <https://eu.idtdna.com/Calc/library-concentration-conversion>

Wednesday 07-02-2024

Cell Culture

HEK

- Decided to wait till Monday to start making virus
- Discarded HEK
- Will reseed HEK Sunday

A20

- Split 1/6
- Keeping going for NSG/BALBc injection
 - Issue with procedure room in AMB, waiting for that to be resolved
- Also plan to incorporate WILDseq

Friday 09-02-2024

Cell Culture

A20

- Split 1/6
- Keeping going for NSG/BALBc injection
 - Issue with procedure room in AMB, waiting for that to be resolved
- Also plan to incorporate WILDseq

Monday 12-02-2024

Cell Culture

Made DMEM-V1

- DMEM for virus production
- Does not have Pen/Strep

Solution	ID code	Volume	% Total volume
DMEM	—	500 mL	78%
FBS	—	128 mL	20%
Glutamax	—	6.5 mL	1%

Thawed HEK

A20

- Split 1/6
- Keeping going for NSG/BALBc injection
 - Planning to inject with Chris tomorrow
- Also plan to incorporate WILDseq once virus is made

Processing scRNAseq

Tuesday 13-02-2024

Cell Culture

HEK

- Morphology looked weird
- 20% FBS is probably too high
- Threw away will start next week

Processing scRNAseq

Wednesday 14-02-2024

Cell Culture

A20

- Split 1/6
- Keeping going for NSG/BALBc injection

- Planning to inject with Chris tomorrow
- Also plan to incorporate WILDseq once virus is made

Processing scRNAseq

Friday 16-02-2024

A20 Xenograft EXP 1

Cell Prep

1. Thawed 2.5mL of Matrigel on ice ~3hours before
2. Transferred 12×10^6 cells to a 15mL eppendorf
3. Pelleted CS
4. Resuspended in 2.5mL of PBS
5. Transferred 400 μ L CS to 6 different 1mL eppendorfs
6. Added 400 μ L of Matrigel, mixed gently and placed back on ice
7. Brought materials over to the AMB for injection

Injection

- Chris Injected
- Study Plan: SP140164

Mouse ID	Earmark	Genotype	Injection Location
TUAD36.2h	NM	NSG	IP
TUAD36.2i	1R	NSG	IP
TUAD36.2k	2R	NSG	IP
TUAD36.2a	NM	NSG	SC
TUAD36.2b	1L	NSG	SC
TUAD36.2c	1R	NSG	SC
TUAP3.1e	1L	BALB/c	IP
TUAP3.1a	NM	BALB/c	IP
TUAP4.1a	2L	BALB/c	IP
TUAP3.1b	1L	BALB/c	SC
TUAP4.1b	2R	BALB/c	SC
TUAP3.1c	1R	BALB/c	SC

- Each mouse received $\sim 1 \times 10^6$ cells
- Mice will be monitored over the coming weeks for tumour development

Monday 19-02-2024

Cell Culture

A20

- Split 1/6
- Used AR-5

Wednesday 21-02-2024

Cell Culture

A20

- Split 1/6
- Used AR-5

HEK

- Seeded New HEK
 - HEK stock are in CM Box 1 -80
 - Used DMEM-V1
1. Thawed vial
 2. Added 1mL DMEM to vial
 3. Transferred to 15mL falcon tube
 4. Slowly added 5mL DMEM
 5. Spun down at 650rpm for 4min
 6. Removed supernant and resuspended vial flicking
 7. Added 1mL DMEM and pipetted multiple times to create single cell suspension
 8. Added to T25 flask
 9. Added 14mL of DMEM
 10. Placed in incubator

Friday 23-02-2024

Cell Culture

A20

- Split 1/6
- Used AR-5

HEK

- Split 1/6
- Used DMEM-V1

Monday 26-02-2024

Cell Culture

A20

- Split into T25
- Froze down 4 vials
- 4×10^6 cells per vial
- Used AR-5

HEK

- Froze down 7 vials
- 6×10^6 cells per vial

WILDseq Virus Production Day 1

- Seeded 2 10cm dishes of HEK cells
 - Coated with Poly-D-Lysine (covered plate, removed liquid, let sit for 10 min, washed with PBS)
 - 4.5×10^6 cells seeded per dish (seeded at 13:00)

March 2024

Wednesday 06-03-2024

- Collected SC tumours
- HEK infection
 - binned A20/HEK

Wednesday 13-03-2024

- Made media
- Seeded in N2-BC4
- Froze down A20
- Collected 3 IP NSG tumours

Thursday 14-03-2024

- Changed media (N2-BC4)

Friday 15-03-2024

Cell Culture

N2-BC4 RTX CDC Assay - version 1

Monday 18-03-2024

Cell Culture - N2 -C4

- Split: seeded 5×10^6 cells

N2-BC4 RTX CDC Assay - version 2

Overview: Trying to optimise RTX CDC assay - Using 0% serum, since there are no complement factors in the media the cell live/dead reading should not change - Next step: add varying amounts of human serum

Friday 22-03-2024

Cell Culture

N2-BC4

- Cells look healthy
- Feeders in flask are also healthy
- Split 5×10^6 back into T175 with 30mL of AR-6

EC50 RTX N2-BC4 22324 - Seeding

- Seeded 96-well plate with N2-BC4 and treated with RTX concentration range with or without 10% Human Serum (HS)
 - Used Heat-Inactivated Serum which does not have functional complement which is why this experiment did not work
- Plate seeding protocol:
 - Diluted cell suspension to seed 20000 cells/well in 50 μ L amounts


Well Number	Cells total	Volume total	CS cells/mL	Flask cell count	Flask Volume	Media Volume
60 wells	1.2x10 ⁶ cells	3 mL	4x10 ⁵	3.06x10 ⁶	400 μ L	2.6 mL

- Made RTX dilutions and added to respective wells in 50 μ L
 - [RTX stock] = 10.3 mg/mL - Drug volumes are added in triplicate - Drug volumes are being added constitute 1/4 of well volume: [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 50 μ L per well ~ minimum of 300 μ L per condition needed (recommend 500 μ L)

Dilution ID	[RTX Well]	Media Volume	Source	[RTX WS]
RTX 1	100 μ g/mL	1 mL	50 μ L stock	400 μ g/mL
RTX 2	50 μ g/mL	500 mL	500 μ L RTX 1	200 μ g/mL
RTX 3	25 μ g/mL	500 mL	500 μ L RTX 2	100 μ g/mL
RTX 4	12.5 μ g/mL	500 mL	500 μ L RTX 3	50 μ g/mL
RTX 5	6.25 μ g/mL	500 mL	500 μ L RTX 4	25 μ g/mL
RTX 6	3.125 μ g/mL	500 mL	500 μ L RTX 5	12.5 μ g/mL
RTX 7	1.56 μ g/mL	500 mL	500 μ L RTX 6	6.25 μ g/mL
RTX 8	.78 μ g/mL	500 mL	500 μ L RTX 7	3.125 μ g/mL
RTX 9	.39 μ g/mL	500 mL	500 μ L RTX 8	1.56 μ g/mL
CNTR	0 μ g/mL	500 mL	-	0 μ g/mL

- Added 10% serum or media control
 - This provides complement factors to bind BL-bound RTX and initiate CDC
 - Serum/media amounts added at 100 μ L/well
 - Serum volume is added in 1:4 ratio (20 μ L serum in 200 μ L final well volume)
 - Serum stock mix is 1.2mL HS : 2.4 mL media / plate
- Plate is incubated for 72 hrs at 37C

EC50_RTX_N2-BC4_250322

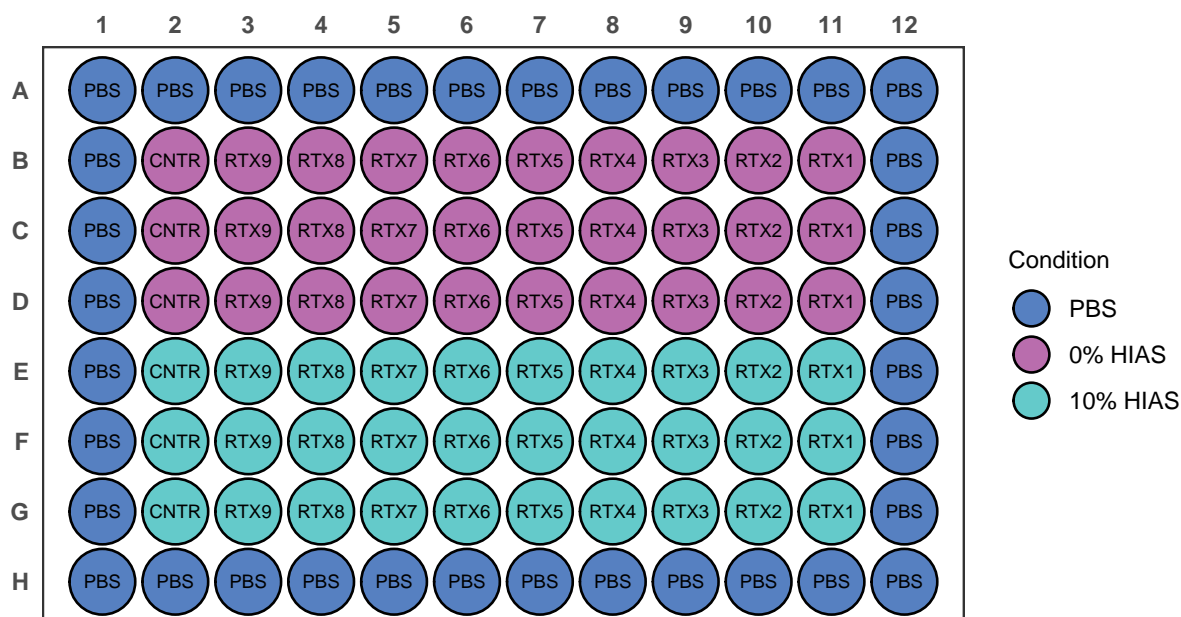

Monday 25-03-2024

EC50 RTX N2-BC4 22324 - Collection

- Collected plate seeded on 22-03-2025
- EC Plate collection protocol:
 - Added 40 μ L Cell Titre Blue (CTB) to each conditioned well
 - 20 μ L CTB/100 μ L of conditioned well recommended by manufacturer
 - Incubated for 1hr at 37C
 - Read on plate reader according to Cell Titre Blue Protocol
- Results: DNW

EC50 RTX N2-BC4 25324 - Seeding (DNW)

- Seeded an EC50 experiment comparing the effects of RTX on N2-BC4 with or without 10% serum
 - Used Heat-Inactivated Serum which does not have functional complement which is why this experiment did not work
- Decreased overall well volume to reduce CTB being used
- Used RTX provided by Jamie
- Plate seeding protocol:
 - Diluted cell suspension to seed 20000 cells/well in 25 μ L amounts

Well Number	Cells total	Volume total	CS cells/mL	Flask cell count	Flask Volume	Media Volume
60 wells	1.2x10 ⁶ cells	1.5 mL	8x10 ⁵	2.92x10 ⁶	410 μ L	1.1 mL

- Made RTX dilutions and added to respective wells in 50 μ L
 - [RTX stock] = 10.3 mg/mL - Drug volumes are added in triplicate - Drug volumes are being added constitute 1/4 of well volume: [RTX working] needs to be 4x [RTX well]

- 6 wells per condition, 25 μ L per well ~ minimum of 150 μ L per condition needed (recommend 200 μ L)


Dilution ID	[RTX Well]	Media Volume	Source	[RTX WS]
RTX 1	100 μ g/mL	400 μ L	15.5 μ L stock	400 μ g/mL
RTX 2	50 μ g/mL	200 mL	200 μ L RTX 1	200 μ g/mL
RTX 3	25 μ g/mL	200 mL	200 μ L RTX 2	100 μ g/mL
RTX 4	12.5 μ g/mL	200 mL	200 μ L RTX 3	50 μ g/mL
RTX 5	6.25 μ g/mL	200 mL	200 μ L RTX 4	25 μ g/mL
RTX 6	3.125 μ g/mL	200 mL	200 μ L RTX 5	12.5 μ g/mL
RTX 7	1.56 μ g/mL	200 mL	200 μ L RTX 6	6.25 μ g/mL
RTX 8	.78 μ g/mL	200 mL	200 μ L RTX 7	3.125 μ g/mL
RTX 9	.39 μ g/mL	200 mL	200 μ L RTX 8	1.56 μ g/mL
CNTR	0 μ g/mL	500 mL	-	0 μ g/mL

3. Added 10% serum or media control

- This provides complement factors to bind BL-bound RTX and initiate CDC
- Serum/media amounts added at 50 μ L/well
- Serum volume is added in 1:4 ratio (10 μ L serum in 100 μ L final well volume)
- Serum stock mix is 600 μ L HS : 1.2 mL media / plate

4. Plate is incubated for 72 hrs at 37C

EC50_RTX_N2-BC4_25324

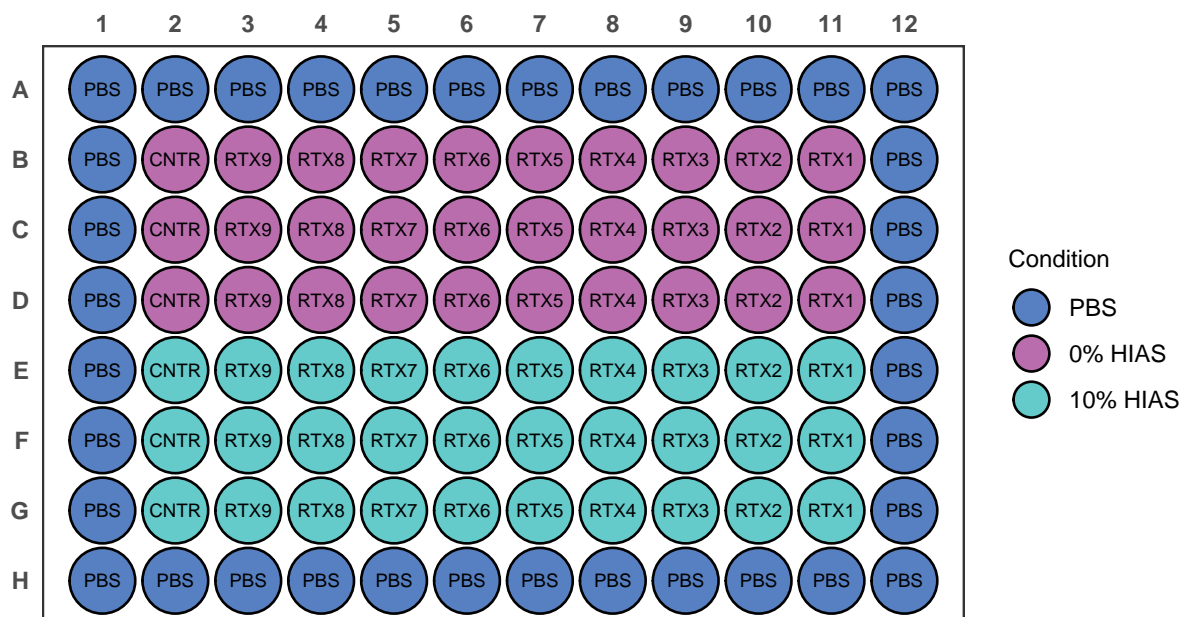

Cell Culture

N2-BC4

- Cells look healthy
- Feeders in flask are also healthy
- Split 5x10⁶ back into T175 with 30mL of AR-6

Wednesday 27-03-2024

Thursday 28-03-2024

Friday 29-03-2024

EC50 RTX N2-BC4/RAMOS-BC5 25324 - Seeding (DNW)

- Seeded an EC50 experiment comparing the effects of RTX on N2-BC4 and RAMOS-BC5 with or without 10% serum
 - Used Heat-Inactivated Serum which does not have functional complement which is why this experiment did not work
- Decreased overall well volume to reduce CTB being used
- Used RTX provided by Jamie
- Plate seeding protocol:
 1. Diluted cell suspension to seed 20000 cells/well in 25 μ L amounts


Cell Line	Well Number	Cells total	Volume total	CS cells/mL	Flask cell count	Flask Volume
N2-BC4	60 wells	1.2x10 ⁶ cells	1.5 mL	8x10 ⁵	2.92x10 ⁶	410 μ L
RAMOS-BC5	60 wells	1.2x10 ⁶ cells	1.5 mL	8x10 ⁵	2.92x10 ⁶	410 μ L

2. Made RTX dilutions and added to respective wells in 50 μ L
 - [RTX stock] = 10.3 mg/mL - Drug volumes are added in triplicate - Drug volumes are being added constitute 1/4 of well volume: [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 25 μ L per well ~ minimum of 150 μ L per condition needed (recommend 200 μ L)

Dilution ID	[RTX Well]	Media Volume	Source	[RTX WS]
RTX 1	100 μ g/mL	400 μ L	15.5 μ L stock	400 μ g/mL
RTX 2	50 μ g/mL	200 mL	200 μ L RTX 1	200 μ g/mL
RTX 3	25 μ g/mL	200 mL	200 μ L RTX 2	100 μ g/mL
RTX 4	12.5 μ g/mL	200 mL	200 μ L RTX 3	50 μ g/mL
RTX 5	6.25 μ g/mL	200 mL	200 μ L RTX 4	25 μ g/mL
RTX 6	3.125 μ g/mL	200 mL	200 μ L RTX 5	12.5 μ g/mL
RTX 7	1.56 μ g/mL	200 mL	200 μ L RTX 6	6.25 μ g/mL
RTX 8	.78 μ g/mL	200 mL	200 μ L RTX 7	3.125 μ g/mL
RTX 9	.39 μ g/mL	200 mL	200 μ L RTX 8	1.56 μ g/mL
CNTR	0 μ g/mL	500 mL	-	0 μ g/mL

3. Added 10% serum or media control
 - This provides complement factors to bind BL-bound RTX and initiate CDC
 - Serum/media amounts added at 50 μ L/well
 - Serum volume is added in 1:4 ratio (10 μ L serum in 100 μ L final well volume)
 - Serum stock mix is 600 μ L HS : 1.2 mL media / plate
4. Plate is incubated for 72 hrs at 37C

EC50_RTX_N2-BC4_250324

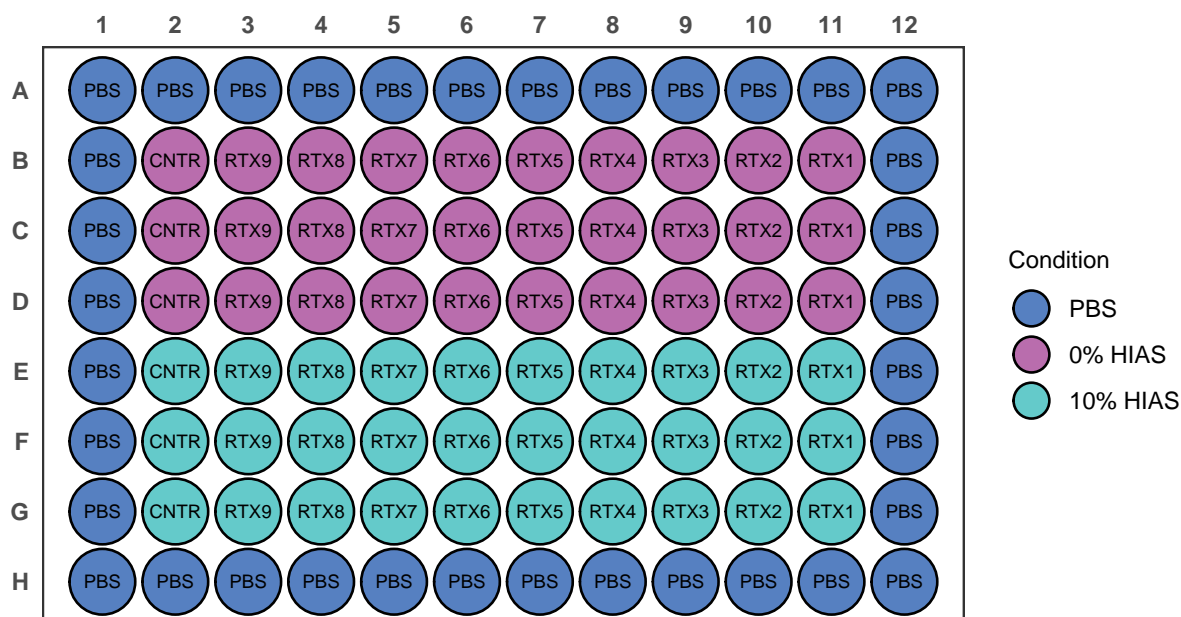

Sunday 31-03-2024

Cell Culture

N2-BC4

- Healthy, split
- Changed FDC:
 1. Thawed vial
 2. Added 1mL DMEM to vial
 3. Transferred to 15mL falcon tube
 4. Slowly added 5mL DMEM
 5. Spun down at 650rpm for 4min
 6. Removed supernant and resuspended vial flicking
 7. Added 1mL DMEM and pipetted multiple times to create single cell suspension
 8. Added to T175 flask
 9. Added 30mL of DMEM
 10. Added 5×10^5 N2 cells from previous flask
 11. Placed in incubator

Ramos-BC5


- Healthy, split
- Count: 2.18×10^6 cells/mL
- Split: 7×10^5 cells in 15 mL

Ramos CDC Testing

- Tested CDC assay adapted from [Ge et al., 2019](#)
- CDC protocol
 1. Seed 250k cells 500uL volumes per well of a 12-well plate

2. Add RTX or media in 250uL volumes per well
 - [RTX] working needs to be 4x final desired well volume
3. Add Serum or media in 250uL volumes per well
 - Serum amount should be 4x final % desired

EC50_RTX_N2-BC4_250331

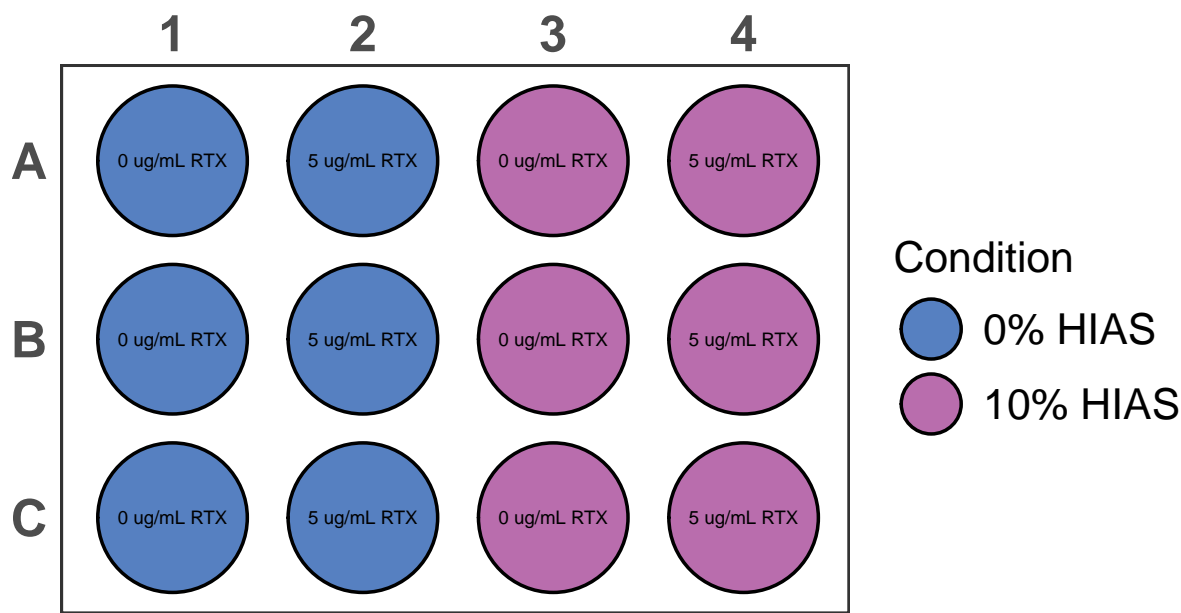

EC50 RTX N2-BC4/RAMOS-BC5 24331 - Seeding (DNW)

- Seeded an EC50 experiment comparing the effects of RTX on N2-BC4 and RAMOS-BC5 with or without 10% serum
- Decreased overall well volume to reduce CTB being used
- Used RTX provided by Jamie
- Only incubated for 24hr
- Plate seeding protocol:
 1. Diluted cell suspension to seed 20000 cells/well in 25μL amounts

Cell Line	Well Number	Cells total	Volume total	CS cells/mL	Flask cell count	Flask Volume
N2-BC4	60 wells	1.2x10 ⁶ cells	1.5 mL	8x10 ⁵	2.92x10 ⁶	410μL
RAMOS-BC5	60 wells	1.2x10 ⁶ cells	1.5 mL	8x10 ⁵	2.92x10 ⁶	410μL

2. Made RTX dilutions and added to respective wells in 50μL
 - [RTX stock] = 10.3 mg/mL - Drug volumes are added in triplicate - Drug volumes are being added constitute 1/4 of well volume: [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 25μL per well ~ minimum of 150μL per condition needed (recommend 200μL)

Dilution ID	[RTX Well]	Media Volume	Source	[RTX WS]
RTX 1	100 μg/mL	400μL	15.5μL stock	400 μg/mL


Dilution ID	[RTX Well]	Media Volume	Source	[RTX WS]
RTX 2	50 $\mu\text{g/mL}$	200 mL	200 μL RTX 1	200 $\mu\text{g/mL}$
RTX 3	25 $\mu\text{g/mL}$	200 mL	200 μL RTX 2	100 $\mu\text{g/mL}$
RTX 4	12.5 $\mu\text{g/mL}$	200 mL	200 μL RTX 3	50 $\mu\text{g/mL}$
RTX 5	6.25 $\mu\text{g/mL}$	200 mL	200 μL RTX 4	25 $\mu\text{g/mL}$
RTX 6	3.125 $\mu\text{g/mL}$	200 mL	200 μL RTX 5	12.5 $\mu\text{g/mL}$
RTX 7	1.56 $\mu\text{g/mL}$	200 mL	200 μL RTX 6	6.25 $\mu\text{g/mL}$
RTX 8	.78 $\mu\text{g/mL}$	200 mL	200 μL RTX 7	3.125 $\mu\text{g/mL}$
RTX 9	.39 $\mu\text{g/mL}$	200 mL	200 μL RTX 8	1.56 $\mu\text{g/mL}$
CNTR	0 $\mu\text{g/mL}$	500 mL	-	0 $\mu\text{g/mL}$

3. Added 10% serum or media control

- This provides complement factors to bind BL-bound RTX and initiate CDC
- Serum/media amounts added at 50 μL /well
- Serum volume is added in 1:4 ratio (10 μL serum in 100 μL final well volume)
- Serum stock mix is 600 μL HS : 1.2 mL media / plate

4. Plate is incubated for 72 hrs at 37C

EC50_RTX_N2-BC4_250331

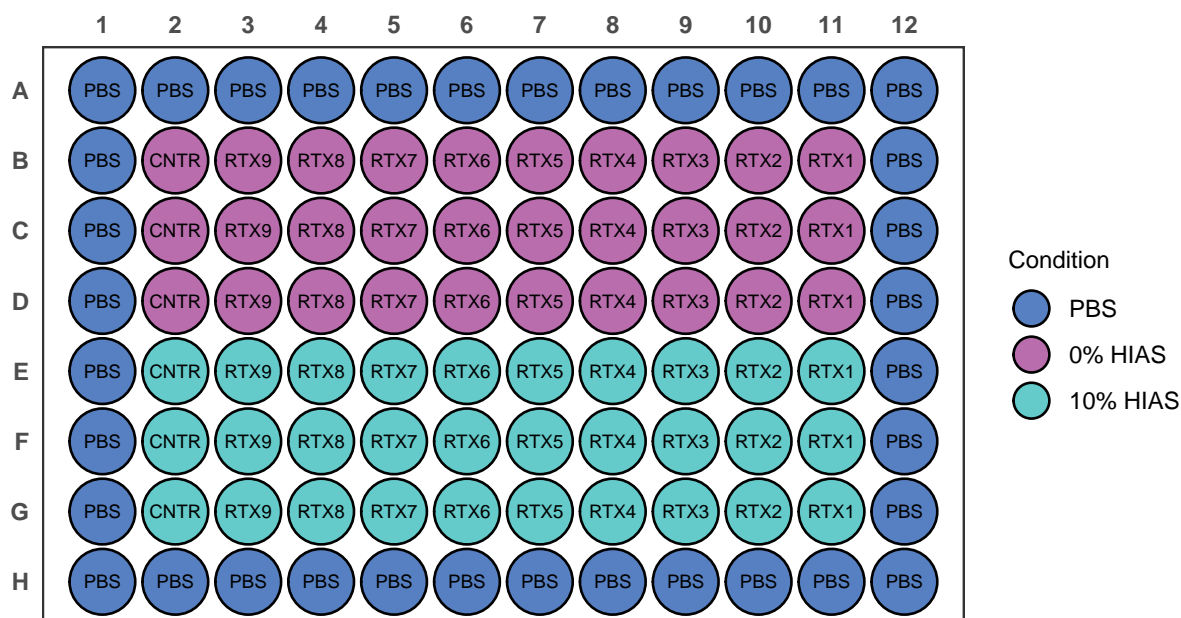

April

Monday 01-04-2024

Cell Culture

N2

- Not ready to split
- Added 5mL media

Ramos

- Count: 1.6×10^6 cells/mL
- Seeded: 5×10^5 cells in 15mL

EC50 Collection: 24329/24331

- Collected plates seeded on 24-03-2024 and 31-3-2024
- EC Plate collection protocol:
 1. Added 40 μ L Cell Titre Blue (CTB) to each conditioned well
 - 20 μ L CTB/100 μ L of conditioned well recommended by manufacturer
 2. Incubated for 1hr at 37C
 3. Read on plate reader according to Cell Titre Blue Protocol
- Results:

EC50 RTX N2-BC4/RAMOS-BC5 24401 - Seeding (DNW)

- Seeded an EC50 experiment comparing the effects of RTX on N2-BC4 and RAMOS-BC5 with or without 10% serum
- Increased [RTX] to 1mg/mL (closer to Sorcha EC50)
- Used RTX provided by Jamie
- Plate seeding protocol:
 1. Diluted cell suspension to seed 20000 cells/well in 25 μ L amounts

Cell Line	Well Number	Cells total	Volume total	CS cells/mL	Flask cell count	Flask Volume
N2-BC4	60 wells	1.2×10^6 cells	1.5 mL	8×10^5	2.92×10^6	410 μ L
RAMOS-BC5	60 wells	1.2×10^6 cells	1.5 mL	8×10^5	2.92×10^6	410 μ L

2. Made RTX dilutions and added to respective wells in 50 μ L
 - [RTX stock] = 10.3 mg/mL - Drug volumes are added in triplicate - Drug volumes are being added constitute 1/4 of well volume: [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 25 μ L per well ~ minimum of 150 μ L per condition needed (recommend 200 μ L)

Dilution ID	[RTX Well]	Media Volume	Source	[RTX WS]
RTX 1	1 mg/mL	400 μ L	15.5 μ L stock	4 mg/mL
RTX 2	500 μ g/mL	200 mL	200 μ L RTX 1	2 μ g/mL
RTX 3	250 μ g/mL	200 mL	200 μ L RTX 2	1 μ g/mL
RTX 4	125 μ g/mL	200 mL	200 μ L RTX 3	5 μ g/mL
RTX 5	62.5 μ g/mL	200 mL	200 μ L RTX 4	2.5 μ g/mL
RTX 6	31.25 μ g/mL	200 mL	200 μ L RTX 5	1.25 μ g/mL
RTX 7	15.6 μ g/mL	200 mL	200 μ L RTX 6	.625 μ g/mL
RTX 8	7.8 μ g/mL	200 mL	200 μ L RTX 7	.3125 μ g/mL
RTX 9	3.9 μ g/mL	200 mL	200 μ L RTX 8	.156 μ g/mL
CNTR	0 μ g/mL	500 mL	-	0 μ g/mL


3. Added 10% serum or media control

- This provides complement factors to bind BL-bound RTX and initiate CDC
- Serum/media amounts added at 50 μ L/well
- Serum volume is added in 1:4 ratio (10 μ L serum in 100 μ L final well volume)

- Serum stock mix is 600 μ L HS : 1.2 mL media / plate

4. Plate is incubated for 72 hrs at 37C

EC50_RTX_N2-BC4_240401

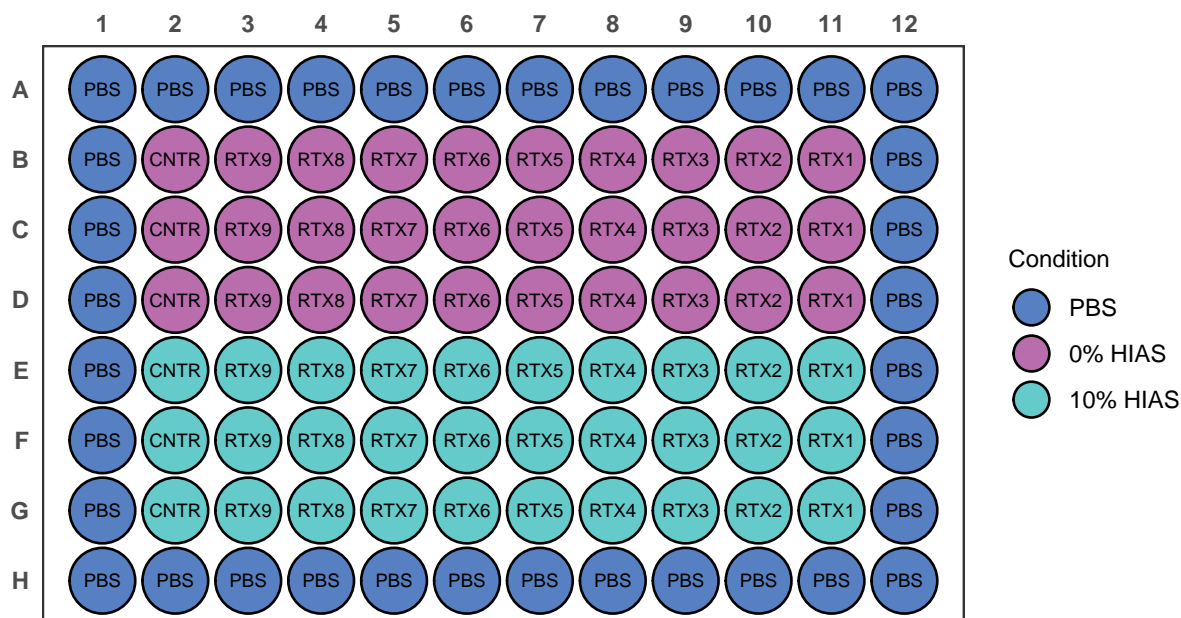


Tuesday 02-04-2024

N2/Ramos CD20 Flow Cytometry

- Got antibody from Chris for CD20
- Cat#: 562873

Cell Stain Protocol:

1. Resuspend cells to 1x10⁶ cells in 1mL
2. Distribute CS into 3 wells of a 96-well plate
3. Spin down @ 1500g for 2min
4. Flick media out
5. Make Cell Stain:
 - 500uL PBS + 2uL Stain -Make this while cells are spinning down
6. Resuspend w/ 100uL in the first well and then resuspend the next 2 wells in the same stain volume
 - Consolidate into 1 well/sample
7. incubate in the fridge for 20min
8. Spin down @ 1500g for 2 min
9. Transfer to FACS tube and run
 - Results:

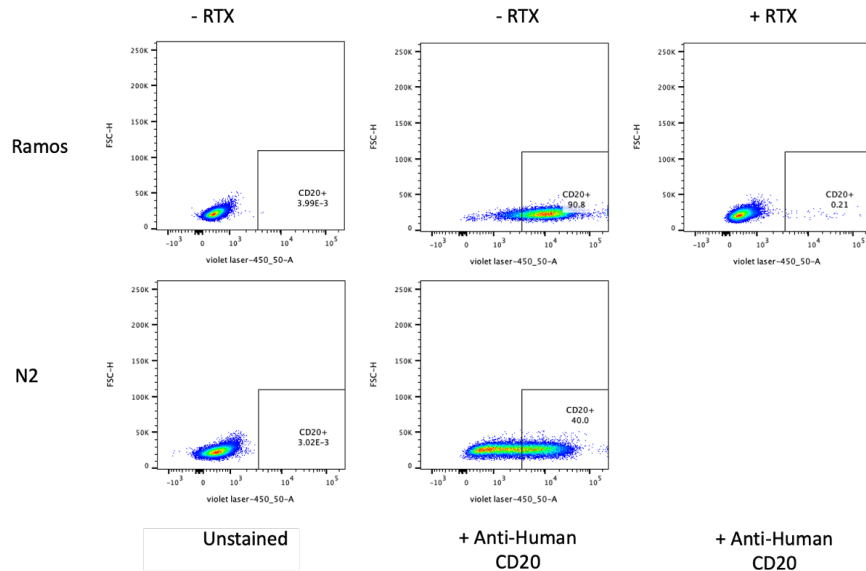

Thursday 04-04-2024

Cell Culture

N2

- Count: 2.1×10^6 cells/mL
- Seeded: 5×10^5 cells in 15mL

Ramos

- Count: 1.9×10^6 cells/mL
- Seeded: 5×10^5 cells in 15mL

EC50 Collection: 240401

- Collected plate seeded on 01-04-2024
- EC Plate collection protocol:
 1. Added $40 \mu\text{L}$ Cell Titre Blue (CTB) to each conditioned well
 - $20 \mu\text{L}$ CTB/ $100 \mu\text{L}$ of conditioned well recommended by manufacturer
 2. Incubated for 1hr at 37°C
 3. Read on plate reader according to Cell Titre Blue Protocol
- Results:DNW

Monday 08-04-2024

Cell Culture

N2

- Count: 1.4×10^6 cells/mL
- Seeded: 3×10^5 cells in 15mL

Ramos

- Count: 2.9×10^6 cells/mL

- Seeded: 3×10^5 cells in 15mL

Tuesday 09-04-2024

Perla Drug Randomizer

- Randomized drugs for Perla Pucci mouse experiment

Drug	Identification
Vehicle	B
Beta 30	F
Beta 10	D
Plo 30	A
Plo 10	E
TESA 4	C
TESA0,4	G

Saturday 14-04-2024

Cell Culture

N2-BC4

- Count: 6.7×10^5 cells/mL
- Seeded: 3×10^5 cells in 15mL

Ramos-BC5

- Count: 6.0×10^5 cells/mL
- Seeded: 3×10^5 cells in 15mL

EC50 RTX N2-BC4/RAMOS-BC5 24414 - Seeding (DNW)

- Seeded an EC50 experiment comparing the effects of RTX on N2-BC4 and RAMOS-BC5 with or without 10% serum
- Increased [RTX] to 1mg/mL (closer to Sorcha EC50)
- Used RTX provided by Jamie
- Looked at cell-intrinsic effects of RTX (no serum)
- Plate seeding protocol:
 1. Diluted cell suspension to seed 20000 cells/well in 25 μ L amounts


Cell Line	Well Number	Cells total	Volume total	CS cells/mL	Flask cell count	Flask Volume
N2-BC4	60 wells	1.2×10^6 cells	1.5 mL	8×10^5	2.92×10^6	410μL
RAMOS-BC5	60 wells	1.2×10^6 cells	1.5 mL	8×10^5	2.92×10^6	410μL

2. Made RTX dilutions and added to respective wells in 50 μ L
 - [RTX stock] = 10.3 mg/mL - Drug volumes are added in triplicate - Drug volumes are being added constitute 1/4 of well volume: [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 25 μ L per well ~ minimum of 150 μ L per condition needed (recommend 200 μ L)

Dilution ID	[RTX Well]	Media Volume	Source	[RTX WS]
RTX 1	1 mg/mL	400 μ L	15.5 μ L stock	4 mg/mL
RTX 2	500 μ g/mL	200 mL	200 μ L RTX 1	2 μ g/mL
RTX 3	250 μ g/mL	200 mL	200 μ L RTX 2	1 μ g/mL
RTX 4	125 μ g/mL	200 mL	200 μ L RTX 3	5 μ g/mL
RTX 5	62.5 μ g/mL	200 mL	200 μ L RTX 4	2.5 μ g/mL
RTX 6	31.25 μ g/mL	200 mL	200 μ L RTX 5	1.25 μ g/mL
RTX 7	15.6 μ g/mL	200 mL	200 μ L RTX 6	.625 μ g/mL
RTX 8	7.8 μ g/mL	200 mL	200 μ L RTX 7	.3125 μ g/mL
RTX 9	3.9 μ g/mL	200 mL	200 μ L RTX 8	.156 μ g/mL
CNTR	0 μ g/mL	500 mL	-	0 μ g/mL

3. Plate is incubated for 72 hrs at 37C

EC50_RTX_N2-BC4_240414

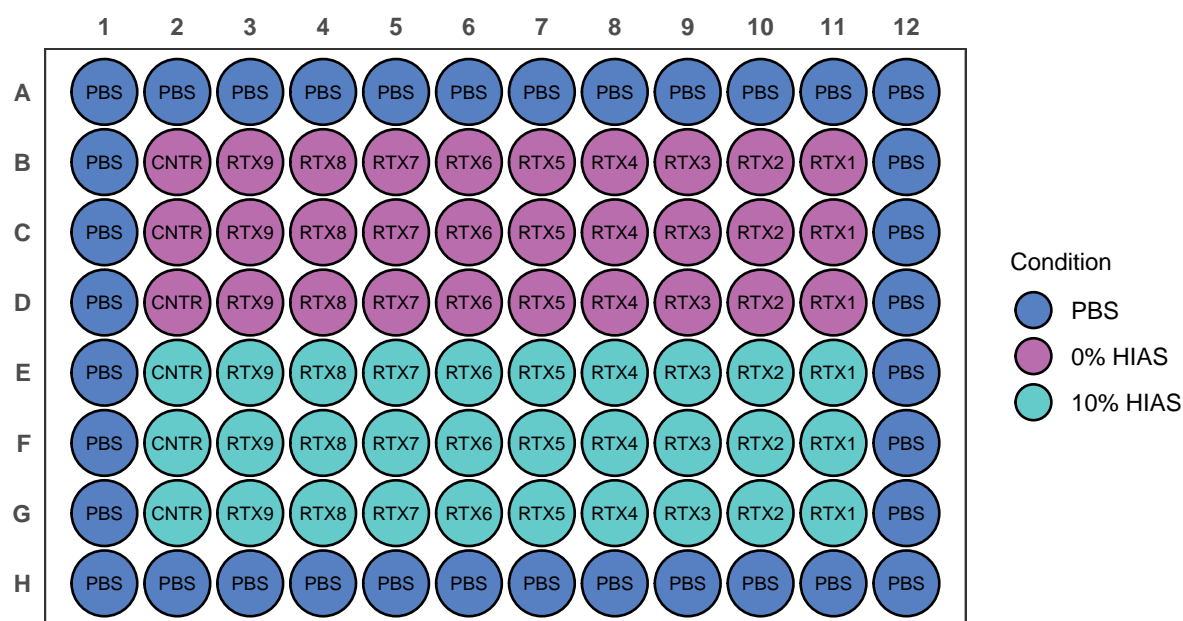

Wednesday 17-04-2024

EC50 Collection: 240414

- Collected plate seeded on 14-04-2024
- EC Plate collection protocol:
 1. Added 40 μ L Cell Titre Blue (CTB) to each conditioned well
 - 20 μ L CTB/100 μ L of conditioned well recommended by manufacturer
 2. Incubated for 1hr at 37C
 3. Read on plate reader according to Cell Titre Blue Protocol
- Results:DNW

Monday 22-04-2024

Cell Culture

N2-BC4

- Seeded: 3×10^5 cells in 15mL

Ramos-BC5

- Seeded: 3×10^5 cells in 15mL

Thawed NK-92

- Obtained from Chris Steele
- Cultured in Advanced RPMI
- Added 200u/mL of IL-2 (#78036.1)
 - Obtained from Emily
 - Product listed as 4.1×10^4 IU/ug
 - Stock diluted to 10ug/mL
 - Added .48 uL stock/mL of media

Wednesday 24-04-2024

Cell Culture

N2-BC4

- Seeded: 3×10^5 cells in 15mL

Ramos-BC5

- Seeded: 3×10^5 cells in 15mL

NK-92

- Looked very unhealthy
- Lots of dead cells/debris in the media
- Spun down, resuspended in 1mL of media
- Split into 4 wells of a 48 well plate with increasing amounts of IL-2

Friday 26-04-2024

Cell Culture

N2-BC4

- Seeded: 3×10^5 cells in 15mL

Ramos-BC5

- Seeded: 3×10^5 cells in 15mL

NK-92

- Cells still look unhealthy
- Tried to consolidate in single well of 24 well plate

- Spun down in eppendorf and resuspended in 1mL of media

Monday 29-04-2024

Cell Culture

N2-BC4

- Seeded: 3×10^5 cells in 15mL

Ramos-BC5

- Seeded: 3×10^5 cells in 15mL

Tuesday 30-04-2024

Cell Culture

N2-BC4

- Count: 3×10^5 cells/mL
- Seeded: 3×10^5 cells in 12mL
- Reseeded FDC cells in T75

Ramos-BC5

- Count: 4×10^5 cells/mL
- Seeded: 4×10^5 cells in 12mL

Human Serum Reciept


- Serum from 3 donors ordered from Cambridge Biosciences
- Serum mixed together, aliquoted, and stored at -80

Product	Product ID	Barcode ID	Volume
Human Serum - Fresh Frozen	SERSKF2SIL10-FSXX	PR24C441891	10 mL
Human Serum - Fresh Frozen	SERSKF2SIL10-FSXX	PR23K435425	10 mL
Human Serum - Fresh Frozen	SERSKF2SIL10-MSXX	PR23D435392	10 mL

Ramos/N2 CDC Testing

- Tested CDC assay adapted from [Ge et al., 2019](#)
 - CDC protocol
1. Seed 250k cells 500uL volumes per well of a 12-well plate
 2. Added RTX or media in 250uL volumes per well
 - [RTX Stock] = 10mg/mL
 - [RTX] working needs to be 4x final desired well volume
 - 1uL RTX stock added to 1mL Media
 3. Add Serum or media in 250uL volumes per well
 - Serum amount should be 4x final % desired

CDC Test N2 + RAMOS 240430

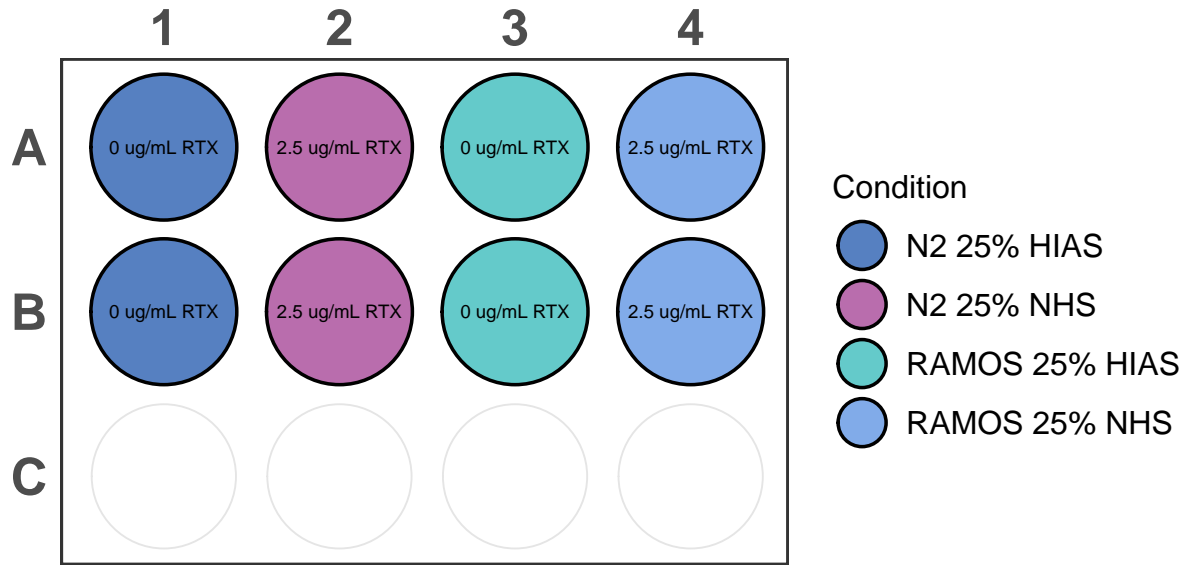

Results: DNW

- No clear differences between +HIAS/+RTX and +NHS/+RTX

- Seemed relatively healthy

- Possible fixes:

* Increase RTX dose * Increase Serum % * Increase the amount of time * Add at the same time instead of pre-incubating with RTX


EC50 RTX RAMOS-BC5 24430 - Seeding

- Seeded an EC50 experiment comparing the effects of RTX on N2-BC4 and RAMOS-BC5 with or without 10% serum
- Increased [RTX] to 1mg/mL (closer to Sorcha EC50)
- Used RTX provided by Jamie
- Plate seeding protocol:
 - Diluted cell suspension to seed 20000 cells/well in 50 μ L amounts
 - Made RTX dilutions and added to respective wells in 25 μ L
 - [RTX stock] = 10.3 mg/mL
 - Drug volumes are added in triplicate
 - Drug volumes are being added constitute 1/4 of well volume: [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 25 μ L per well ~ minimum of 150 μ L per condition needed (recommend 200 μ L)

Dilution ID	Well [RTX] (μ g/mL)	RTX Source	Source Volume (μ L)	Media Volume (μ L)	Working Stock [RTX] (μ L/mL)
RTX 1	257.500000	Stock	60	540	1030.0000
RTX 2	128.750000	RTX 1	300	300	515.0000
RTX 3	64.375000	RTX 2	300	300	257.5000
RTX 4	32.187500	RTX 3	300	300	128.7500
RTX 5	16.093750	RTX 4	300	300	64.3750
RTX 6	8.046875	RTX 5	300	300	32.1875
CNTR	0.000000	-	-	600	0.0000

3. Plate is incubated for 72 hrs at 37C

EC50 CDC Test N2 + RAMOS 240430

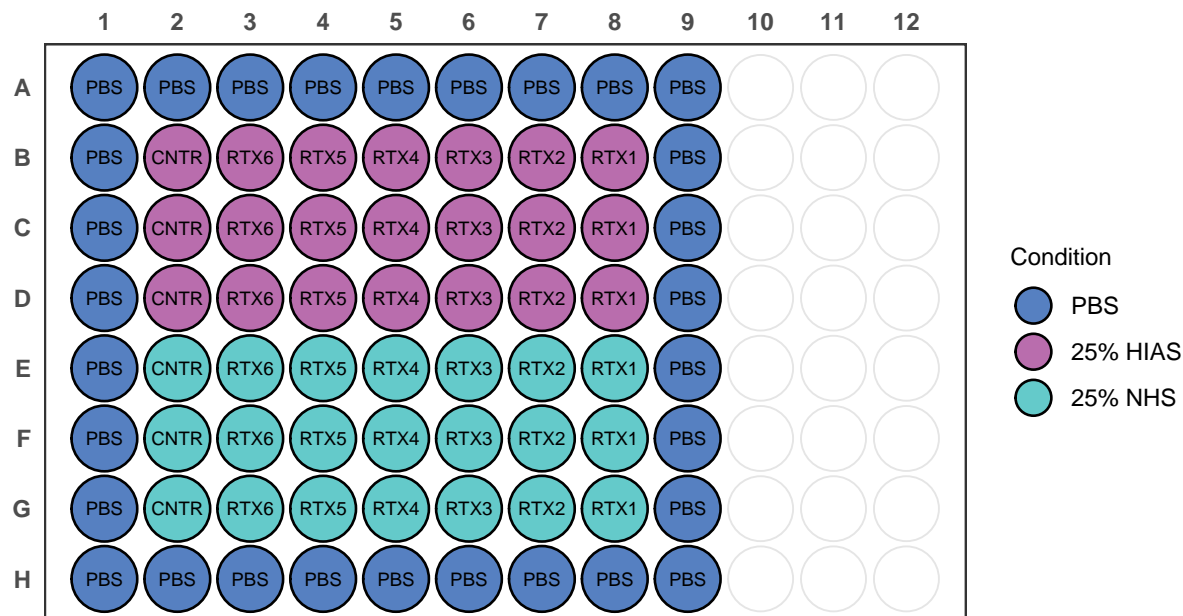

May

Wednesday 01-05-2024

Cell Culture

N2-BC4

- Count: 3×10^5 cells/mL
- Seeded: 3×10^5 cells in 12mL
- Reseeded FDC cells in T75

Ramos-BC5

- Count: 4×10^5 cells/mL
- Seeded: 4×10^5 cells in 12mL

June

Monday 10-06-2024

Cell Culture

- Split cells

Cell Line Name	Count	Seeding Density
Ramos BC 1	3.3×10^6 cells/mL	6.5×10^4 cells/mL
Ramos BC 3	3.3×10^6 cells/mL	6.5×10^4 cells/mL
Ramos BC 5	3.6×10^6 cells/mL	6.5×10^4 cells/mL

EC50 RTX RAMOS-BC5 240610 - Seeding

- Seeded an EC50 experiment comparing the effects of RTX on N2-BC4 and RAMOS-BC5 with or without 10% serum
- Seeded 3 plates each with a different RAMOS barcode population


Plate ID	Cell Line
Plate 1	RAMOS BC 1
Plate 2	RAMOS BC 3
Plate 3	RAMOS BC 5

- Increased [RTX] to >1mg/mL (closer to Sorcha EC50)
- Used RTX provided by Jamie
- Plate seeding protocol:
 1. Diluted cell suspension to seed 10000 cells/well in 50 μ L amounts
 2. Made RTX dilutions and added to respective wells in 25 μ L
 - [RTX stock] = 10.3 mg/mL
 - Drug volumes are added in triplicate
 - Drug volumes are being added constitute 1/4 of well volume: [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 25 μ L per well ~ minimum of 150 μ L per condition needed (recommend 200 μ L)

Dilution ID	[RTX Well]	Media Volume	Source	[RTX WS]
RTX 1	1.66 mg/mL	1200 μ L	800 μ L stock	6.66 mg/mL
RTX 2	500 μ g/mL	1 mL	1 mL RTX 1	3.33 μ g/mL
RTX 3	250 μ g/mL	1 mL	1 mL RTX 2	1.66 μ g/mL
RTX 4	125 μ g/mL	1 mL	1 mL RTX 3	0.833 μ g/mL
RTX 5	62.5 μ g/mL	1 mL	1 mL RTX 4	0.416 μ g/mL
RTX 6	31.25 μ g/mL	1 mL	1 mL RTX 5	0.208 μ g/mL
RTX 7	15.6 μ g/mL	1 mL	1 mL RTX 6	0.104 μ g/mL
RTX 8	7.8 μ g/mL	1 mL	1 mL RTX 7	0.052 μ g/mL
RTX 9	3.9 μ g/mL	1 mL	1 mL RTX 8	0.026 μ g/mL
CNTR	0 μ g/mL	500 mL	-	0 μ g/mL

3. Added HIAS/NHS to indicated wells
 - 25 μ L/well
 - Final well volume = 25% Serum (HIAS/NHS)
4. Plate is incubated for 48 hrs at 37C

EC50 CDC Test N2 + RAMOS 240430

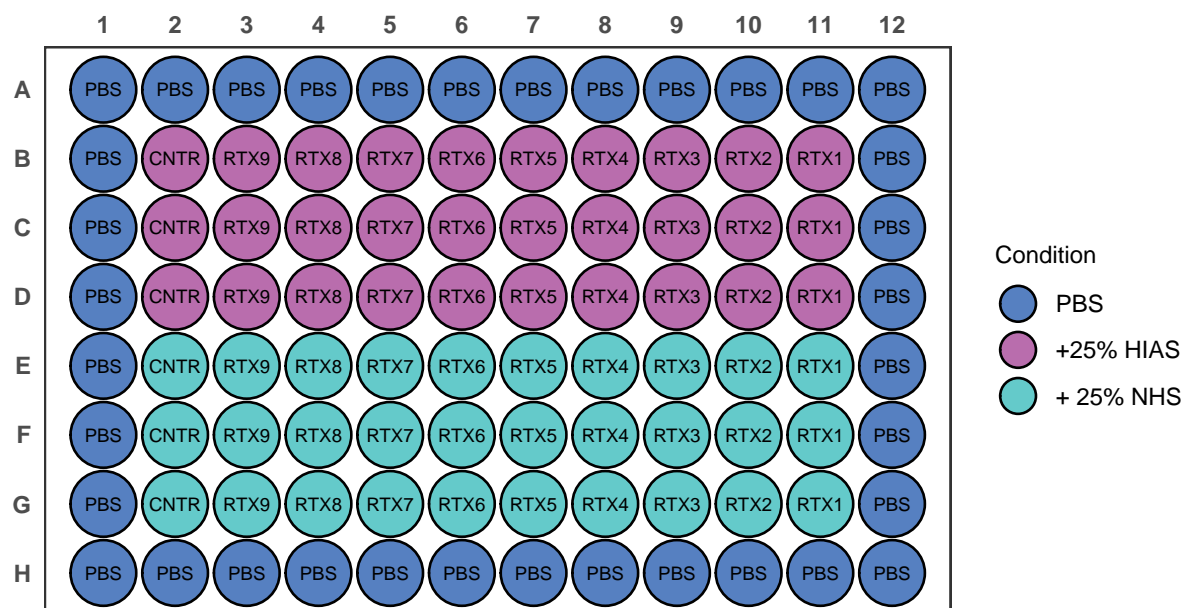

Wednesday 13-06-2024


Cell Culture

- Split cells

Cell Line Name	Count	Seeding Density	Flask	Flask Volume
Ramos BC 1	3.3×10^6 cells/mL	6.5×10^4 cells/mL	T75	20mL

EC50 Collection: EC50_RTX_N2-BC4_10064

- Collected plates seeded on 10-06-2024
- EC Plate collection protocol:
 1. Added $40\mu\text{L}$ Cell Titre Blue (CTB) to each conditioned well
 - $20\mu\text{L}$ CTB/ $100\mu\text{L}$ of conditioned well recommended by manufacturer
 2. Incubated for 1hr at 37°C
 3. Read on plate reader according to Cell Titre Blue Protocol
- Results:

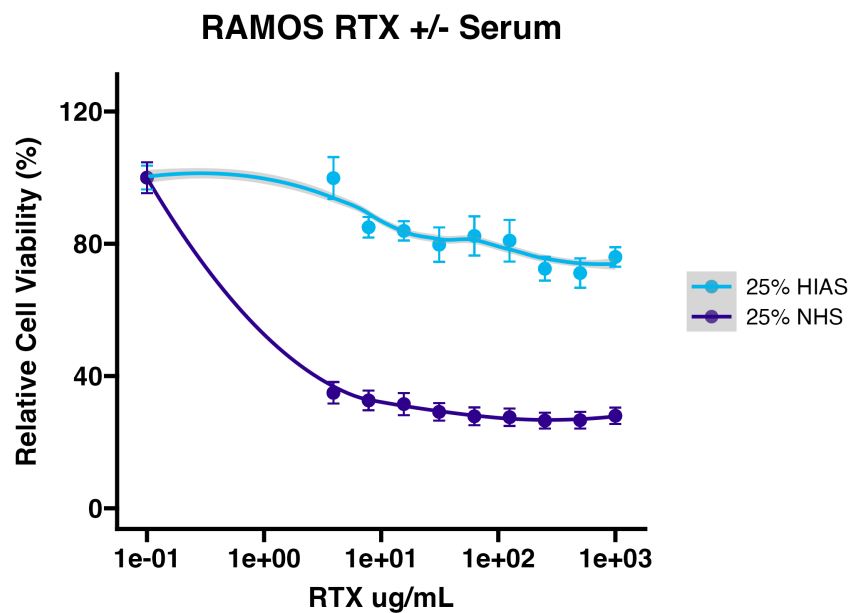

Thursday 13-06-2024

Cryopreservation - RAMOS-BC 1/3/5

- Cryopreserved 3 vials of the following cell lines:
- Ramos-BC 1
- Ramos-BC 3
- Ramos-BC 5
- $\sim 2 \times 10^6$ cells/vial
- Freezing media: FBS + 10% DMSO

Protocol

- 1.

Monday 24-06-2024

Cell Culture

- Transferred RAMOS BC 3 to RPMI-1640

Making RPMI

- RPMI_1

Solution	ID code	Volume	% Total volume
RPMI 1640	—	500 mL	89%
FBS	—	56 mL	10%
Glutamax	—	5.6 mL	1%

Splitting Ramos BC 3

Cell Line Name	Count	Seeding Density	Flask	Flask Volume
Ramos BC 1	2x10 ⁶ cells/mL	1x10 ⁵ cells/mL	T75	20mL

1. Transferred CS to 50mL flask
2. Spun down RAMOS BC 3
3. Removed media and resuspended in 5mL of RPMI_1

Monday 24-06-2024

Cell Culture

- Split RAMOS
- Spun down flask and counted
 - Count: 2x10⁶ cells/mL
 - Added 1mL of CS in 19mL media (1x10⁵cells/mL)

CDC Tests

- Set up 6 well plates and incubated for various amounts of time to look at RTX-CDC in an expanded format
- Sample Plate setup:

EC50 CDC Test RAMOS 240624

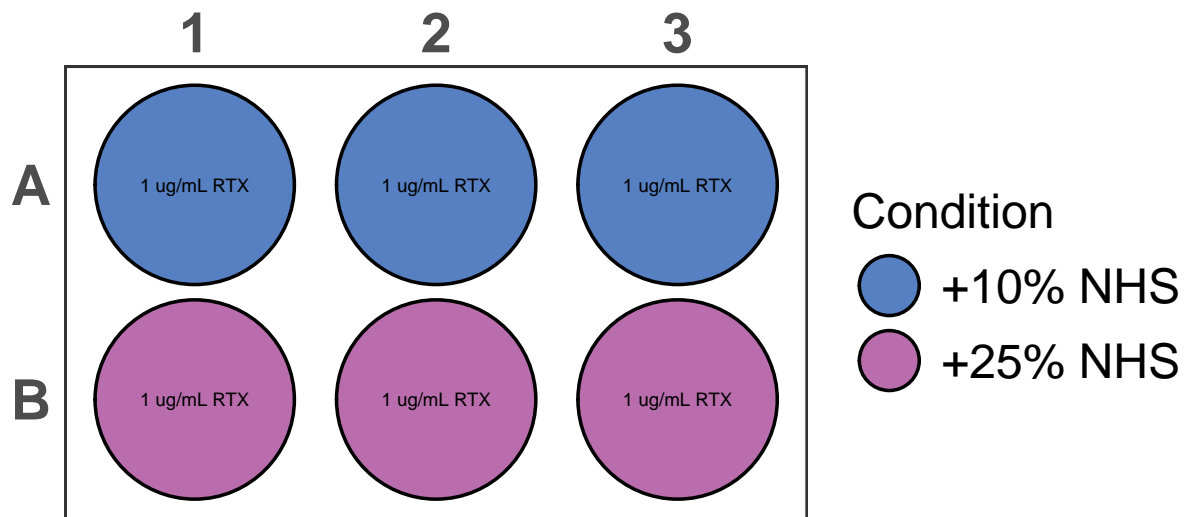


Plate Set Up: 1. Make CS

July

Wednesday 10-07-2024

Cell Culture

- Split RAMOS BC 1
- Count: 1.65x10⁶
- Added 1.5mL CS into 18.5media

RAMOS - RTX In Vitro CDC Pressure

- Began RTX CDC In Vitro experiments
- Seeded RAMOS BC 1 into 2x 6 well plates
- Froze down 3 vials of RAMOS BC 1
 - These can be considered Day 0/Baseline pools
 - Marked with black mark on top of caps

Freezing Down Cells Protocol


1. Counted cells in suspension
2. Took volume of cell suspension such that each vial would contain at least 2×10^6 cells
3. Spun down CS @ 300 ref for 5min
4. Discarded supernatant and resuspended in freezing media such that 1ml of freezing media contains 2×10^6 cells - Freezing media: FBS + 10% DMSO
5. Added 1mL CS in freezing media per cryovial
6. Cryovials were immediately put into freezing caddy and placed in -80 freezer
7. After 24hrs vials removed from freezing caddy and placed in Liquid Nitrogen
 - Location: Tank 1, Rack 5, Box 6

RAMOS RTX DP - Dose 1


- Began RTX CDC In Vitro dosing
- Seeded RAMOS BC 1 into 2x 6 well plates

Plate Layout


RAMOS RTX DP1 Control 240724

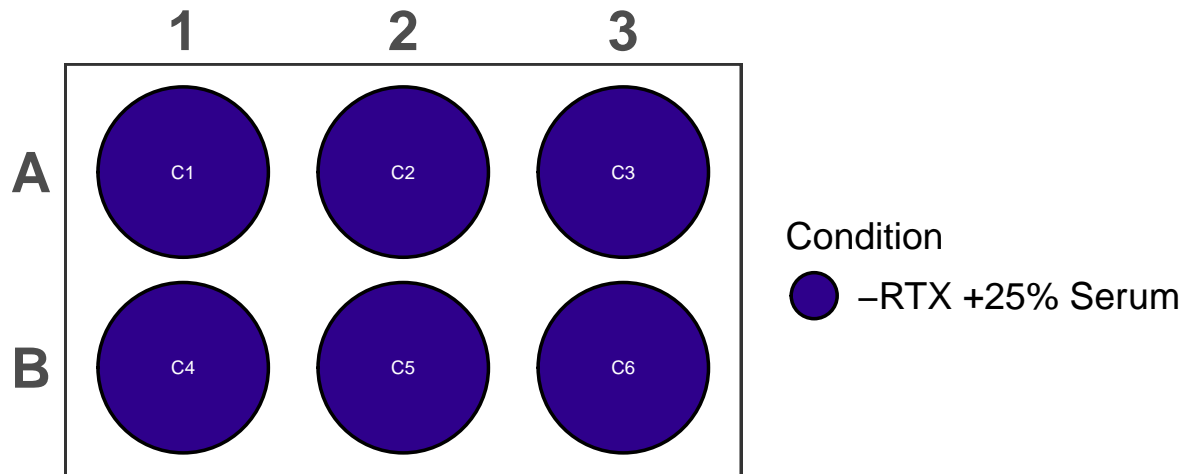

RAMOS RTX DP1 Rituximab 240724

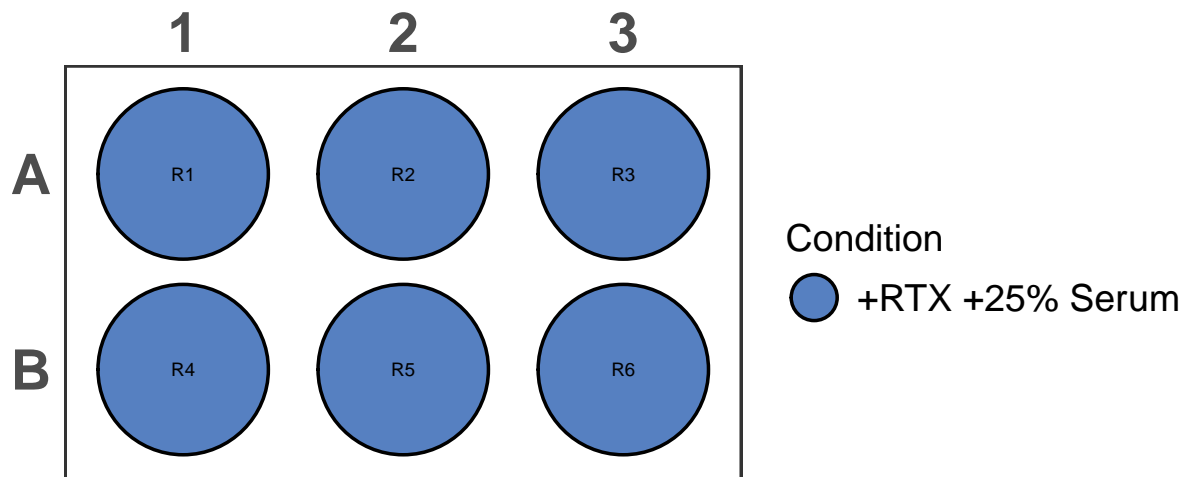

RAMOS BC 1 RTX DP – Dose 1 Plate 1

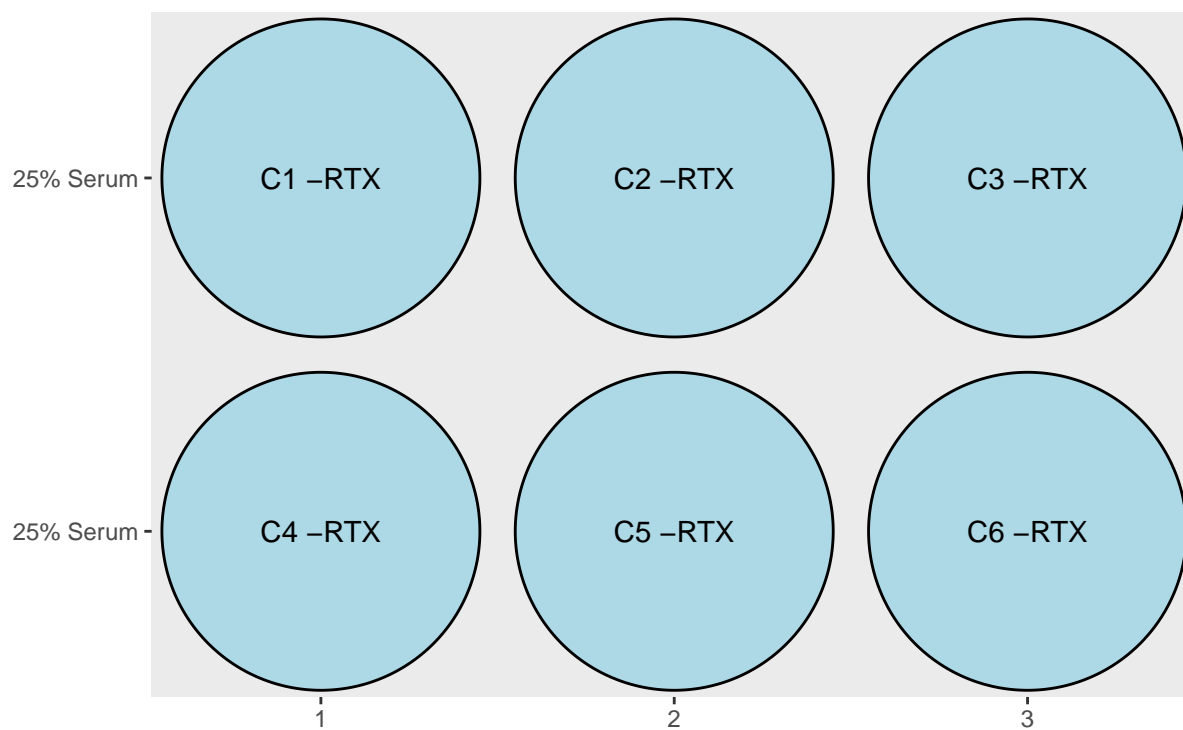

RAMOS BC 1 RTX DP – Dose 1 Plate 2

