

Lab Notebook 2024

Callum Malcolm

Contents

LN Repository	4
Rack 1	4
January 2024	6
February 2024	6
March 2024	6
Monday 24-03-2024	6
MR - CD20 Flow EpiDrug Pretreatment	6
Protocol	6
Results	6
Tuesday 25-03-2024	6
RTX EC50_250324 RBL2 - Seeding	6
Plate seeding protocol:	6
RBL2 RTX DP - Dose 1	7
Dosing Protocol	7
Plate Layout	8
RBL2 RTX DP - Sample Collection Baseline	9
Collection Protocol	9
Thursday 27-03-2024	9
EC50_250325 Collection - RBL2 RTX CDC Baseline	9
Results:	9
RBL2 RTX DP - Dose 1 Collection	9
Collection Protocol	9
Friday 28-03-2024	10
RAMOS-RTXDP SplitSeq: Qubit	10
SS Qubit Results	10
Monday 31-03-2024	10
RAMOS-RTXDP SplitSeq: Library Pooling and Submission	10
RTX EC50_250331 RBL2 - Seeding	10
Plate seeding protocol:	11
April 2025	12
Tuesday 01-04-2024	12
RTX EC50_250401 RBL2 - Seeding	12
Plate seeding protocol:	12

RBL2 RTX DP - Dose 1	13
Dosing Protocol	13
Plate Layout	14
Wednesday 02-04-2024	15
Media Prep	15
Advanced RPMI for PDX	15
EC50_250331 Collection - RBL2 RTX CDC Baseline	15
Results:	15
RTX EC50_250402 RBL2 - Seeding	16
Plate seeding protocol:	16
Thursday 03-04-2024	17
RBL2 RTX DP - Dose 1 Collection	17
Collection Protocol: RTX Treated	17
EC50_250401 Collection - RBL2 RTX CDC Baseline	17
Results:	18
Friday 04-04-2024	18
RBL2 RTX DP - Dose 1 Culture	18
RBL2 RTX-Treated	18
RBL2 Control	18
EC50_250402 Collection - RBL2 RTX CDC Baseline	18
Results:	19
Monday 07-04-2024	19
RBL2 RTX DP - Dose 1 Culture	19
RBL2 Baseline	19
RBL2 RTX-Treated	19
RBL2 Control	19
RTX EC50_250407 RBL2 - Seeding	19
Plate seeding protocol:	20
Tuesday 08-04-2024	21
Cell Culture	21
RBL2 Baseline	21
RBL2 RTX-Treated	21
RBL2 Control	21
Wednesday 09-04-2024	21
EC50_250407 Collection - RBL2 RTX CDC Baseline	21
Results:	22
Friday 11-04-2024	22
Cell Culture	22
RBL2 Baseline + RBL2-Cx1	22
Monday 14-04-2024	22
Cell Culture	22
RBL2 Baseline + RBL2-Cx1	22
Monday 21-04-2024	22
Cell Culture	22
RBL2 Baseline + RBL2-Cx1	22
RTX EC50_250421 RBL2 - Seeding	23

Plate seeding protocol:	23
Tuesday 22-04-2024	24
RTX EC50_250422 RBL2 - Seeding	24
Plate seeding protocol:	24
Wednesday 23-04-2024	25
EC50_250421 Collection - RBL2 RTX CDC Baseline	25
Results:	25
Thursday 23-04-2024	26
Cell Culture	26
RBL2 Barcoded	26
EC50_250422 Collection - RBL2 RTX CDC Baseline	26
Results:	27
Monday 28-04-2024	27
Cell Culture	27
RBL2 Barcoded	27
EC50_250430 Collection - RBL2 RTX CDC Baseline	27
Results:	28

LN Repository

Rack 1

- Location: Tank 2, Rack 1, Row H (Bottom)

Location	Cap ID	Description	Date
1	Grey	Empty - Marker	-
2	Ramos BC 1	Ramos RTX CDC Baseline	12/06/2024
3	Ramos BC 1	Ramos RTX CDC Baseline	12/06/2024
4	Ramos BC 1	Ramos RTX CDC Baseline	12/06/2024
5	C4 DP2	Ramos RTX CDC C4-DP2	-
6	C5 DP2	Ramos RTX CDC C5-DP2	-
7	C1 DP2	Ramos RTX CDC C1-DP2	-
8	R3 DP2	Ramos RTX CDC R3-DP2	-
9	C6 DP2	Ramos RTX CDC C6-DP2	-
10	C3 DP2	Ramos RTX CDC C3-DP2	-
11	C2 DP2	Ramos RTX CDC C2-DP2	-
12	-	-	-
13	-	-	-
14	-	-	-
15	-	-	-
16	-	-	-
17	-	-	-
18	-	-	-
19	-	-	-
20	-	-	-
21	-	-	-
22	-	-	-
23	-	-	-
24	-	-	-
25	-	-	-
26	-	-	-
27	-	-	-
28	-	-	-
29	-	-	-
30	-	-	-
31	Ramos BC 1	Ramos Barcode Pool 1	12/06/2024
32	Ramos BC 1	Ramos Barcode Pool 1	04/01/2023
33	Ramos BC 1	Ramos Barcode Pool 1 *	16/07/2024
34	Ramos BC 3	Ramos Barcode Pool 3	12/06/2024
35	Ramos BC 3	Ramos Barcode Pool 3	04/06/2024
36	Ramos BC 3	Ramos Barcode Pool 3	12/06/2024
37	Ramos BC 5	Ramos BC Pool 5	12/06/2024
38	Ramos BC 5	Ramos Barcode Pool 5	04/16/2024
39	Ramos BC 6	Ramos Barcode Pool 6	04/16/2024
40	-	-	-
41	-	-	-
42	-	-	-
43	-	-	-
44	-	-	-
45	-	-	-
46	-	-	-
47	-	-	-

Location	Cap ID	Description	Date
48	-	-	-
49	-	-	-
50	-	-	-
51	RBL1	RBL1 PDX	31/07/2023
52	RBL1	RBL1 PDX	31/07/2024
53	RBL1 PDX	RBL1 PDX	31/07/2024
54	BLLW	BLLW PDX Pool	31/07/2024
55	BLLW	BLLW PDX Pool	31/07/2024
56	N4	N4 PDX pool	07/11/2023
57	N4	N4 PDX pool	07/11/2023
58	N2 BC	N2 Barcoded pool	11/05/2023
59	N2 BC 5	N2 barcode pool 5	29/04/2024
60	A20	A20 Cell Pool	13/10/2024
61	A20	A20 Stock	13/10/2024
62	-	-	-
63	-	-	-
64	-	-	-
65	-	-	-
66	-	-	-
67	-	-	-
68	-	-	-
69	-	-	-
70	-	-	-
71	-	-	-
72	-	-	-
73	-	-	-
74	-	-	-
75	-	-	-
76	-	-	-
77	-	-	-
78	-	-	-
79	-	-	-
80	-	-	-
81	-	-	-
82	-	-	-
83	-	-	-
84	-	-	-
85	-	-	-
86	-	-	-
87	-	-	-
88	-	-	-
89	-	-	-
90	-	-	-
91	NA	NA	NA
92	NA	NA	NA
93	NA	NA	NA
94	A20 ME	B-IP-724-1L	-
95	A20 ME	B-IP-723 NM	-
96	A20 ME	B-IP-723-2L	-
97	A20 ME	723-2R	-
98	A20 ME	723-1L	-
99	A20 ME	710 NM - A20 Mouse Experiment	-

Location	Cap ID	Description	Date
100	NA	NA	NA

January 2024

February 2024

March 2024

Monday 24-03-2024

MR - CD20 Flow EpiDrug Pretreatment

Protocol

1. 1×10^6 cells from each treatment group were split into 3 wells of a 96-well plate
2. Plate spun at 300xg 5 min and supernatant discarded
3. Cells were resuspended in FACS Staining buffer
 - 1ul CD20 (BD Cat# 562873) per 250uL
 - 12 uL in 3000uL
4. Plate incubated in the dark for 20min at 4C
5. Plate spun at 300xg 5 min and supernatant discarded
6. Samples were resuspended in 400uL PBS and transferred to FACS tubes

Results

- No clear alteration of CD20 expression due to Epigenetic drug pretreatment

Tuesday 25-03-2024

RTX EC50_250324 RBL2 - Seeding

- Seeded an EC50 experiment comparing the effects of RTX on Baseline RBL2 in the presence of 25% NHS
- Seeded 1 plates with the same RBL2 Baseline population
- Used Rixathon (Catalogue#:)

Plate seeding protocol:

1. Diluted cell suspension to seed 10000 cells/well in 50 μ L amounts

Plate	Cell Line	Cell Count	Required Cell total	Required Volume total	CS cells/mL	Stock Volume (uL)	Media Volume (mL)
Plate 1	RBL2	5.65×10^5	7.00×10^5	3.5	1.61×10^5	1238.9	2.2611

2. Made RTX dilutions and added to respective wells in 25 μ L
 - [RTX stock] = 10.3 mg/mL
 - Drug volumes are added in triplicate
 - Drug volumes are being added constitute 1/4 of well volume:

- [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 25 μ L per well ~ minimum of 150 μ L per condition needed (recommend 200 μ L)


Dilution ID	Well [RTX] (μ g/mL)	RTX Source	Source Volume (μ L)	Media Volume (μ L)	Working Stock [RTX] (μ g/mL)
RTX 1	20.0	Stock	4.6500000000000004	595.4	79.8
RTX 2	10.0	RTX 1	300	300.0	39.9
RTX 3	5.0	RTX 2	300	300.0	20.0
RTX 4	2.5	RTX 3	300	300.0	10.0
RTX 5	1.2	RTX 4	300	300.0	5.0
RTX 6	0.6	RTX 5	300	300.0	2.5
RTX 7	0.3	RTX 6	300	300.0	1.2
RTX 8	0.2	RTX 7	300	300.0	0.6
RTX 9	0.1	RTX 8	300	600.0	0.3
CNTR	0.0	-	-	1000.0	0.0

3. Added HIAS/NHS to indicated wells

- 25 μ L/well
- Final well volume = 25% Serum (HIAS/NHS)

4. Plate is incubated for 48 hrs at 37C

EC50_250324 CDC Test RBL2

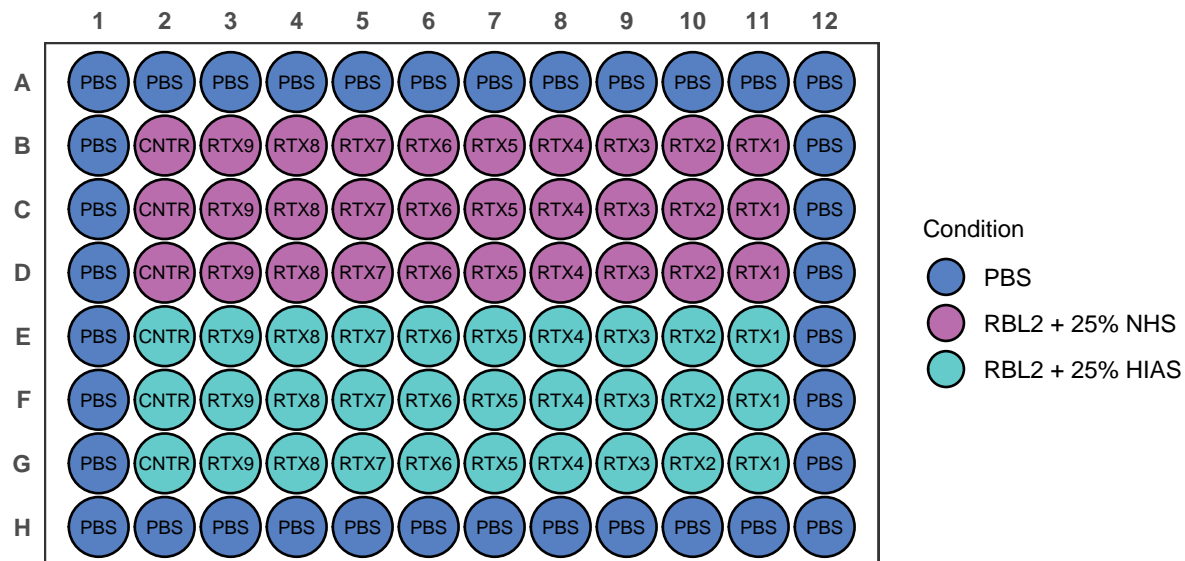

RBL2 RTX DP - Dose 1

- Began RTX CDC In Vitro dosing
- Seeded RBL2 into 2x 6 well plates


Dosing Protocol

1. Count CS and dilute to 2x10⁵ cells in 1 mL
 - If cell count is below either re-culture or add required CS amount, spin down, and resuspend in 1mL
2. Add 1mL of cell suspension containing 2x10⁵ cells to respective wells of 6-well plate


3. Made RTX dilutions and added to respective Rx wells in 500 μ L
 - [RTX stock] = 10.3 mg/mL
 - Drug volumes are being added constitute 1/4 of well volume:
 - [RTX working] needs to be 4x [RTX well]
 - 6 wells per RTX dosing, 500 μ L per well ~ minimum of 3000 μ L per condition needed (recommend 3500 μ L)
 - 1.4 uL RTX stock in 3.5mL media
 - 500 μ L media added to Cx wells
4. Added NHS to all wells
 - 500 μ L/well
 - Final well volume = 25% Serum (NHS)
4. Plates incubated for 24 hrs at 37C

Plate Layout

RBL2 RTX DP1 Control 250325

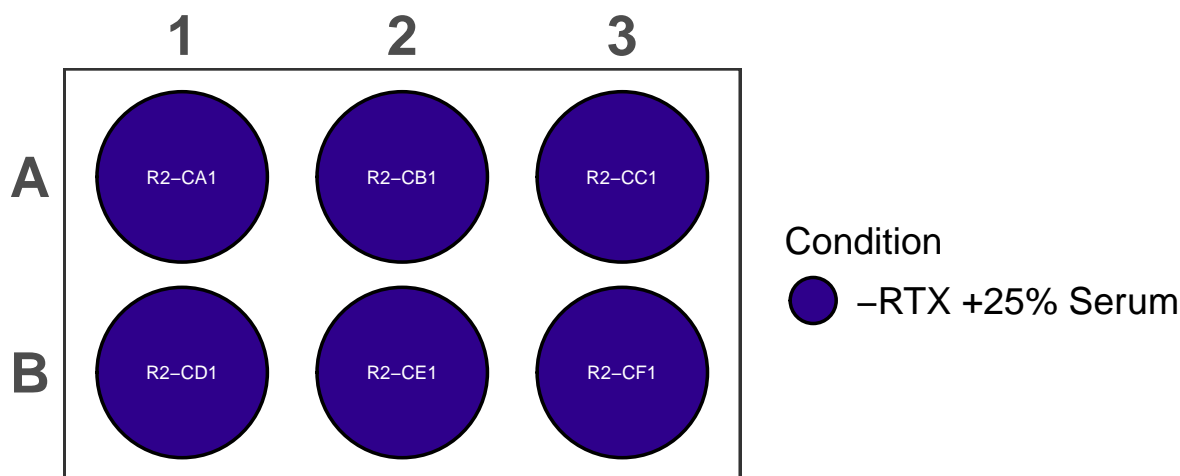

RBL2 RTX DP1 RR-RBL2 250325

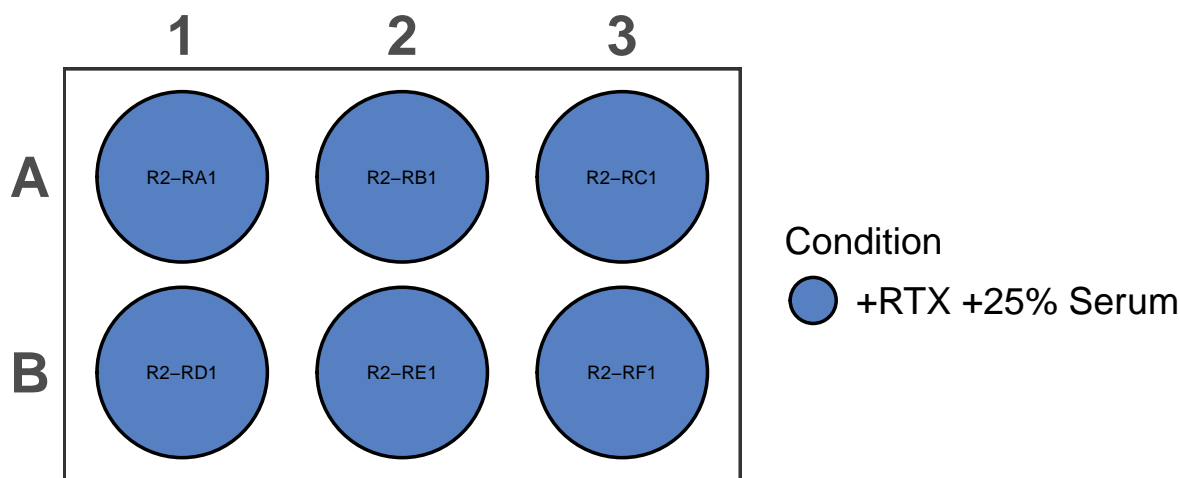

RBL2 RTX DP - Sample Collection Baseline

Collection Protocol

1. 1×10^6 cells added to eppendorf tubes
2. Spun down @ max speed @ 4C
3. Supernatant removed
4. Pellets flash frozen
5. Pellets saved at -80C


Thursday 27-03-2024

EC50_250325 Collection - RBL2 RTX CDC Baseline

- Collected plates seeded on 25-03-2025
- EC Plate collection protocol:
 1. Added 20 μ L Cell Titre Blue (CTB) to each conditioned well
 - 20 μ L CTB/100 μ L of conditioned well recommended by manufacturer
 2. Incubated for 2hr at 37C
 3. Read on plate reader according to Cell Titre Blue Protocol

Results:

- HIAS impact on RBL2 CDC similar to RAMOS experiments
- NHS impact on RTX CDC is consistent with previous experiments
 - Cell viability was only reduced to 50%
 - Will repeat with higher doses of RTX (100ug/ml) and more careful cell counting

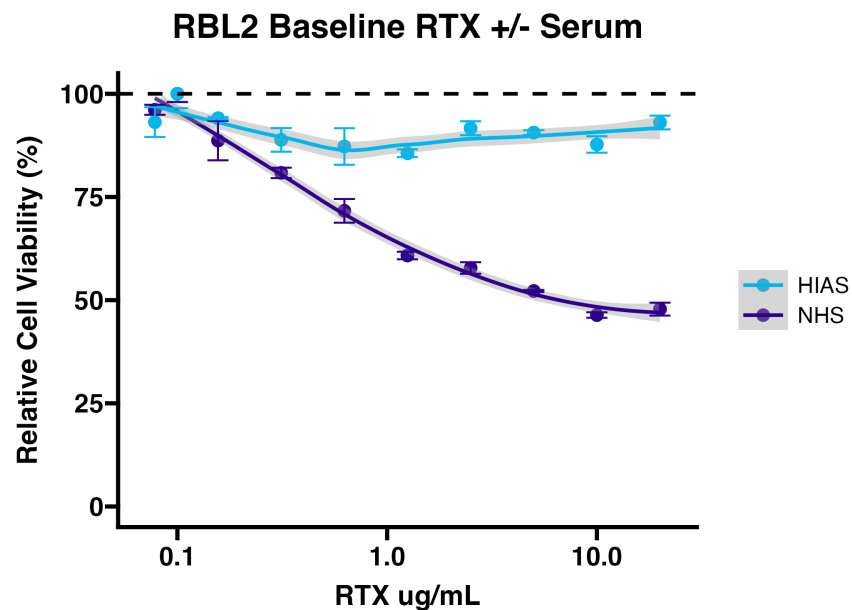

RBL2 RTX DP - Dose 1 Collection

Collection Protocol

1. Well volume transferred to 15ml Eppendorf
2. Eppendorfs supn down at 200 rcf for 6 min
3. Supernatant discarded and resuspended in 2ml

4. Collected cells added to individual wells of 6-well plates

Friday 28-03-2024

RAMOS-RTXDP SplitSeq: Qubit

- Needed to resubmit library due to low concentration

SS Qubit Results

Sample ID	DNA ng/uL	DNA nM
1	20.4	30.9
2	20.8	31.5
3	6.9	10.4
4	24.4	37.0
5	21.0	31.8
6	12.4	18.8
7	19.4	29.4
8	12.3	18.6

- CRUK-CI Genomics core tends to ask for 30-40uL of 10-20nM Library Pool

Monday 31-03-2024

RAMOS-RTXDP SplitSeq: Library Pooling and Submission

- Library Pooled for submission
- Submission ID: SLX-24264

Library Concentration (ng/μl)	Library Concentration (nM)	Library Volume (μl)	10 mM Tris-HCl, pH 8.5 (μl)	Pooling Volume (μl)
Library 1	30.9	10.3	2.2	4.9
Library 2	31.5	9.9	2.6	5.0
Library 3	10.6	0.0	0.0	11.8
Library 4	37.0	8.5	4.0	5.0
Library 5	31.8	9.8	2.7	5.0
Library 6	18.8	0.0	0.0	6.6
Library 7	29.4	10.6	1.9	5.0
Library 8	18.6	0.0	0.0	6.7

Total Pool Volume	Total Pool Concentration (ng/uL)
50	12.5

RTX EC50_250331 RBL2 - Seeding

- Seeded an EC50 experiment comparing the effects of RTX on Baseline RBL2 in the presence of 25% NHS
- Seeded 1 plates with the RBL2 Baseline population
- Used Rixathon (Catalogue#:)

Plate seeding protocol:

1. Diluted cell suspension to seed 10000 cells/well in 50 μ L amounts


Plate	Cell Line	Cell Count	Required Cell total	Required Volume total	CS cells/mL	Stock Volume (uL)	Media Volume (mL)
Plate 1	RBL2	1.26×10^6	7.00×10^5	3.5	3.60×10^5	555.5	2.9445

2. Made RTX dilutions and added to respective wells in 25 μ L
 - [RTX stock] = 10.3 mg/mL
 - Drug volumes are added in triplicate
 - Drug volumes are being added constitute 1/4 of well volume:
 - [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 25 μ L per well ~ minimum of 150 μ L per condition needed (recommend 200 μ L)

Dilution ID	Well [RTX] (μ g/mL)	RTX Source	Source Volume (μ L)	Media Volume (μ L)	Working Stock [RTX] (μ g/mL)
RTX 1	100.0	Stock	23.3	576.6	400.1
RTX 2	50.0	RTX 1	299.95	300.0	200.0
RTX 3	25.0	RTX 2	299.95	300.0	100.0
RTX 4	12.5	RTX 3	299.95	300.0	50.0
RTX 5	6.3	RTX 4	299.95	300.0	25.0
RTX 6	3.1	RTX 5	299.95	300.0	12.5
RTX 7	1.6	RTX 6	299.95	300.0	6.3
RTX 8	0.8	RTX 7	299.95	300.0	3.1
RTX 9	0.4	RTX 8	299.95	599.9	1.6
CNTR	0.0	-	-	1000.0	0.0

3. Added HIAS/NHS to indicated wells
 - 25 μ L/well
 - Final well volume = 25% Serum (HIAS/NHS)
4. Plate is incubated for 48 hrs at 37C

EC50_250331 CDC Test RBL2

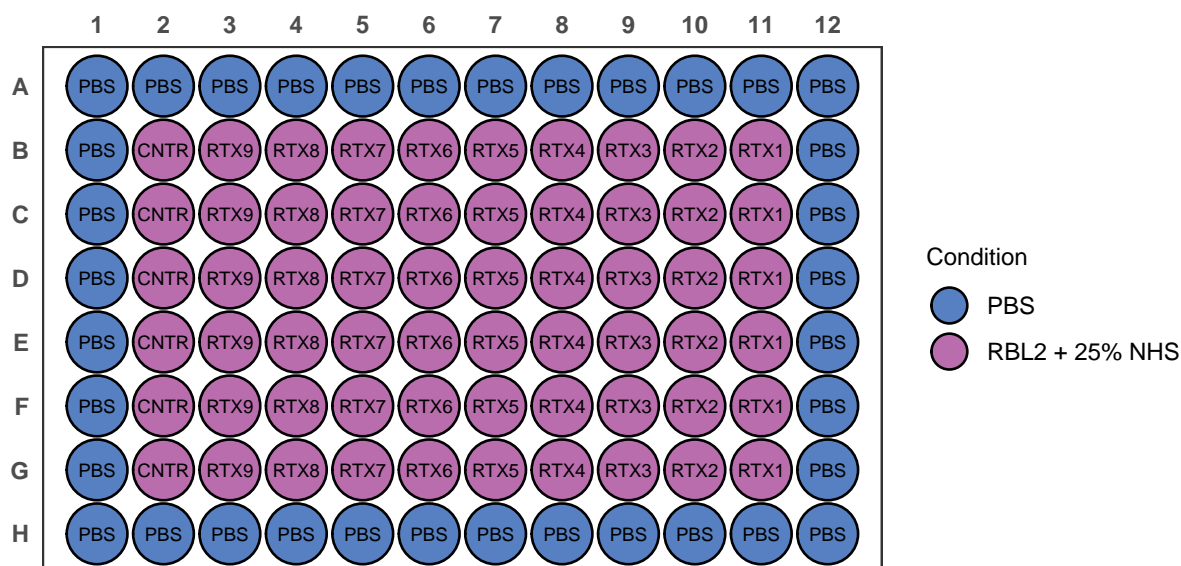

April 2025

Tuesday 01-04-2024

RTX EC50_250401 RBL2 - Seeding

- Seeded an EC50 experiment comparing the effects of RTX on Baseline RBL2 in the presence of 25% NHS
- Seeded 1 plates with the RBL2 Baseline population
- Used Rixathon (Catalogue#:)

Plate seeding protocol:

1. Diluted cell suspension to seed 10000 cells/well in 50 μ L amounts

Plate	Cell Line	Cell Count	Required Cell total	Required Volume total	CS cells/mL	Stock Volume (uL)	Media Volume (mL)
Plate 1	RBL2	9.70×10^5	7.00×10^5	3.5	2.77×10^5	721.6	2.7784

2. Made RTX dilutions and added to respective wells in 25 μ L
 - [RTX stock] = 10.3 mg/mL
 - Drug volumes are added in triplicate
 - Drug volumes are being added constitute 1/4 of well volume:
 - [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 25 μ L per well ~ minimum of 150 μ L per condition needed (recommend 200 μ L)


Dilution ID	Well [RTX] ($\mu\text{g/mL}$)	RTX Source	Source Volume (μL)	Media Volume (μL)	Working Stock [RTX] ($\mu\text{g/mL}$)
RTX 1	100.0	Stock	23.3	576.6	400.1
RTX 2	50.0	RTX 1	299.95	300.0	200.0
RTX 3	25.0	RTX 2	299.95	300.0	100.0
RTX 4	12.5	RTX 3	299.95	300.0	50.0
RTX 5	6.3	RTX 4	299.95	300.0	25.0
RTX 6	3.1	RTX 5	299.95	300.0	12.5
RTX 7	1.6	RTX 6	299.95	300.0	6.3
RTX 8	0.8	RTX 7	299.95	300.0	3.1
RTX 9	0.4	RTX 8	299.95	599.9	1.6
CNTR	0.0	-	-	1000.0	0.0

3. Added HIAS/NHS to indicated wells

- 25 μL /well
- Final well volume = 25% Serum (HIAS/NHS)

4. Plate is incubated for 48 hrs at 37C

EC50_250331 CDC Test RBL2

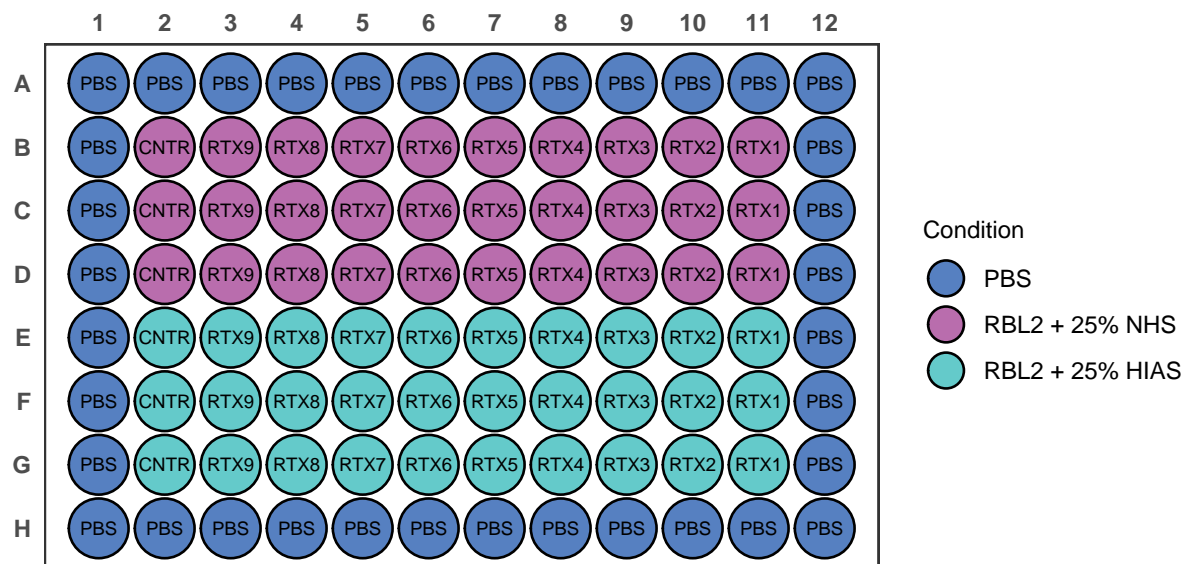

RBL2 RTX DP - Dose 1

- Restarted RTX CDC In Vitro dosing
- Concerns with initial seeding/dosing
 - RBL2 appear to be more resistant to RTX dosing
 - Plan is to start with 25 $\mu\text{g/mL}$ RTX
 - 25 -> 50 -> 100 -> 200
- Seeded RBL2 into 2x 6 well plates


Dosing Protocol

- Count CS and dilute to 5×10^5 cells in 1 mL
 - If cell count is below either re-culture or add required CS amount, spin down, and resuspend in 1mL
- Add 1mL of cell suspension containing 5×10^5 cells to respective wells of 6-well plate


3. Made RTX dilutions and added to respective Rx wells in 500 μ L
 - [RTX stock] = 10.3 mg/mL
 - Drug volumes are being added constitute 1/4 of well volume:
 - [RTX working] needs to be 4x [RTX well]
 - 6 wells per RTX dosing, 500 μ L per well ~ minimum of 3000 μ L per condition needed (recommend 3500 μ L)
 - 25ug/mL [RTX]
 - 34 μ L RTX stock in 3.5mL media
 - 100ug/mL WS
 - 500 μ L media added to Cx wells
4. Added NHS to all wells
 - 500 μ L/well
 - Final well volume = 25% Serum (NHS)
4. Plates incubated for 24 hrs at 37C

Plate Layout

RBL2 RTX DP1 Control 250401

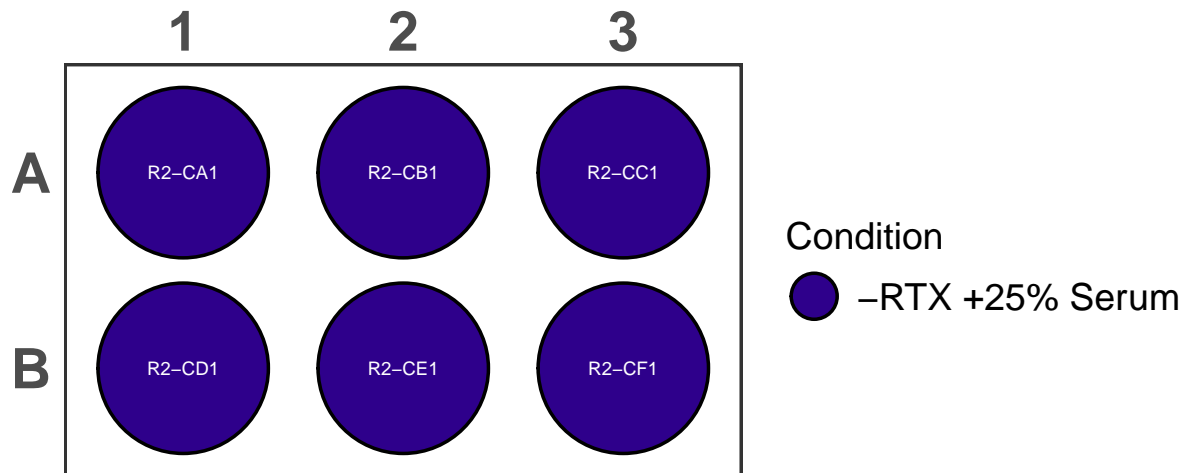

RBL2 RTX DP1 RR-RBL2 250401

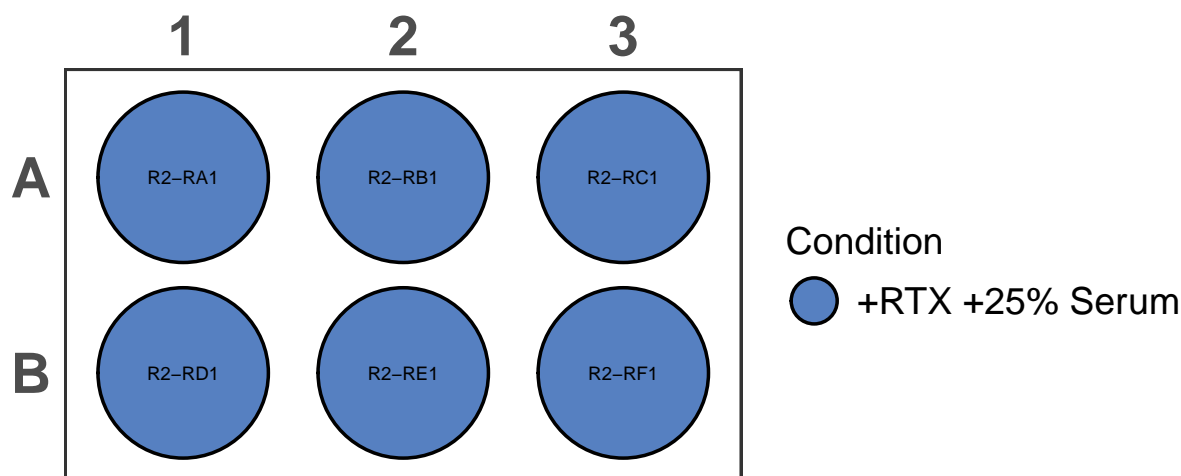

Wednesday 02-04-2024

Media Prep

Advanced RPMI for PDX

- 20% FBS ##### Recipe


Solution	ID code	Volume	% Total volume
Advanced RMPI	—	500 mL	78%
Glutamax	—	6.5 mL	20%
Pen-Strep	—	6.5 mL	1%
FBS	—	128 mL	20%

EC50_250331 Collection - RBL2 RTX CDC Baseline

- Collected plates seeded on 25-03-2025
- EC Plate collection protocol:
 1. Added 20 μ L Cell Titre Blue (CTB) to each conditioned well
 - 20 μ L CTB/100 μ L of conditioned well recommended by manufacturer
 2. Incubated for 2hr at 37C
 3. Read on plate reader according to Cell Titre Blue Protocol

Results:

- Forgot to seed HIAS
- NHS impact on RTX CDC is consistent with previous experiments
 - Cell viability still wasnt reduced below 50%
 - Will repeat with higher doses of RTX (1000ug/ml) and increased cell/well amount (4×10^5)

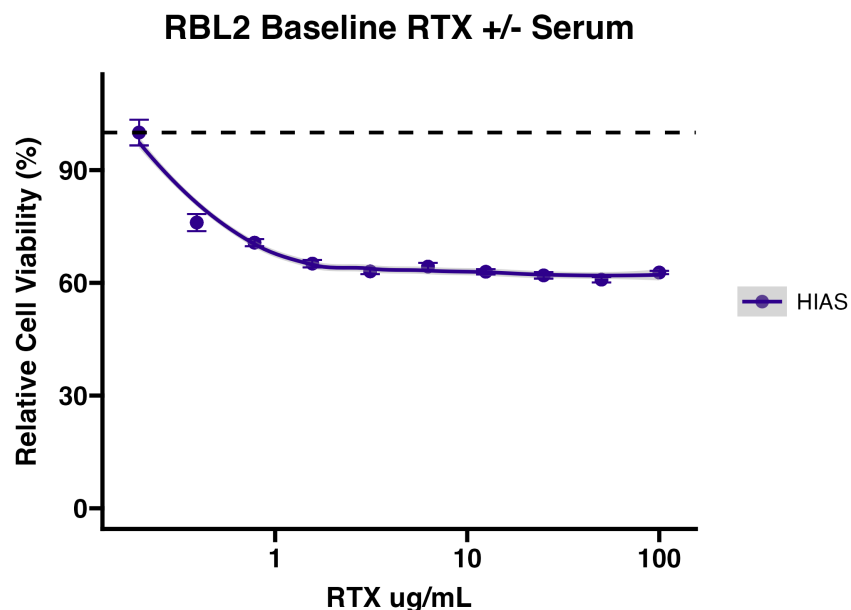

RTX EC50_250402 RBL2 - Seeding

- Seeded an EC50 experiment comparing the effects of RTX on Baseline RBL2 in the presence of 25% NHS
- Seeded 1 plates with the RBL2 Baseline population
- Used Rixathon (Catalogue#:)

Plate seeding protocol:

- Diluted cell suspension to seed 40000 cells/well in 50 μ L amounts

Plate	Cell Line	Cell Count	Required Cell total	Required Volume total	CS cells/mL	Stock Volume (uL)	Media Volume (mL)
Plate 1	RBL2	9.70×10^5	7.00×10^5	3.5	2.77×10^5	721.6	2.7784

- Made RTX dilutions and added to respective wells in 25 μ L
 - [RTX stock] = 10.3 mg/mL
 - Drug volumes are added in triplicate
 - Drug volumes are being added constitute 1/4 of well volume:
 - [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 25 μ L per well ~ minimum of 150 μ L per condition needed (recommend 200 μ L)

Dilution ID	Well [RTX] (μ g/mL)	RTX Source	Source Volume (μ L)	Media Volume (μ L)	Working Stock [RTX] (μ g/mL)
RTX 1	100.0	Stock	23.3	576.6	400.1
RTX 2	50.0	RTX 1	299.95	300.0	200.0
RTX 3	25.0	RTX 2	299.95	300.0	100.0
RTX 4	12.5	RTX 3	299.95	300.0	50.0
RTX 5	6.3	RTX 4	299.95	300.0	25.0


Dilution ID	Well [RTX] ($\mu\text{g/mL}$)	RTX Source	Source Volume (μL)	Media Volume (μL)	Working Stock [RTX] ($\mu\text{g/mL}$)
RTX 6	3.1	RTX 5	299.95	300.0	12.5
RTX 7	1.6	RTX 6	299.95	300.0	6.3
RTX 8	0.8	RTX 7	299.95	300.0	3.1
RTX 9	0.4	RTX 8	299.95	599.9	1.6
CNTR	0.0	-	-	1000.0	0.0

3. Added HIAS/NHS to indicated wells

- 25 μL /well
- Final well volume = 25% Serum (HIAS/NHS)

4. Plate is incubated for 48 hrs at 37C

EC50_250331 CDC Test RBL2

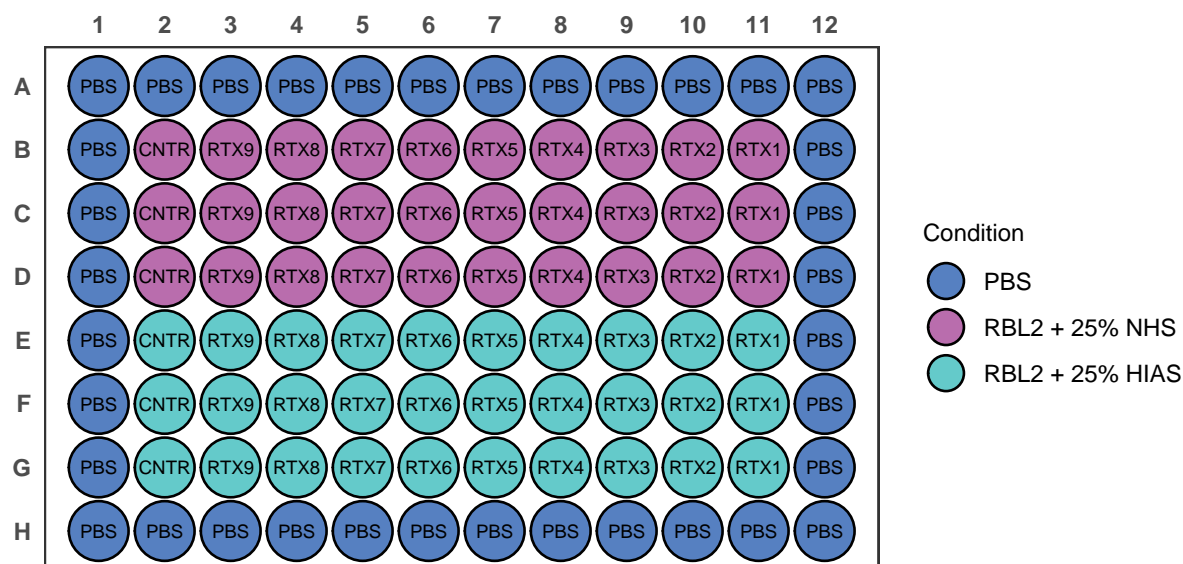

Thursday 03-04-2024

RBL2 RTX DP - Dose 1 Collection

Collection Protocol: RTX Treated

1. Well volume transferred to 15ml Eppendorf
 2. Eppendorfs supn down at 200 rcf for 6 min
 3. Supernatant discarded and resuspended in 1ml
 4. Collected RTX-treated cells added to individual wells of 24-well plates
- Control wells expanded to T25s in 4ml media


EC50_250401 Collection - RBL2 RTX CDC Baseline

- Collected plates seeded on 25-03-2025
- EC Plate collection protocol:
 1. Added 20 μL Cell Titre Blue (CTB) to each conditioned well
 - 20 μL CTB/100 μL of conditioned well recommended by manufacturer
 2. Incubated for 2hr at 37C

3. Read on plate reader according to Cell Titre Blue Protocol

Results:

- HIAS Results as expected
 - Need to seed one more as reference
- NHS impact on RTX CDC is consistent with previous experiments
 - Cell viability still wasn't reduced below 50%
 - Will repeat with higher doses of RTX (1000ug/ml) and increased cell/well amount (4×10^5)

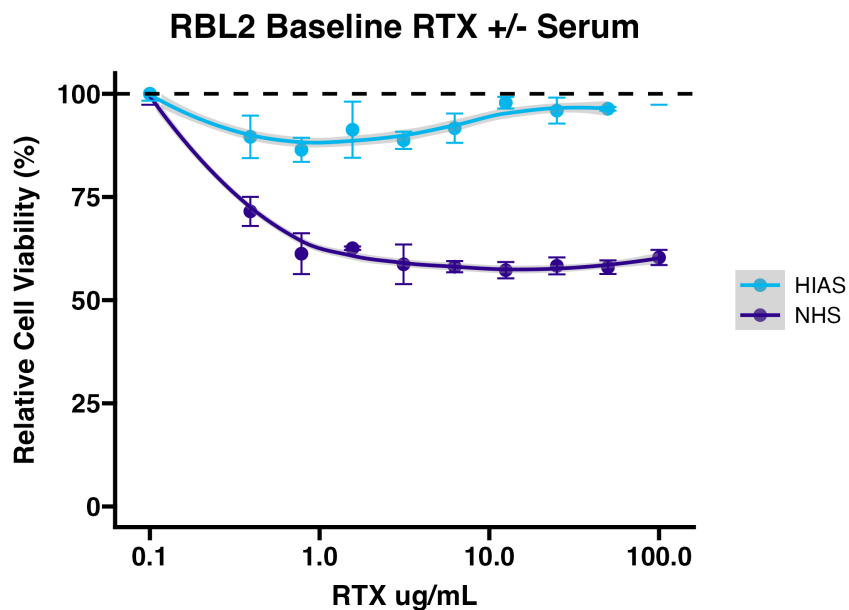

Friday 04-04-2024

RBL2 RTX DP - Dose 1 Culture

RBL2 RTX-Treated

- Expanded to T25
- Transferred 1mL of CS and added 3mL of fresh media
- Cells looked relatively healthy, good growth from yesterday

RBL2 Control


- Spun down and resuspended in 4 ml of media
- Look fine

EC50_250402 Collection - RBL2 RTX CDC Baseline

- Collected plates seeded on 02-04-2025
- EC Plate collection protocol:
 1. Added 20μL Cell Titre Blue (CTB) to each conditioned well
 - 20μL CTB/100μL of conditioned well recommended by manufacturer
 2. Incubated for 2hr at 37C
 3. Read on plate reader according to Cell Titre Blue Protocol

Results:

- Seeding conditions here are wildly different
 - Down to 30% viability at ~4ug/mL
 - Seems like the issue was cell seeding
 - Will repeat 07-04-2025 with 20ug/mL but 4×10^5 cells

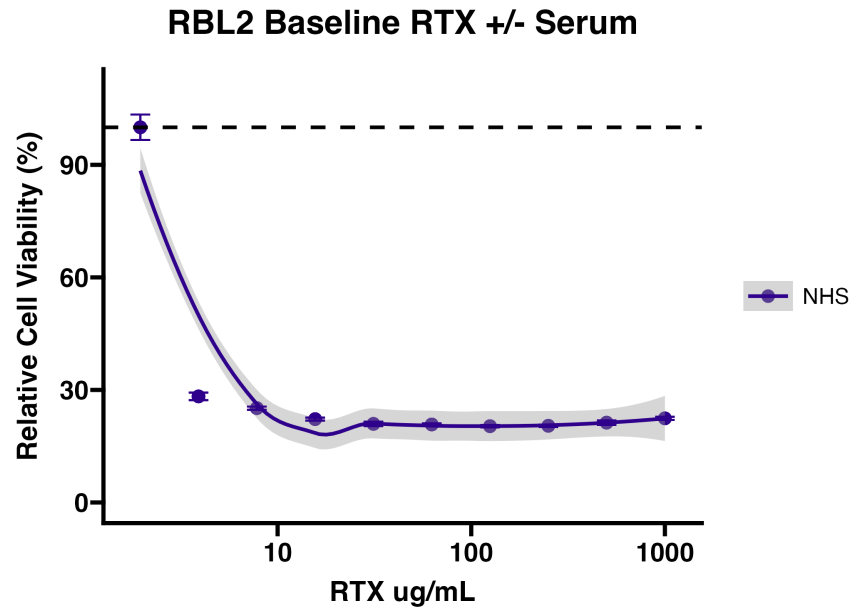

Monday 07-04-2024

RBL2 RTX DP - Dose 1 Culture

RBL2 Baseline

- Healthy
- Use for EC50
- Split
- Would like to discontinue the line once we have firmly established optimal dose-response conditions for RTX

RBL2 RTX-Treated

- Cells look rough
- Leave for one more day to see if they're just growing slowly

RBL2 Control

- Look good
- Expand to T75

RTX EC50_250407 RBL2 - Seeding

- Seeded an EC50 experiment comparing the effects of RTX on Baseline RBL2 in the presence of 25% NHS or 25% HIAS
 - Looking to establish optimal EC50 range
- Seeded 1 plates with the RBL2 Baseline population based on previous results:

- 40k cells per well
- Starting at 100ug/mL dose
- Used Rixathon (Catalogue#:)

Plate seeding protocol:

1. Diluted cell suspension to seed 40000 cells/well in 50 μ L amounts


Plate	Cell Line	Cell Count	Required Cell total	Required Volume total	CS cells/mL	Stock Volume (uL)	Media Volume (mL)
Plate 1	RBL2	1.63×10^6	2.80×10^6	4	4.08×10^5	1717.7	2.2823

2. Made RTX dilutions and added to respective wells in 25 μ L
 - [RTX stock] = 10.3 mg/mL
 - Drug volumes are added in triplicate
 - Drug volumes are being added constitute 1/4 of well volume:
 - [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 25 μ L per well ~ minimum of 150 μ L per condition needed (recommend 200 μ L)

Dilution ID	Well [RTX] (μ g/mL)	RTX Source	Source Volume (μ L)	Media Volume (μ L)	Working Stock [RTX] (μ g/mL)
RTX 1	100.0	Stock	19.417475728155338	480.6	400.0
RTX 2	50.0	RTX 1	250	250.0	200.0
RTX 3	25.0	RTX 2	250	250.0	100.0
RTX 4	12.5	RTX 3	250	250.0	50.0
RTX 5	6.2	RTX 4	250	250.0	25.0
RTX 6	3.1	RTX 5	250	250.0	12.5
RTX 7	1.6	RTX 6	250	250.0	6.2
RTX 8	0.8	RTX 7	250	250.0	3.1
RTX 9	0.4	RTX 8	250	500.0	1.6
CNTR	0.0	-	-	1000.0	0.0

3. Added HIAS/NHS to indicated wells
 - 25 μ L/well
 - Final well volume = 25% Serum (HIAS/NHS)
4. Plate is incubated for 48 hrs at 37C

EC50_250407 CDC Test RBL2

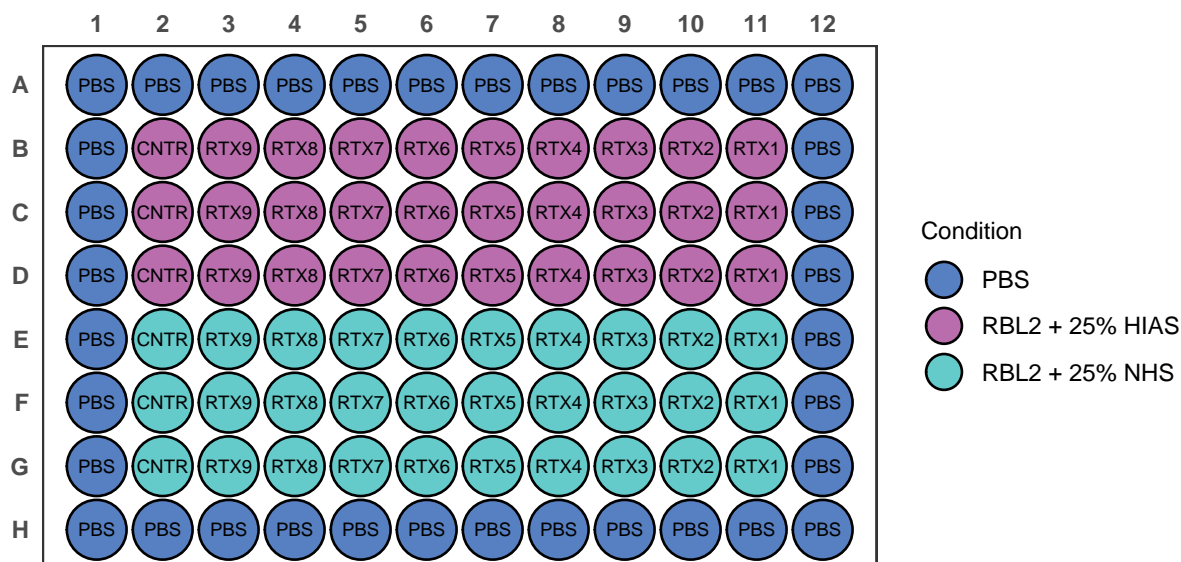

Tuesday 08-04-2024

Cell Culture

RBL2 Baseline

- Healthy
- Use for EC50
 - Seed tomorrow based on results
- Split 1/4
- Would like to discontinue the line once we have firmly established optimal dose-response conditions for RTX

RBL2 RTX-Treated

- Cells look rough
- Spun down and resuspend in 12 well plates
- Expanded too quickly

RBL2 Control

- Look good
- Leave for another day


Wednesday 09-04-2024

EC50_250407 Collection - RBL2 RTX CDC Baseline

- Collected plates seeded on 07-04-2025
- EC Plate collection protocol:
 1. Added 20 μ L Cell Titre Blue (CTB) to each conditioned well
 - 20 μ L CTB/100 μ L of conditioned well recommended by manufacturer
 2. Incubated for 2hr at 37C
 3. Read on plate reader according to Cell Titre Blue Protocol

Results:

- Seeding conditions seem to be optimized
 - HIAS conditions seem weird
 1. Seems to be due to seeding error in the first row
 - RTX+/NHS+ pushed down to acceptable levels at highest dose
 1. Could potentially go higher? (100ug/mL??)
 2. Would prefer to limit to 20ug/mL for resource maintenance

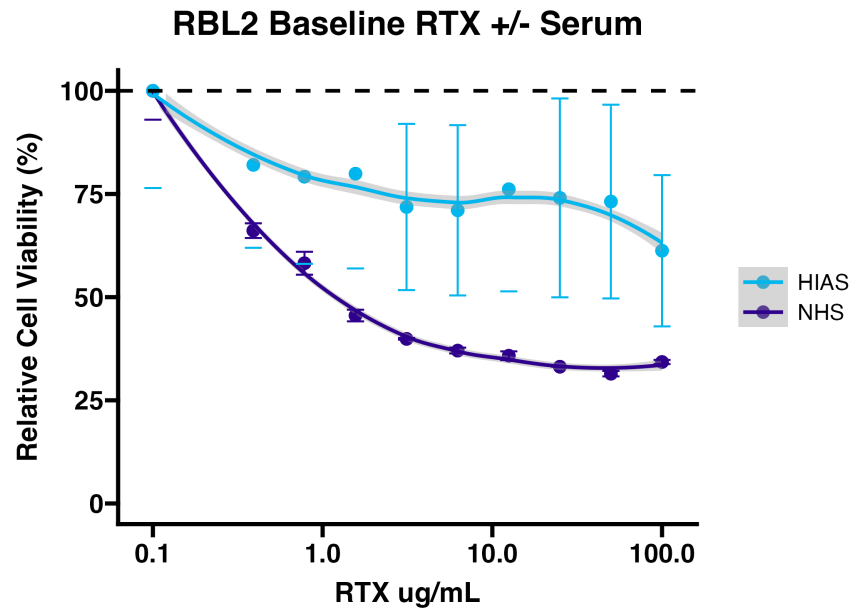

Friday 11-04-2024

Cell Culture

RBL2 Baseline + RBL2-Cx1

- Split 1/10 ### RBL2-Rx1
- Add 1ml media to 6-well

Monday 14-04-2024

Cell Culture

RBL2 Baseline + RBL2-Cx1

- Split 1/10 ### RBL2-Rx1
- Expanded to T25

Monday 21-04-2024

Cell Culture

RBL2 Baseline + RBL2-Cx1

- Split 1/10 ### RBL2-Rx1
- Expanded to T75

RTX EC50_250421 RBL2 - Seeding

- Seeded an EC50 experiment comparing the effects of RTX on Baseline RBL2 in the presence of 25% NHS or 25% HIAS
 - Looking to establish optimal EC50 range
- Seeded 1 plates with the RBL2 Baseline population based on previous results:
 - 40k cells per well
 - Starting at 100ug/mL dose
- Used Rixathon (Catalogue#:)

Plate seeding protocol:

1. Diluted cell suspension to seed 40000 cells/well in 50 μ L amounts
 - Seeding was incorrect


Plate	Cell Line	Cell Count	Required Cell total	Required Volume total	CS cells/mL	Stock Volume (uL)	Media Volume (mL)
Plate 1	RBL2	2.80×10^6	3.20×10^6	4	8.00×10^5	1142.8	2.8572

2. Made RTX dilutions and added to respective wells in 25 μ L
 - [RTX stock] = 10.3 mg/mL
 - Drug volumes are added in triplicate
 - Drug volumes are being added constitute 1/4 of well volume:
 - [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 25 μ L per well ~ minimum of 150 μ L per condition needed (recommend 200 μ L) ** Incorrect dosing

Dilution ID	Well [RTX] (μ g/mL)	RTX Source	Source Volume (μ L)	Media Volume (μ L)	Working Stock [RTX] (μ g/mL)
RTX 1	100.0	Stock	19.417475728155338	480.6	400.0
RTX 2	50.0	RTX 1	250	250.0	200.0
RTX 3	25.0	RTX 2	250	250.0	100.0
RTX 4	12.5	RTX 3	250	250.0	50.0
RTX 5	6.2	RTX 4	250	250.0	25.0
RTX 6	3.1	RTX 5	250	250.0	12.5
RTX 7	1.6	RTX 6	250	250.0	6.2
RTX 8	0.8	RTX 7	250	250.0	3.1
RTX 9	0.4	RTX 8	250	500.0	1.6
CNTR	0.0	-	-	1000.0	0.0

3. Added HIAS/NHS to indicated wells
 - 25 μ L/well
 - Final well volume = 25% Serum (HIAS/NHS)
4. Plate is incubated for 48 hrs at 37C

EC50_250421 CDC Test RBL2

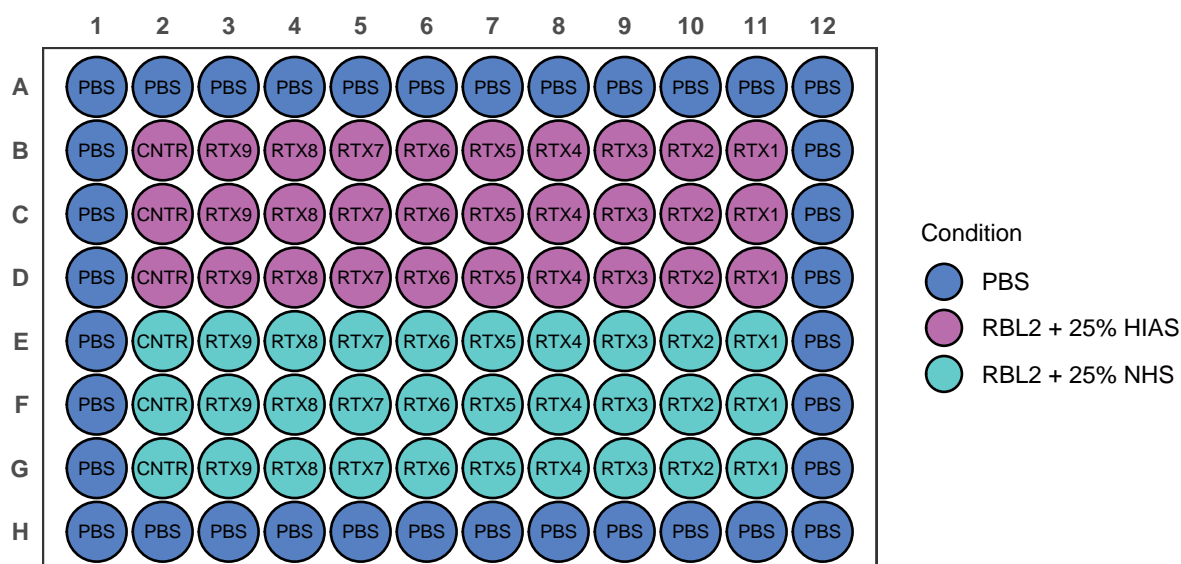

Tuesday 22-04-2024

RTX EC50_250422 RBL2 - Seeding

- Seeded an EC50 experiment comparing the effects of RTX on Baseline RBL2 in the presence of 25% NHS or 25% HIAS
 - Looking to establish optimal EC50 range
- Seeded 1 plates with the RBL2 Baseline population based on previous results:
 - 40k cells per well
 - Starting at 100ug/mL dose
- Used Rixathon (Catalogue#:)

Plate seeding protocol:

- Diluted cell suspension to seed 40000 cells/well in 50μL amounts
 - Seeding was incorrect

Plate	Cell Line	Cell Count	Required Cell total	Required Volume total	CS cells/mL	Stock Volume (uL)	Media Volume (mL)
Plate 1	RBL2	8.00×10^5	3.20×10^6	4	8.00×10^5	4000	0

- Made RTX dilutions and added to respective wells in 25μL
 - [RTX stock] = 10.3 mg/mL
 - Drug volumes are added in triplicate
 - Drug volumes are being added constitute 1/4 of well volume:
 - [RTX working] needs to be 4x [RTX well]
 - 6 wells per condition, 25μL per well ~ minimum of 150μL per condition needed (recommend 200μL) ** Incorrect dosing


Dilution ID	Well [RTX] ($\mu\text{g/mL}$)	RTX Source	Source Volume (μL)	Media Volume (μL)	Working Stock [RTX] ($\mu\text{g/mL}$)
RTX 1	100.0	Stock	19.417475728155338	480.6	400.0
RTX 2	50.0	RTX 1	250	250.0	200.0
RTX 3	25.0	RTX 2	250	250.0	100.0
RTX 4	12.5	RTX 3	250	250.0	50.0
RTX 5	6.2	RTX 4	250	250.0	25.0
RTX 6	3.1	RTX 5	250	250.0	12.5
RTX 7	1.6	RTX 6	250	250.0	6.2
RTX 8	0.8	RTX 7	250	250.0	3.1
RTX 9	0.4	RTX 8	250	500.0	1.6
CNTR	0.0	-	-	1000.0	0.0

3. Added HIAS/NHS to indicated wells

- 25 μL /well
- Final well volume = 25% Serum (HIAS/NHS)

4. Plate is incubated for 48 hrs at 37C

EC50_250422 CDC Test RBL2

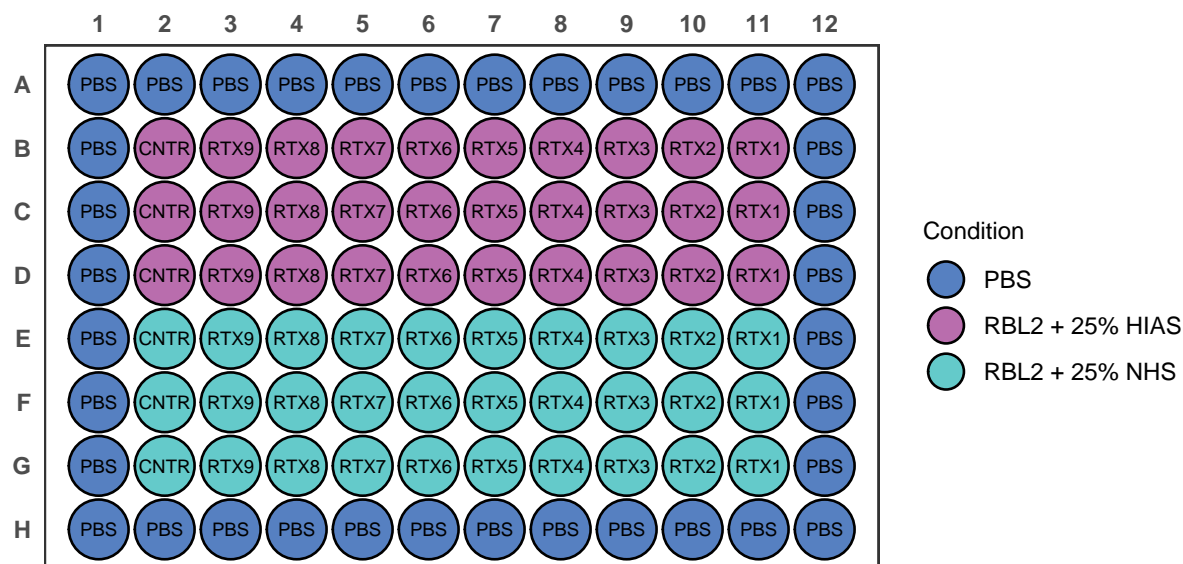

Wednesday 23-04-2024


EC50_250421 Collection - RBL2 RTX CDC Baseline

- Collected plates seeded on 21-04-2025
- EC Plate collection protocol:
 1. Added 20 μL Cell Titre Blue (CTB) to each conditioned well
 - 20 μL CTB/100 μL of conditioned well recommended by manufacturer
 2. Incubated for 2hr at 37C
 3. Read on plate reader according to Cell Titre Blue Protocol

Results:

- Seeding error

- Seeded incorrect number of cells (undercounted - overseeded)
- HIAS condition is lower than ideal
 1. Seeding error?
 2. HIAS quality?
- Still struggling to crack 50% Cell viability reduction
 1. In this case likely due to seeding
 2. EC50_250422 looks similar however
 3. Seed less cells?

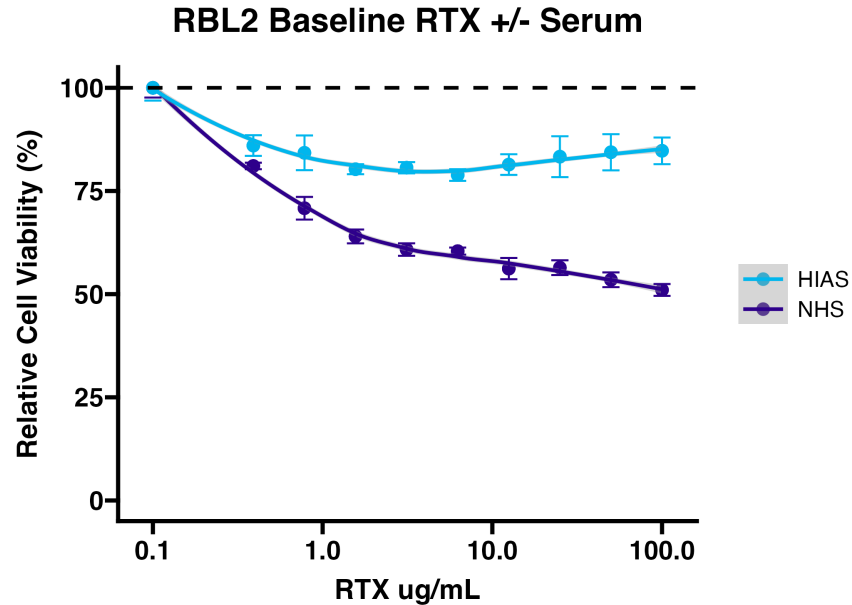

- Conclusions:
 - 40k might be too much
 1. Seed 20k on top and 30k on bottom?
 2. Thaw new cells

Thursday 23-04-2024

Cell Culture

RBL2 Barcoded

- Thawed RBL2 ### NK-92
- Made NK92 media #### Recipe

Solution	ID code	Volume	% Total volume
RMPI 1640 + L-Glutamine	—	500 mL	78%
FBS	—	50 mL	20%
NaPyr	—	5.5 mL	1%
FBS	—	128 mL	20%

- Thawed NK-92


EC50_250422 Collection - RBL2 RTX CDC Baseline

- Collected plates seeded on 22-04-2025

- EC Plate collection protocol:
 1. Added 20 μ L Cell Titre Blue (CTB) to each conditioned well
 - 20 μ L CTB/100 μ L of conditioned well recommended by manufacturer
 2. Incubated for 2hr at 37C
 3. Read on plate reader according to Cell Titre Blue Protocol

Results:

- Seeding error
 - Seeded incorrect number of cells (undercounted - overseeded)
 - HIAS condition is lower than ideal
 1. Seeding error?
 2. HIAS quality?
 - Still struggling to crack 50% Cell viability reduction
 1. In this case likely due to seeding
 2. EC50_250422 looks similar however
 3. Seed less cells?

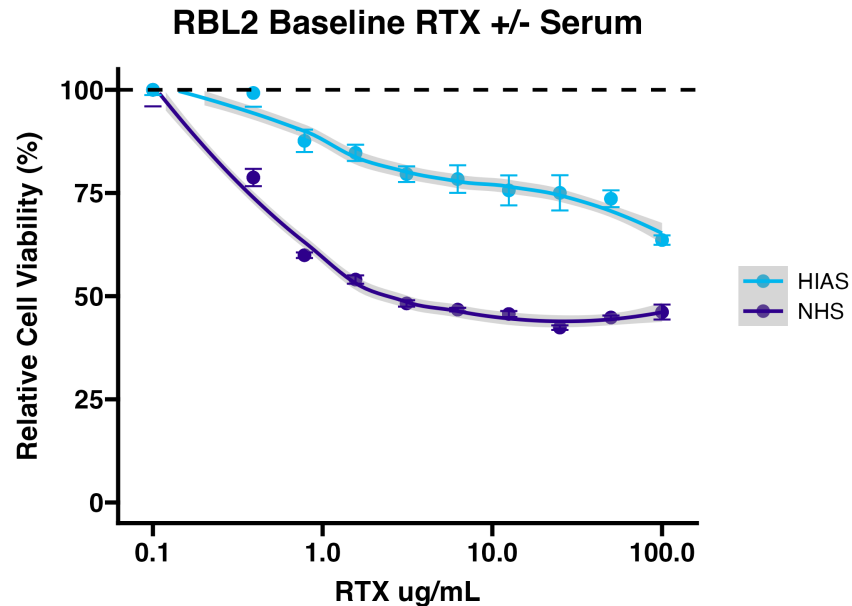

- Conclusions:
 - Re-thaw cells and restart

Monday 28-04-2024

Cell Culture

RBL2 Barcoded

- Split RBL2 ### NK-92
- Add media


EC50_250430 Collection - RBL2 RTX CDC Baseline


- Collected plates seeded on 22-04-2025
- EC Plate collection protocol:
 1. Added 20 μ L Cell Titre Blue (CTB) to each conditioned well

- 20 μ L CTB/100 μ L of conditioned well recommended by manufacturer
- 2. Incubated for 2hr at 37C
- 3. Read on plate reader according to Cell Titre Blue Protocol

Results:

- Seeding error
 - Seeded incorrect number of cells (undercounted - overseeded)
 - HIAS condition is lower than ideal
 1. Seeding error?
 2. HIAS quality?
 - Still struggling to crack 50% Cell viability reduction
 1. In this case likely due to seeding
 2. EC50_250422 looks similar however
 3. Seed less cells?

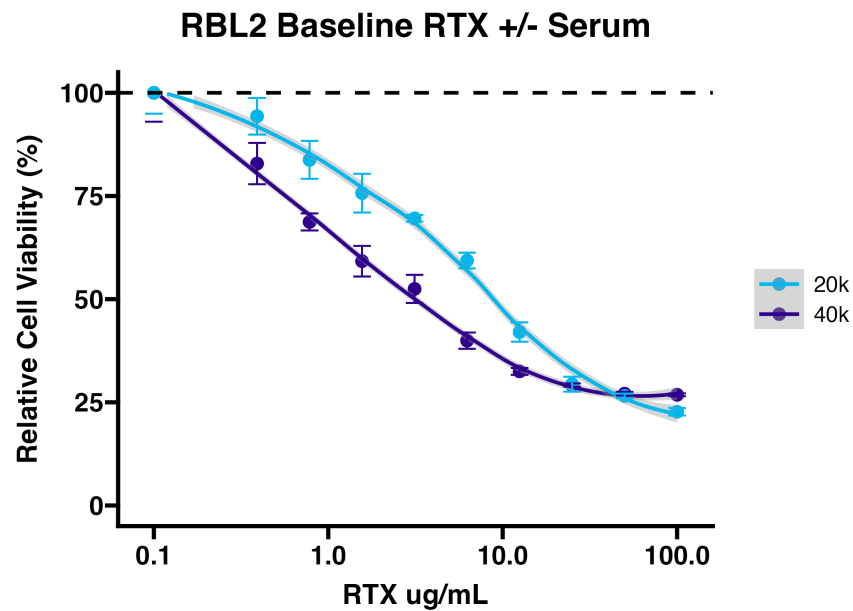

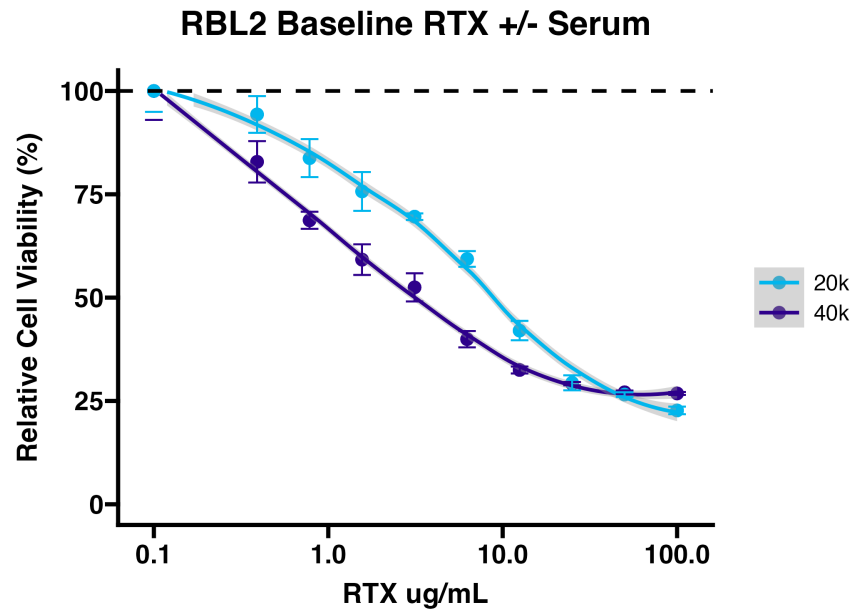

- Conclusions:
 - Re-thaw cells and restart