

# Turner Lab Notebook

Callum Malcolm

## Contents

<b>MARCH 2023</b>	<b>2</b>
<b>Wednesday 15-03-2023</b>	<b>2</b>
Cell Culture . . . . .	2
N2 Growth Curve . . . . .	2
<b>Thursday 16-03-2023</b>	<b>2</b>
N2/BLLW Growth Curve Day 1 . . . . .	2
<b>Friday 17-03-2023</b>	<b>3</b>
N2/BLLW Growth Curve Day 2 . . . . .	3
Split Cells . . . . .	3
<b>Saturday 18-03-2023</b>	<b>3</b>
N2/BLLW Growth Curve Day 3 . . . . .	3
<b>Sunday 18-03-2023</b>	<b>4</b>
N2/BLLW Growth Curve Day 4 . . . . .	4
<b>Monday 20-03-2023</b>	<b>4</b>
N2 Growth Curve Day 5 . . . . .	4
Split Cells . . . . .	4
Seed BLLW . . . . .	4
<b>Tuesday 21-03-2023</b>	<b>4</b>
<b>Friday 24-03-2023</b>	<b>4</b>
Seed Cells - N2/RBL1/JIRE . . . . .	4
<b>Monday 27-03-2023</b>	<b>4</b>
Split Cells . . . . .	4
<b>Wednesday 29-03-2023</b>	<b>5</b>
Split Cells . . . . .	5
<b>Friday 31-03-2023</b>	<b>5</b>
Split Cells . . . . .	5
<b>APRIL 2023</b>	<b>5</b>
<b>Monday 03-04-2023</b>	<b>5</b>
Seed Cells - N2/RBL1/JIRE . . . . .	5
Seed EC50 . . . . .	5

<b>Wednesday 05-04-2023</b>	<b>5</b>
Seed Cells - N2/RBL1/JIRE . . . . .	5
<b>Thursday 06-04-2023</b>	<b>5</b>
Collect EC50 Plate . . . . .	5
<b>Tuesday 11-04-2023</b>	<b>6</b>
Seed EC50_230411 . . . . .	6
<b>Wednesday 12-04-2023</b>	<b>7</b>
Cell Culture . . . . .	7
Side Population - N2_230412 . . . . .	7
<b>Tuesday 02-05-2023</b>	<b>7</b>
Seeding Cells - N2/RBL1/RAMOS . . . . .	7
<b>Wednesday 03-05-2023</b>	<b>7</b>
Split Cells . . . . .	7
<b>Thursday 04-05-2023</b>	<b>7</b>
<b>Monday 19-06-2023</b>	<b>7</b>
Split Cells . . . . .	7
Seed EC50 . . . . .	8
<b>Saturday 20-05-2023</b>	<b>8</b>
<b>JUNE 2023</b>	<b>8</b>
<b>Tuesday 20-06-2023</b>	<b>8</b>
RNA Extraction for SP Test 3 Sample . . . . .	8
Tapestation_230620 . . . . .	10
Reverse Transcription . . . . .	10
PCR Amplification of cDNA . . . . .	11
<b>Wednesday 21-06-2023</b>	<b>12</b>
Tapestation_230619_SPTTest3_cDNA . . . . .	12
<b>Thursday 22-06-2023</b>	<b>12</b>
Collect EC50 Plate . . . . .	12
<b>Friday 23-06-23</b>	<b>12</b>
Split Cells . . . . .	12
Seed EC50 . . . . .	12
<b>Monday 26-06-23</b>	<b>13</b>
Collect EC50 Plate . . . . .	13
<b>JULY 2023</b>	<b>14</b>
<b>Monday 24-07-2023</b>	<b>14</b>
Seeding Cells - RAMOS/FDC . . . . .	14
Made New Media: AR-1 . . . . .	14
<b>Tuesday 25-07-2023</b>	<b>15</b>
Seed Cells - N2-BC4/RBL1/BLLW/RAMOS-BC1 . . . . .	15

<b>Wednesday 26-07-2023</b>	<b>15</b>
Cell Culture . . . . .	15
<b>Thursday 27-07-2023</b>	<b>15</b>
Split Cells - RAMOS . . . . .	15
Lucy Drug Randomizer . . . . .	16
<b>Friday 28-07-2023</b>	<b>16</b>
<b>OCTOBER 2023</b>	<b>16</b>
<b>Friday 06-10-2023</b>	<b>16</b>
Cell Culture . . . . .	16
Made New Media: AR-2 . . . . .	16
<b>Monday 09-10-2023</b>	<b>16</b>
Cell Culture . . . . .	16
Seeded A20 . . . . .	16
Seeded FDC . . . . .	17
<b>Tuesday 10-10-2023</b>	<b>17</b>
Cell Culture . . . . .	17
Seeded N4 . . . . .	17
Seeded RAMOS-BC . . . . .	17
Checked A20 . . . . .	17
<b>Wednesday 11-10-2023</b>	<b>17</b>
Cell Culture . . . . .	17
A20 Culture . . . . .	17
RAMOS-BC4 Culture . . . . .	17
N4 Culture . . . . .	17
<b>Thursday 12-10-2023</b>	<b>18</b>
Cell Culture . . . . .	18
Ramos-BC4 . . . . .	18
A20 . . . . .	18
N4 . . . . .	18
<b>Friday 13-10-2023</b>	<b>18</b>
Cell Culture . . . . .	18
A20 . . . . .	18
EC50_231013 . . . . .	18
EC50 Plate ID: EC50_231013_Plate1 . . . . .	18
EC50 Plate ID: EC50_231013_Plate2 . . . . .	19
<b>Monday 16-10-2023</b>	<b>20</b>
Cell Culture . . . . .	20
A20 . . . . .	20
N4 . . . . .	20
Collect EC50_231013 . . . . .	20
Plate ID: EC50_231013_Plate1 . . . . .	20
Plate ID: EC50_231013_Plate2 . . . . .	22
<b>Thursday 19-10-2023</b>	<b>22</b>
Cell Culture . . . . .	22
A20 . . . . .	22

N4 . . . . .	22
<b>Friday 20-10-2023</b>	<b>22</b>
Cell Culture . . . . .	22
A20 . . . . .	22
N4 . . . . .	22
<b>Sunday 22-10-2023</b>	<b>23</b>
Cell Culture . . . . .	23
RAMOS Barcode . . . . .	23
<b>Monday 23-10-2023</b>	<b>23</b>
Cell Culture . . . . .	23
N4 . . . . .	23
A20 . . . . .	23
Ramos-BC 1 . . . . .	23
Ramos-BC 5 . . . . .	23
Ramos-BC 6 . . . . .	23
Seed Feeders (FDC) . . . . .	23
<b>Tuesday 23-10-2023</b>	<b>23</b>
Cell Culture . . . . .	23
N4 . . . . .	23
A20 . . . . .	24
Ramos-BC 1 . . . . .	24
Ramos-BC 5 . . . . .	24
Ramos-BC 6 . . . . .	24
FDC . . . . .	24
Transduce A20/N4 with WILDseq Library . . . . .	24
<b>Wednesday 25-10-23</b>	<b>24</b>
Attempt 1: A20/N4 Barcoding - Media change . . . . .	24
Made New Media: AR-3 . . . . .	25
<b>Thursday 26-10-23</b>	<b>25</b>
Attempt 1: A20/N4 Barcoding - Sort . . . . .	25
<b>Friday 27-10-23</b>	<b>25</b>
Cell Culture . . . . .	25
A20 . . . . .	25
N4 . . . . .	25
Ramos-BC 1 . . . . .	25
Ramos-BC 5 . . . . .	25
Ramos-BC 6 . . . . .	25
EC50_231027 . . . . .	26
Overview . . . . .	26
EC50 Protocol: . . . . .	26
EC50 Plate ID: EC50_231027_Plate1 . . . . .	26
EC50 Plate ID: EC50_231027_Plate2 . . . . .	27
<b>Monday 27-10-23</b>	<b>28</b>
Cell Culture . . . . .	28
A20 . . . . .	28
N4 . . . . .	28
Ramos-BC 1 . . . . .	28

Ramos-BC 5 . . . . .	29
Ramos-BC 6 . . . . .	29
Seed Feeders . . . . .	29
Collect EC50_231027 . . . . .	29
Plate ID: EC50_231027_Plate . . . . .	29
<b>Tuesday 31-10-2023</b>	<b>29</b>
Attempt 2: Transduce A20/N4 with WILDseq Library . . . . .	29
Drug Pressure Experiment 1: RAMOS Day 1 . . . . .	30
Overview . . . . .	30
Day 1 . . . . .	30
Seeded RAMOS BC . . . . .	30
Collected RAMOS Baseline . . . . .	31
<b>NOVEMBER 2023</b>	<b>31</b>
<b>Wednesday 01-11-2023</b>	<b>31</b>
Drug Pressure Experiment 1: RAMOS Day 2 . . . . .	31
Attempt 2: Barcoding A20/N4 . . . . .	32
<b>Thursday 02-11-2023</b>	<b>32</b>
Cell Culture . . . . .	32
A20 . . . . .	32
N4 . . . . .	32
Ramos-BC 1 . . . . .	32
Ramos-BC 5 . . . . .	32
Ramos-BC 6 . . . . .	32
Seeded Feeders . . . . .	32
Drug Pressure Experiment 1: RAMOS Day 3 . . . . .	32
<b>Friday 03-11-2023</b>	<b>33</b>
Cell Culture . . . . .	33
A20 . . . . .	33
N4 . . . . .	33
Ramos-BC 1 . . . . .	33
Ramos-BC 5 . . . . .	33
Ramos-BC 6 . . . . .	33
Drug Pressure Experiment 1: RAMOS Day 4 . . . . .	33
<b>Monday 06-11-2023</b>	<b>33</b>
Cell Culture . . . . .	33
A20 . . . . .	33
N4.1 . . . . .	33
N4.2 . . . . .	34
Ramos-BC 1 . . . . .	34
Ramos-BC 5 . . . . .	34
Ramos-BC 6 . . . . .	34
Seed Feeders . . . . .	34
Drug Pressure Experiment 1: RAMOS Day 7 . . . . .	34
<b>Tuesday 07/11/2023</b>	<b>34</b>
Made New Media: AR-4 . . . . .	34
Cell Culture . . . . .	34
A20 . . . . .	34
N4.1 . . . . .	35

N4.2 . . . . .	35
Drug Pressure Experiment 1: RAMOS Day 5 . . . . .	35
Attempt 3: Transduce N4 with WILDseq Library . . . . .	35
<b>Wednesday 08/11/2023</b>	<b>35</b>
Attempt 3: Barcoding N4 . . . . .	35
<b>Thursday 09/11/2023</b>	<b>35</b>
Attempt 3: Barcoding N4 . . . . .	35
A20 . . . . .	36
<b>Friday 10/11/2023</b>	<b>36</b>
Attempt 3: Barcoding N4 . . . . .	36
Cell Culture . . . . .	36
RAMOS-BC1 . . . . .	36
<b>Monday 13/11/2023</b>	<b>36</b>
Start RAMOS-Drug Pressure v2 . . . . .	36
Cell Culture . . . . .	38
N4-BC Attempt 3 . . . . .	38
<b>Thursday 13/11/2023</b>	<b>38</b>
RAMOS-DPv2 . . . . .	38
<b>Friday 17/11/2023</b>	<b>38</b>
Cell Culture . . . . .	38
<b>Monday 20/11/2023</b>	<b>38</b>
WILDseq - Mouse Experiment 1 . . . . .	38
Overview . . . . .	38

**MARCH 2023**

**Wednesday 15-03-2023**

### Cell Culture


- Given N2/RBL1/BLLW PDX line from Chris Steele

### N2 Growth Curve

Looking to compare N2 growth on and off feeders

1. Seeded 2mL of N2/BLLW cell suspension at  $0.25 \times 10^6$  cells/mL

## 96 Well Plate Map

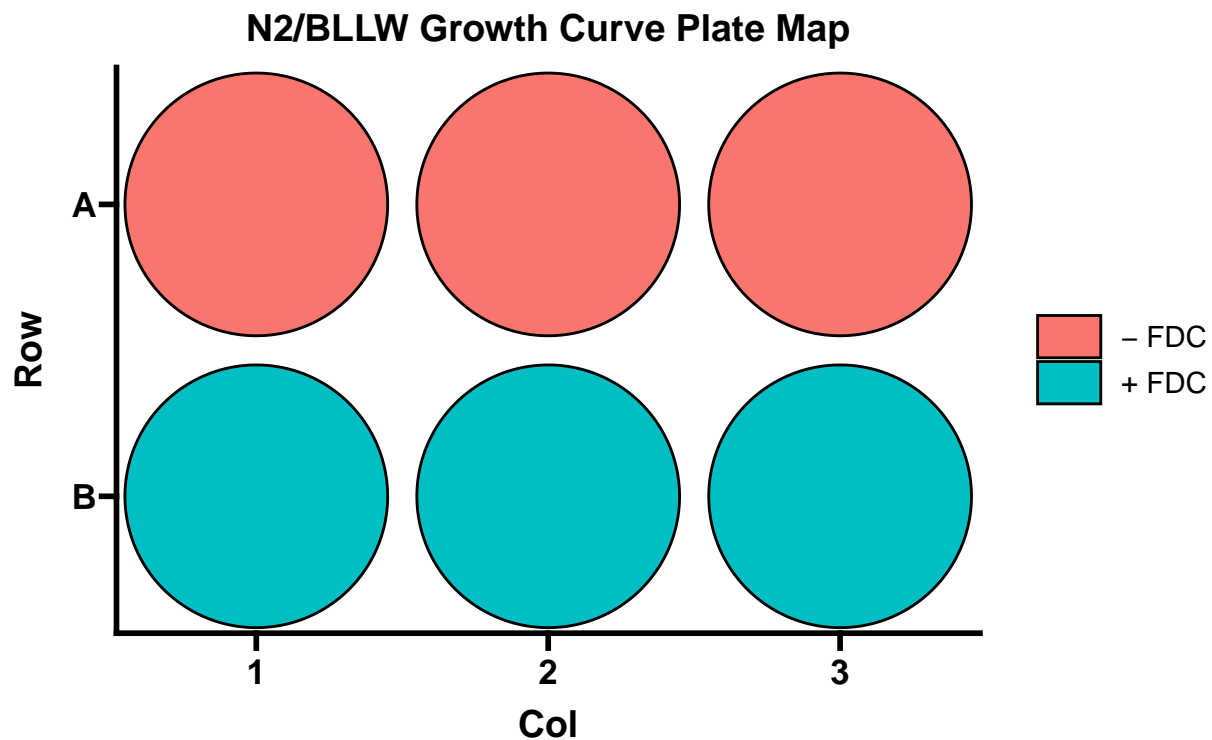

**Thursday 16-03-2023**

### N2/BLLW Growth Curve Day 1

- Cells gently mixed
- 30  $\mu$ L cell suspension taken and counted
- Only counted 1 square and multiplied by 4

Cell Line	FDC Status	Count	Cell Number
N2	+ FDC	25	2000000
N2	+ FDC	26	2080000
N2	+ FDC	12	960000
N2	- FDC	15	1200000
N2	- FDC	14	1120000
N2	- FDC	6	480000
BLLW	+ FDC	6	480000
BLLW	+ FDC	6	480000
BLLW	+ FDC	5	400000
BLLW	- FDC	3	240000
BLLW	- FDC	3	240000
BLLW	- FDC	6	480000

## Friday 17-03-2023

### N2/BLLW Growth Curve Day 2

- Cells gently mixed
- 30  $\mu$ L cell suspension taken and counted
- Only counted 1 square and multiplied by 4

Cell Line	FDC Status	Count	Cell Number
N2	+ FDC	36	2880000
N2	+ FDC	28	2240000
N2	+ FDC	38	3040000
N2	- FDC	61	4880000
N2	- FDC	38	3040000
N2	- FDC	16	1280000
BLLW	+ FDC	3	240000
BLLW	+ FDC	5	400000
BLLW	+ FDC	5	400000
BLLW	- FDC	1	80000
BLLW	- FDC	0	0
BLLW	- FDC	1	80000

### Split Cells

Cell Line	Cells/mL	Cells/mL Seeded	CS volume	Final Volume	Seeded Cells/mL
N2	$8.8 \times 10^5$	$3 \times 10^5$	6.8 mL	12 mL	$3 \times 10^5$
RBL1	$8.6 \times 10^5$	$3 \times 10^5$	11.6 mL	12 mL	$5 \times 10^5$

## Saturday 18-03-2023

### N2/BLLW Growth Curve Day 3

- Counted by Chris



## Sunday 18-03-2023

### N2/BLLW Growth Curve Day 4

- Counted by Chris

## Monday 20-03-2023

### N2 Growth Curve Day 5

- BLLW were all dead, experiment stopped
- Cells gently mixed
- 30  $\mu$ L cell suspension taken and counted
- Only counted 1 square and multiplied by 4

Cell Line	FDC Status	Count	Cell Number
N2	+ FDC	146	11680000
N2	+ FDC	112	8960000
N2	+ FDC	132	10560000
N2	- FDC	112	8960000
N2	- FDC	155	12400000
N2	- FDC	141	11280000

### Split Cells

- Split cell lines N2/RBL1/JIRE
- Dirty split 1:2

### Seed BLLW

- Given BLLW aliquot from Chris Steel
- Thawed and resuspended in 12 mL media
- Added to T75 pre-seeded with FDC cells (from Chris Steel)

## Tuesday 21-03-2023

### Large Infection in TC

- N2 Growth Curve stopped

## Friday 24-03-2023

### Seed Cells - N2/RBL1/JIRE

- Vials given by Chris Steel
- Thawed and resuspended in 12 mL media
- Added to T75 pre-seeded with FDC cells (from Chris Steel)

## Monday 27-03-2023

### Split Cells

- Split cell lines N2/RBL1/JIRE
- Dirty split 1:3

## **Wednesday 29-03-2023**

### **Split Cells**

- Split by Chris Steel.

## **Friday 31-03-2023**

### **Split Cells**

- Split cell lines N2/RBL1/JIRE
- Dirty split 1:4

## **APRIL 2023**

## **Monday 03-04-2023**

### **Large Infection in TC**

### **Seed Cells - N2/RBL1/JIRE**

- Vials given by Chris Steel
- Thawed and resuspended in 12 mL media
- Added to T75 pre-seeded with FDC cells (from Chris Steel)

### **Seed EC50**

DNW - N2/JIRE/RBL1  
- Vin/Dactolisib

## **Wednesday 05-04-2023**

### **Seed Cells - N2/RBL1/JIRE**

- Vials given by Chris Steel
- Thawed and resuspended in 12 mL media
- Added to T75 pre-seeded with FDC cells (from Chris Steel)

## **Thursday 06-04-2023**


### **Collect EC50 Plate**

- Plate ID: EC50\_230619\_Plate1  
- DNW

Tuesday 11-04-2023

Seed EC50\_230411

## 96 Well Plate Map

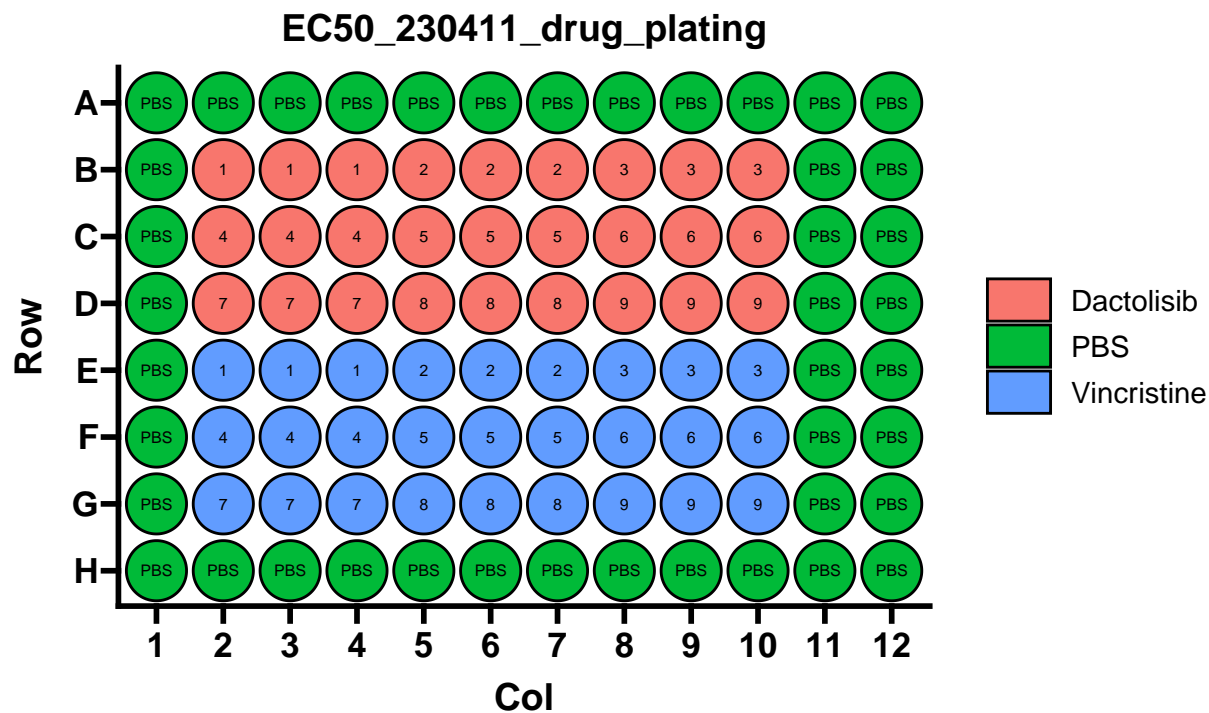

### Drug Dilutions

Drug	Dilution	Concentration	Units	Starting_Concentration
Dactolisib	1	20	uM	20000
Dactolisib	2	2	uM	2000
Dactolisib	3	0.2	uM	200
Dactolisib	4	0.02	uM	20
Dactolisib	5	0.002	uM	2
Dactolisib	6	0.0002	uM	0.2
Dactolisib	7	2e-05	uM	0.02
Dactolisib	8	2e-06	uM	0.002
Dactolisib	9	2e-07	uM	0.0002

Drug	Dilution	Concentration	Units	Starting_Concentration
Vincristine	1	30	nM	30000
Vincristine	2	3	nM	3000
Vincristine	3	0.3	nM	300
Vincristine	4	0.03	nM	30
Vincristine	5	0.003	nM	3
Vincristine	6	0.0003	nM	0.3
Vincristine	7	3e-05	nM	0.03
Vincristine	8	3e-06	nM	0.003
Vincristine	9	3e-07	nM	0.0003

- 3 plates seeded
  - Plate 1: N2

- Plate 2: JIRE
- Plate 3: RBL1

## Wednesday 12-04-2023

### Cell Culture

- Thaw RBL1

### Side Population - N2\_230412

- Hoescht
  - 2.5 / 2.75 / 3
- DNW

# MAY 2023

## Tuesday 02-05-2023

### Seeding Cells - N2/RBL1/RAMOS

1. Thawed N2, RBL1, and RAMOS into T75 flasks  
Note: N2 and RBL1 flasks seeded with feeders

## Wednesday 03-05-2023

### Split Cells

Cell Line	Cells/mL	Cells/mL Seeded	CS volume	Final Volume	Last day split	P number
N2	8.93x10 <sup>5</sup>	3x10 <sup>5</sup> .	6.71 mL	12 mL	Tues 02-05-23	2
RBL1	7.73x10 <sup>5</sup>	3x10 <sup>5</sup> .	7.76 mL	12 mL	Tues 02-05-23	2

2. Expanded N2-Barcodes to T175
3. Seeded feeders into 96-well
  - a. Diluted feeders to 3.2x10<sup>4</sup> cells/mL
  - b. Add 100μL cell suspension to each well

## Thursday 04-05-2023

## Monday 19-06-2023

### Split Cells


- Split N2/RBL1/JIRE
- Added 10 mL media to BLLW
  - Lots of dead cells observed

## Seed EC50

- Plate ID: EC50\_230619\_Plate1
- Cyclophosphamide (Stock: 20mM) EC50 for N2 and RBL1

## 96 Well Plate Map

### EC50\_230619\_plate1

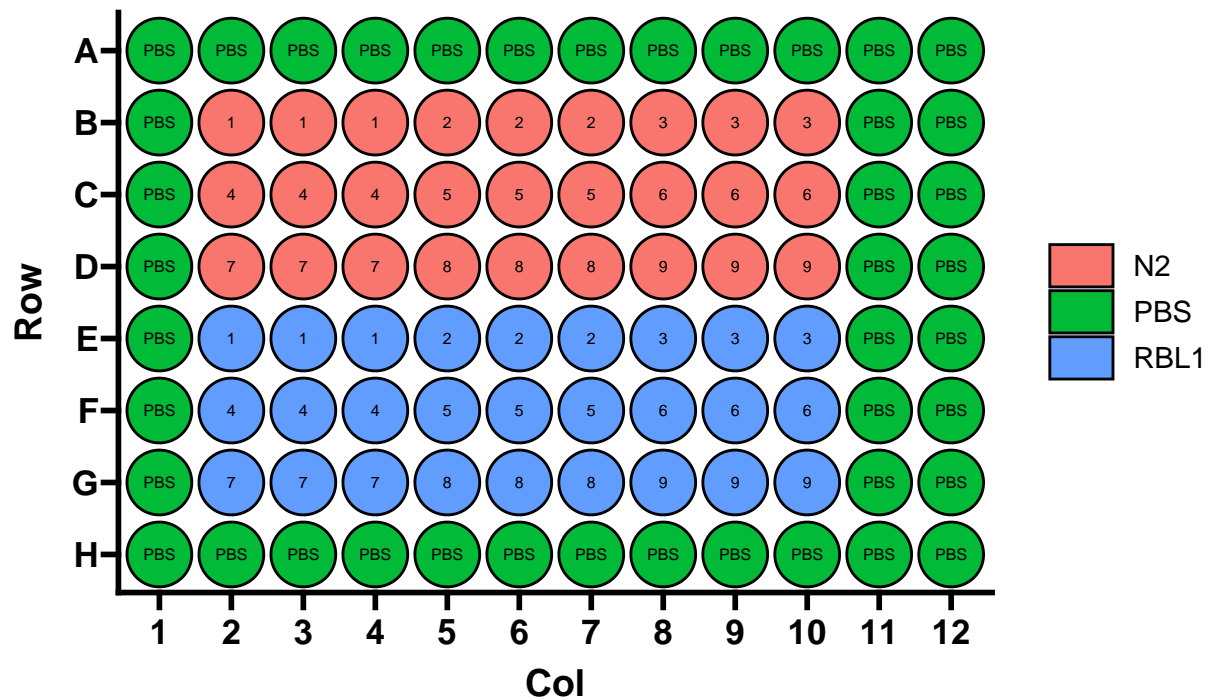

## Drug Dilutions

Drug	Dilution	Concentration	Units	Starting_Concentration
Cyclophosphamide	1	20	uM	20000
Cyclophosphamide	2	2	uM	2000
Cyclophosphamide	3	0.2	uM	200
Cyclophosphamide	4	0.02	uM	20
Cyclophosphamide	5	0.002	uM	2
Cyclophosphamide	6	0.0002	uM	0.2
Cyclophosphamide	7	2e-05	uM	0.02
Cyclophosphamide	8	2e-06	uM	0.002
Cyclophosphamide	9	2e-07	uM	0.0002

Saturday 20-05-2023

JUNE 2023

Tuesday 20-06-2023

## RNA Extraction for SP Test 3 Sample

- Sample collected:

- Low number of cells so used Pico-Pure Kit
- 1 sample (n = 1)

*Note before starting:*

- Pre-heat block to 42C
- Pre-make DNase treatment master mix (1.1n)

Component	Volume	MM volume
DNase I Stock	5 $\mu$ L	5.5 $\mu$ L
DNA Digestion Buffer	35 $\mu$ L	38.5 $\mu$ L

### **RNA Extraction Pico-Pure Protocol**

1. Extract cells with 100  $\mu$ L of Extraction Buffer (XB)
  - Resuspend the cell pellet gently by pipetting
  - DO NOT VORTEX
2. Incubate at 42C for 30min
3. Centrifuge sample at 3000xg for two minutes
4. Collect supernatant
  - Can stop here and freeze RNA at 70C
5. Pre-condition the RNA Purification Column
  - Add 250  $\mu$ L Conditioning Buffer onto the purification column filter membrane
  - Incubate RNA Purification Column with CB for 5min at RT
  - Centrifuge purification column at 16000xg for 1 min
6. Pipette 100  $\mu$ L of 70% Ethanol (EtOH) into cell extract
  - Mix well by pipetting
  - DO NOT VORTEX
7. Add cell extract to column
8. Centrifuge column for 2 minutes at 100 x g and then 16000 x g for 30s
9. Pipette 100  $\mu$ L Wash Buffer 1 (W1) into the purification column and centrifuge for 1min at 8000 x g
10. Add 40  $\mu$ L of DNase treatment master mix to purification tube
  - Add directly onto membrane
11. Incubate at RT for 15 min
12. Pipette 40  $\mu$ L Wash Buffer 1 (W1) onto purification column membrane and centrifuge at 8000 x g for 15 seconds
13. Pipette 100  $\mu$ L Wash Buffer 1 (W2) into purification column membrane and centrifuge for one minute at 8000 x g
14. Pipette 100  $\mu$ L Wash Buffer 1 (W2) into purification column membrane and centrifuge for two minutes at 16000 x g
  - Check if any wash buffer remains in column and centrifuge again at 16000 x g for one minute to clear

all the liquid

15. Transfer purification column to new 0.5mL tube
  - Insert purification column/0.5mL assembly into 1.5mL tube for centrifugation
16. Pipette EB directly onto membrane of the purification column
  - Gently touch tip to surface of membrane to ensure maximum absorption of EB into membrane
  - Use 11  $\mu\text{L}$  to 30  $\mu\text{L}$  EB
17. Incubate purification column for one minute at room temperature
18. Centrifuge column for one minute at 1000 x g to distribute EB in column
19. Centrifuge column for one minute at 16000 x g to elute RNA
20. Started PCR and completed Tapestation

## Tapestation\_230620

### High-Sensitivity RNA ScreenTape Protocol

#### Results

- Inconclusive, proceeded with PCR

## Reverse Transcription

### Notes Before Starting:

- Label PCR strip tubes
- Thaw RT reagents
- Uses WILDseq specific RT primer
- Dilute RT primer to 2 $\mu\text{M}$ 
  - Do this in a PCR hood
  - **IMPORTANT** Do not mix primers and carefully note which RT primer is added to which sample
- Make RT and ExoI master mixes (MM) before starting
  - MM amounts are 1.1n

### Protocol:

*Note:* Used “rt” protocol in callum folder on thermocycler for heating steps

1. In PCR strip tubes diluted 1-5ug RNA in a total volume of 10  $\mu\text{L}$  of RNase/DNase-free water
2. Added 1 ul of RT primer
3. Added 1 ul of dNTPs
4. Denatured at 65 C for 5 mins in PCR machine
5. Added straight to ice for 2 mins
6. Spun to get liquid to bottom of tube
7. Prepared master mix of RT enzyme and buffers (Made prior to starting)

Component	Volume	MM volume
5x SSIV Buffer	4 $\mu\text{L}$	4.4 $\mu\text{L}$

Component	Volume	MM volume
SSIV RT	1 $\mu$ L	1.1 $\mu$ L
100 mM DTT	1 $\mu$ L	1.1 $\mu$ L
RNAse Out	1 $\mu$ L	1.1 $\mu$ L

8. Added 7  $\mu$ L of RT MM prepared above to each sample and mix
  - Spun briefly to get liquid to bottom of the tube
9. In PCR Machine: incubated at 53C for 10 mins followed by 80C for 10 mins
10. Added 3  $\mu$ L of ExoI MM to each sample

Component	Volume	MM volume
Thermolabile Exonuclease I	1 $\mu$ L	1.1 $\mu$ L
NEBuffer r3.1*	2 $\mu$ L	2.2 $\mu$ L

- Most PCR buffers are compatible
11. In PCR Machine: heat at 37C for 4 mins followed by 80C for 1 min
    - Spin briefly to get liquid to bottom of the tube
  12. Add 1  $\mu$ L of RNAse H and incubate at 37C for 20 mins
  13. Dilute final cDNA 1 in 2 with DNase/RNAse-free water

## PCR Amplification of cDNA

1. Prepare mastermix of below reagents at 1.1x the number of samples

Component	Volume	MM volume
Q5 Master Mix	25 $\mu$ L	
P5_R1_fwd (100 $\mu$ M)	2.5 $\mu$ L	
P7_rev	2.5 $\mu$ L	
DNase/RNAse H20	13.5 $\mu$ L	
Q5 polymerase	0.5 $\mu$ L	
<b>Total Master Mix volume</b>	<b>40 <math>\mu</math>L</b>	
cDNA from above	10 $\mu$ L	

2. Perform PCR using the following parameters:

Steps	Time
Step 1: 98C	30s
<b>Step 2: 98C</b>	<b>10s</b>
<b>Step 3: 61C</b>	<b>30s</b>
<b>Step 4: 72C</b>	<b>30s</b>
Step 5: 72C	2 min
Step 6: 12C	Hold

**20-25 cycles of steps 2-4**

Ran overnight and left at 12C



**Wednesday 21-06-2023**

**Tapestation\_230619\_SPTest3\_cDNA**

High-Sensitivity DNA ScreenTape Protocol



**Results**

- \* Did not work, bands present but faint and not at 200bp
- \* Potential solutions:
  - Repeat PCR with increased cycles
  - Repeat RNA extraction with trizol protocol (Anna)

**Thursday 22-06-2023**

**Collect EC50 Plate**

- Plate ID: EC50\_230619\_Plate1

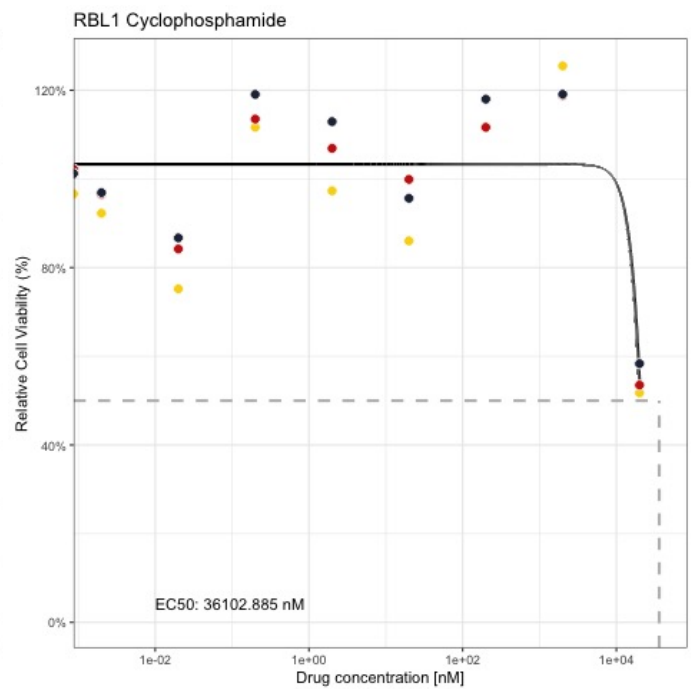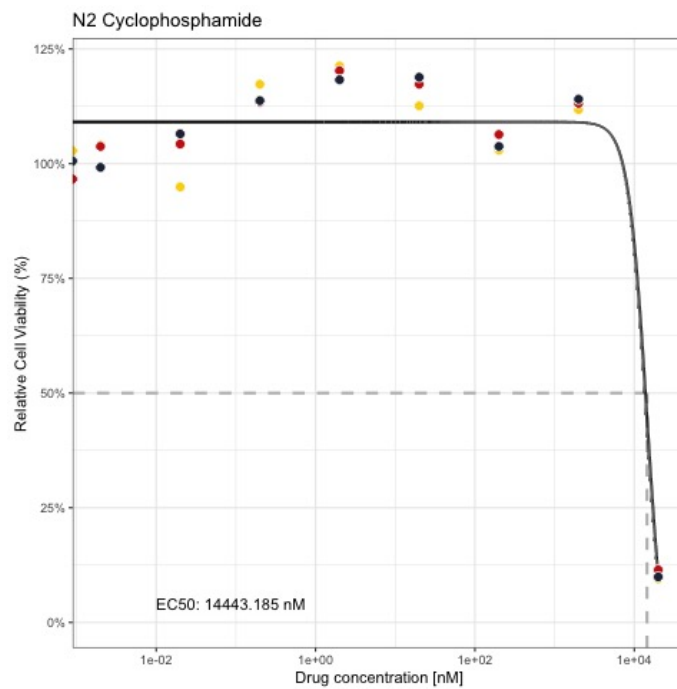

**Friday 23-06-23**

**Split Cells**


1. Threw out BLLW (all dead)
2. Split JIRE 1/2
3. Split N2/RBL1 1/4

**Seed EC50**

- EC50 Plate ID: EC50\_230623\_Plate1

# 96 Well Plate Map

## EC50\_230623\_plate1

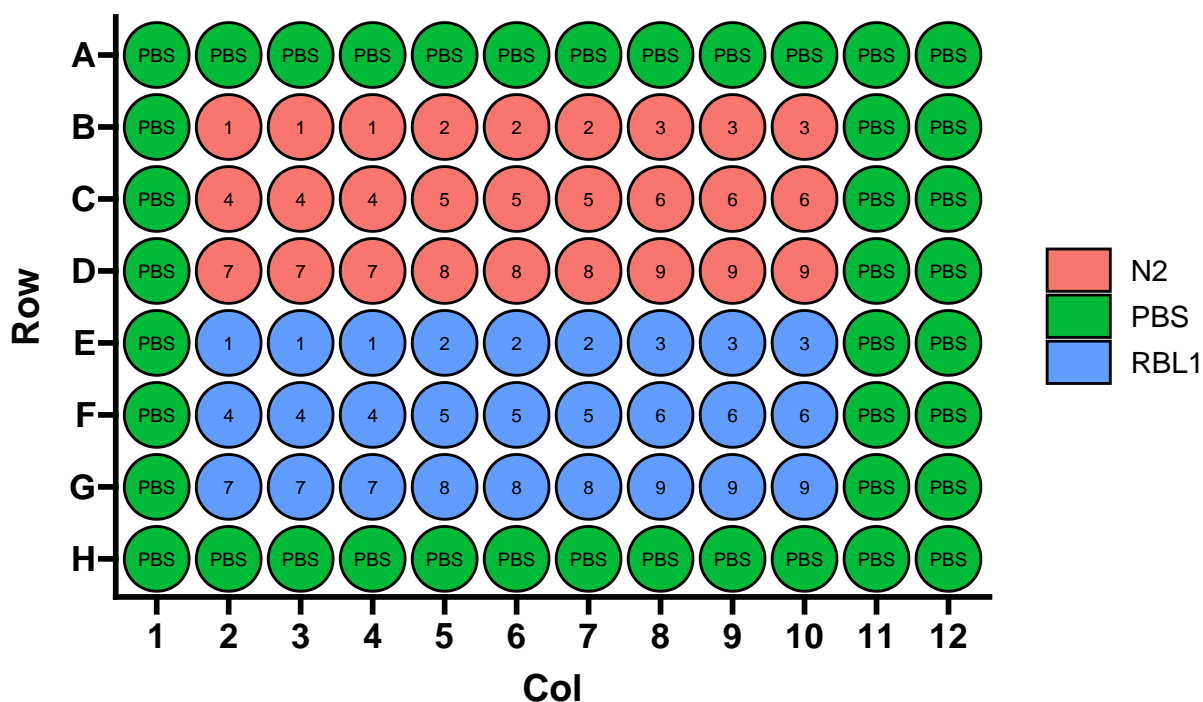


### Drug Dilutions

Drug	Dilution	Concentration	Units	Starting_Concentration
Cyclophosphamide	1	20	uM	20000
Cyclophosphamide	2	13.3	uM	13300
Cyclophosphamide	3	8.89	uM	8890
Cyclophosphamide	4	5.93	uM	5930
Cyclophosphamide	5	3.95	uM	3950
Cyclophosphamide	6	2.63	uM	2630
Cyclophosphamide	7	1.76	uM	1760
Cyclophosphamide	8	1.17	uM	1170
Cyclophosphamide	9	0.78	uM	780

Monday 26-06-23

### Collect EC50 Plate

- Plate ID: EC50\_230623\_Plate1

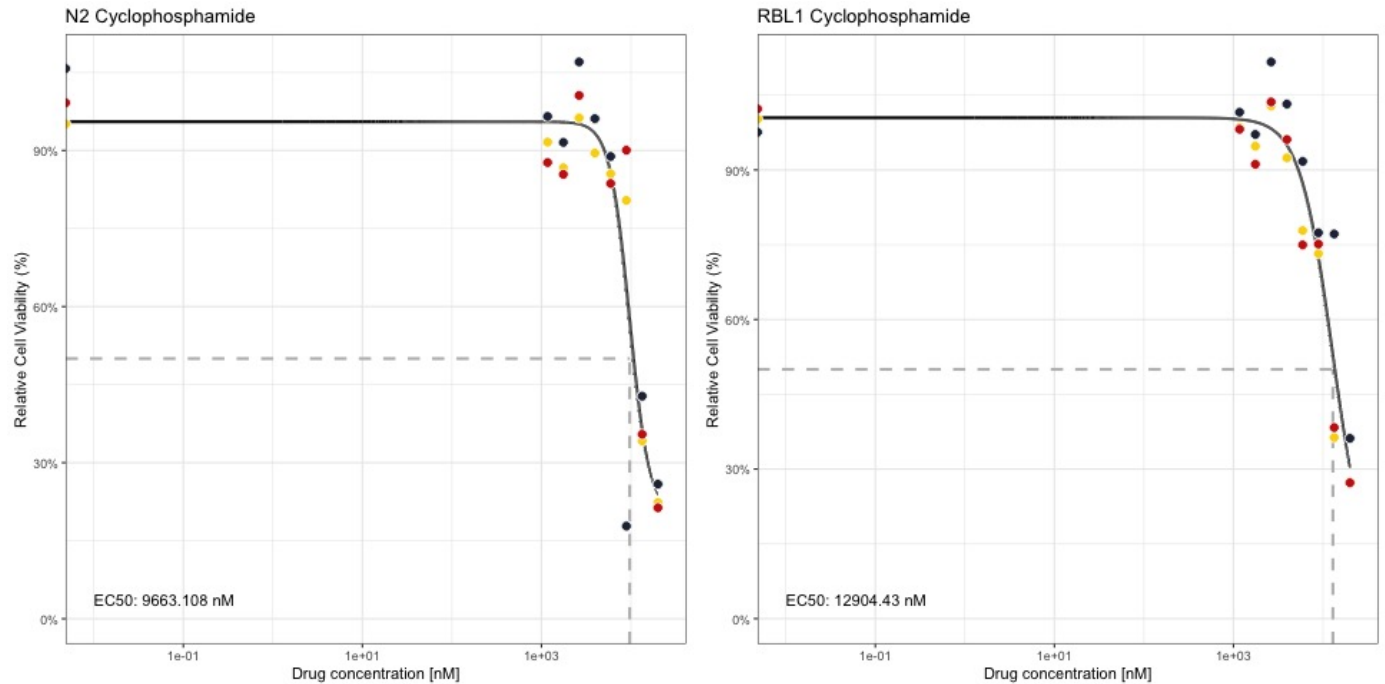

Break for writing MRes Report: 27/06/2023 - 24/07/2023

**JULY 2023**

**Monday 24-07-2023**

**Seeding Cells - RAMOS/FDC**

**Made New Media: AR-1**

1. Made new media: AR-1

Solution	ID code	Volume	% Total volume
Advanced RMPI	—	500 mL	78%
Glutamax	—	6.5 mL	20%
Pen-Strep	—	6.5 mL	1%
FBS	—	128 mL	1%

2. Seeded Feeders

- 2 vial of FDC cells thawed
- Resuspended in 12 mL of media
- 6 mL of FDC cell suspension were added to 2 separate T75 flasks (Seeding RBL1/BLLW on 25-07-2023)
- Remaining 12 mL of Media were added to T175 (Seeding N2 on 25-07-2023)

3. Seeded Ramos WILDseq Barcode Pool 2

- Thawed a vial of RAMOS pool 2 (previously barcoded DATE)
- Resuspended in 5 mL
- Spun down 90\* for 5min

- Resuspended in 12 mL of AR-1 media
- Seeded into T75

### Experimental Plans

1. Barcode BLLW/RBL1
2. Set up Drug pressure experiment for RAMOS barcoded cells
  - EC50's for cyclophosphamide and methotrexate need to be completed
  - Drugs to test (EC50):

Drug	Concentration	Replicates
DMSO	1%	3
Doxorubicin	200nM	3
Methotrexate	100nM	3
Vincristine	1nM	3

3. Do Side Population Experiment for N2 Barcoded Cells
  - Grow up sorted side pop to avoid issues with low cell numbers

## Tuesday 25-07-2023

### Seed Cells - N2-BC4/RBL1/BLLW/RAMOS-BC1

- Got vials from Chris
- Thawed, resuspended in 6 mL AR-1
- Spun down at 90\* for 5 min
- Discarded supernatant and resuspended pellet in 12 mL of AR-1
- Added RBL1/BLLW to T75's pre-seeded with FDC (24-07-2023)
- Added RAMOS-BC1 to new T75
  - RAMOS-BC2 from 24-07-2023 was dead
  - Was taken from -80, New RAMOS vial from LN
- Added N2 to T125 seeded with FDC (24-07-2023)

## Wednesday 26-07-2023

### Cell Culture

- Checked cells, all flasks looked healthy
- Added media
  - N2 - 20ml of AR-1
  - RAMOS - 10 mL of AR-1
  - RBL1 - 10ml of AR-1
  - BLLW - 10ml of AR-1

## Thursday 27-07-2023

### Split Cells - RAMOS

- Split Ramos 1:2
- Added 10ml of AR-1 to N2
- Plan to split RBL1/BLLW tomorrow

## Lucy Drug Randomizer

- Randomized drugs for Lucy Hare mouse experiment

Drug	Identification
Alectinib	A
Erda	B
Vehicle	C
Alectinib + Erda	D

## Friday 28-07-2023

Drug	Dilution	Concentration	Units	Starting Concentration
Cyclophosphamide	1	20	uM	20000
Cyclophosphamide	2	13.3	uM	13300
Cyclophosphamide	3	8.89	uM	8890
Cyclophosphamide	4	5.93	uM	5930
Cyclophosphamide	5	3.95	uM	3950
Cyclophosphamide	6	2.63	uM	2630
Cyclophosphamide	7	1.76	uM	1760
Cyclophosphamide	8	1.17	uM	1170
Cyclophosphamide	9	0.78	uM	780

## OCTOBER 2023

### Friday 06-10-2023

#### Cell Culture

##### Made New Media: AR-2

- Made new media: AR-2

Solution	ID code	Volume	% Total volume
Advanced RMPI	—	500 mL	78%
Glutamax	—	6.5 mL	20%
Pen-Strep	—	6.5 mL	1%
FBS	—	128 mL	1%

## Monday 09-10-2023

#### Cell Culture

##### Seeded A20

- Got A20 from Swetha box: 2nd rack in shared -80
- AR-2

### **Seeded FDC**

- Seeded feeders, taken from feeder box, irradiated by Chris
- Will seed N4 PDX tomorrow

## **Tuesday 10-10-2023**

### **Cell Culture**

#### **Seeded N4**

- Given by Chris -AR-2
- Seeded onto FDC cells from 09-10
- Future experiments:
  - Barcode PDX
  - EC50: Vincristine/Doxorubicin/Methotrexate/Cyclophosphamide (Active)

#### **Seeded RAMOS-BC**

- Seeded barcoded Ramos cells - Pool #4
  - From Chris
  - LN tank 1 - Rack 5 - box 2nd from Bottom
- AR-2
- Future Experiments:
  - Combination Drug pressure: Vin/Dox/Meth

#### **Checked A20**

- Look healthy
- ~40% confluency

## **Wednesday 11-10-2023**

### **Cell Culture**

#### **A20 Culture**

- ~70% confluent
- Split A20 in half
  - Half in 1xT75 and other half in 1xT175
  - T175 to be frozen down

#### **RAMOS-BC4 Culture**

- ~60% confluent
- Spun down at 400xg for 5 min and resuspended in 20ml AR-2

#### **N4 Culture**

- Look healthy, attached to feeders with some floaters
- Will add media tomorrow
- Plans:
  - Freeze down aliquots (2-3x10<sup>6</sup> per vial needed, ideally freeze down 3x vials)
  - Plan barcoding experiment

## Thursday 12-10-2023

### Cell Culture

#### Ramos-BC4

- Infected, will need to grab another vial
- Is the Stock infected? Will thaw 2 different pools this time

#### A20

- 1 x T75 + 1 x T175 both about ~60% confluent
- Media is slightly turbid, Jamie did a visual check and cells seem clean (A20 have a tendency towards slight turbidity - Jamie)
- Will add 10mL of AR-2 and observe till tomorrow
- If clean:
  - Freeze down 3 vials from T175
  - Set up EC50s with left over cells (Vin/Dox/Meth)
  - Split

#### N4

- Look healthy, media clean
- About ~40-50% confluent
- Will add 10mL of AR-2 and observe on Friday

## Friday 13-10-2023

### Cell Culture

#### A20

- Freeze down (10% DMSO in FBS)
- A20 count:  $2.91 \times 10^6$ /mL
  - 2mL + 3mL Media
- 3 vials frozen and put in LN Tank 6, Rack 5, Box 6

#### EC50\_231013

##### EC50 Plate ID: EC50\_231013\_Plate1

- Doxorubicin top/Vincristine bottom

Doxorubicin Stock: 10mM DMSO


Drug	Dilution	Concentration	Units	Starting_Concentration
Doxorubicin	1	300	nM	300000
Doxorubicin	2	200	nM	200000
Doxorubicin	3	133	nM	133000
Doxorubicin	4	88.9	nM	88900
Doxorubicin	5	59.3	nM	59300
Doxorubicin	6	39.5	nM	39500
Doxorubicin	7	26.3	nM	26300
Doxorubicin	8	17.6	nM	17600
Doxorubicin	9	11.7	nM	11700

Vincristine stock: 30 $\mu$ M

Drug	Dilution	Concentration	Units	Starting_Concentration
Vincristine	1	10	nM	10000
Vincristine	2	7.14	nM	7140
Vincristine	3	5.1	nM	5100
Vincristine	4	3.64	nM	3640
Vincristine	5	2.6	nM	2600
Vincristine	6	1.86	nM	1860
Vincristine	7	1.33	nM	1330
Vincristine	8	0.949	nM	949
Vincristine	9	0.678	nM	678

## 96 Well Plate Map

### EC50\_231013\_plate1

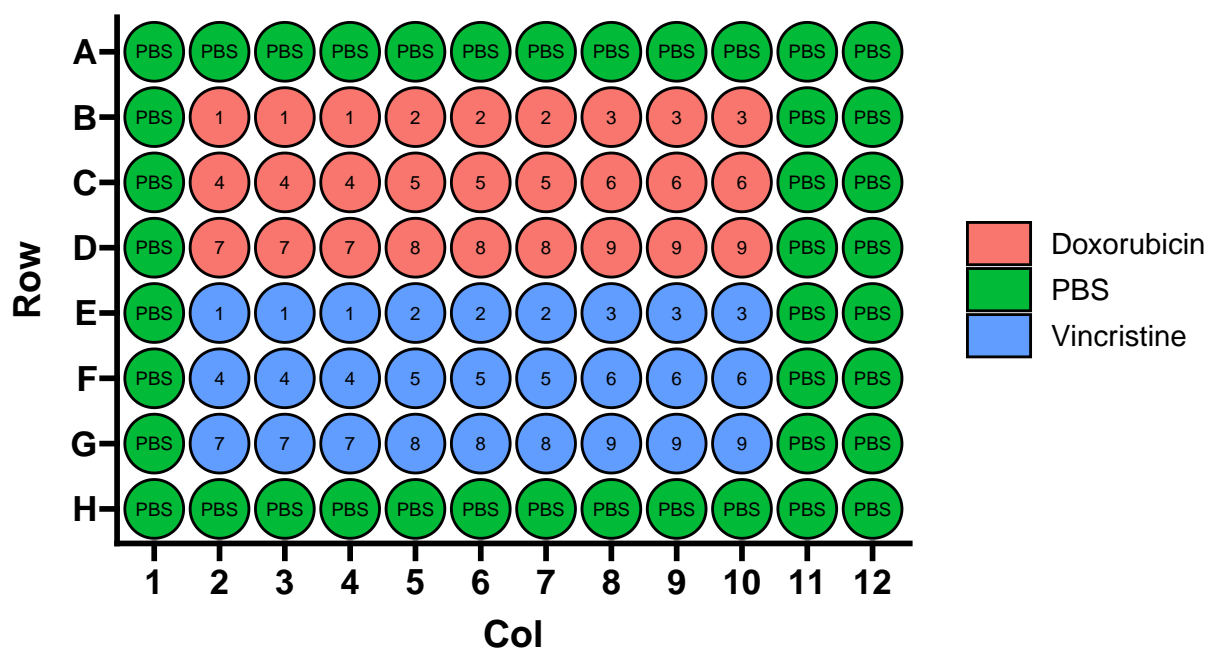

### EC50 Plate ID: EC50\_231013\_Plate2

- Methotrexate top/Cyclophosphamide bottom

Methotrexate Stock: 10mM DMSO

Drug	Dilution	Concentration	Units	Starting_Concentration
Methotrexate	1	500	nM	500000
Methotrexate	2	250	nM	250000
Methotrexate	3	125	nM	125000
Methotrexate	4	62.5	nM	62500
Methotrexate	5	31.2	nM	31200
Methotrexate	6	15.6	nM	15600
Methotrexate	7	7.81	nM	7810
Methotrexate	8	3.91	nM	3910
Methotrexate	9	1.95	nM	1950


Cyclophosphamide stock: 20 mM DMSO



Drug	Dilution	Concentration	Units	Starting_Concentration
Cyclophosphamide	1	20	uM	20000
Cyclophosphamide	2	14.3	uM	14300
Cyclophosphamide	3	10.2	uM	10200
Cyclophosphamide	4	7.29	uM	7290
Cyclophosphamide	5	5.21	uM	5210
Cyclophosphamide	6	3.72	uM	3720
Cyclophosphamide	7	2.66	uM	2660
Cyclophosphamide	8	1.9	uM	1900
Cyclophosphamide	9	1.36	uM	1360

## 96 Well Plate Map

### EC50\_231013\_plate2

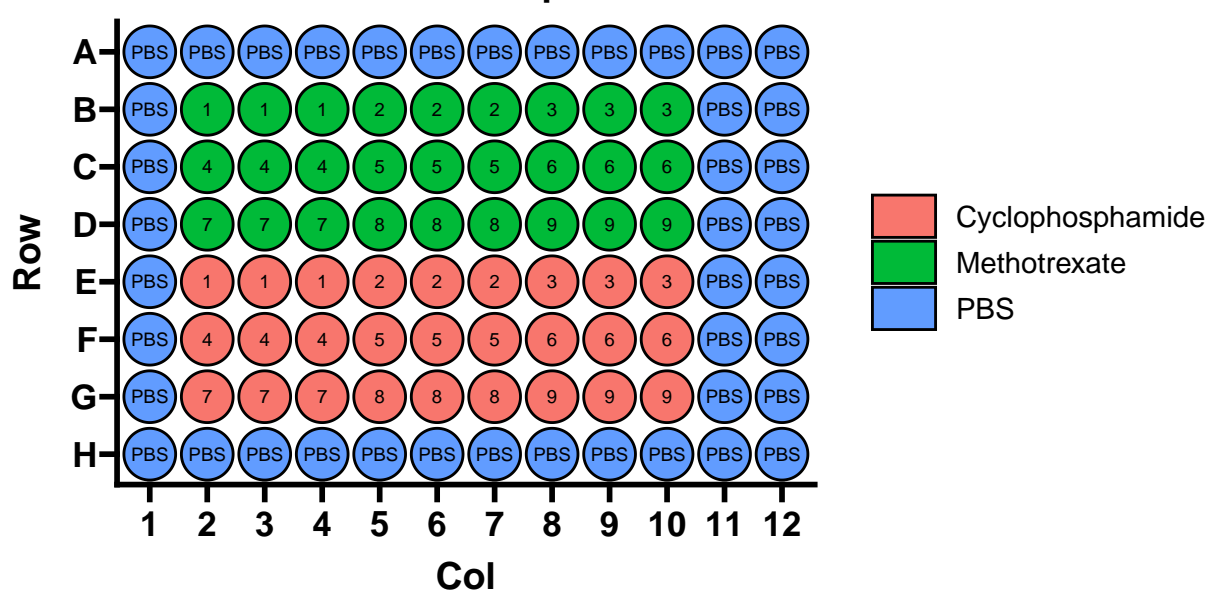

Monday 16-10-2023

Cell Culture

A20

- Split 1/5
- AR-2

N4

- Split 1/5
- AR-2

Collect EC50\_231013

Plate ID: EC50\_231013\_Plate1

- Vincristine DNW

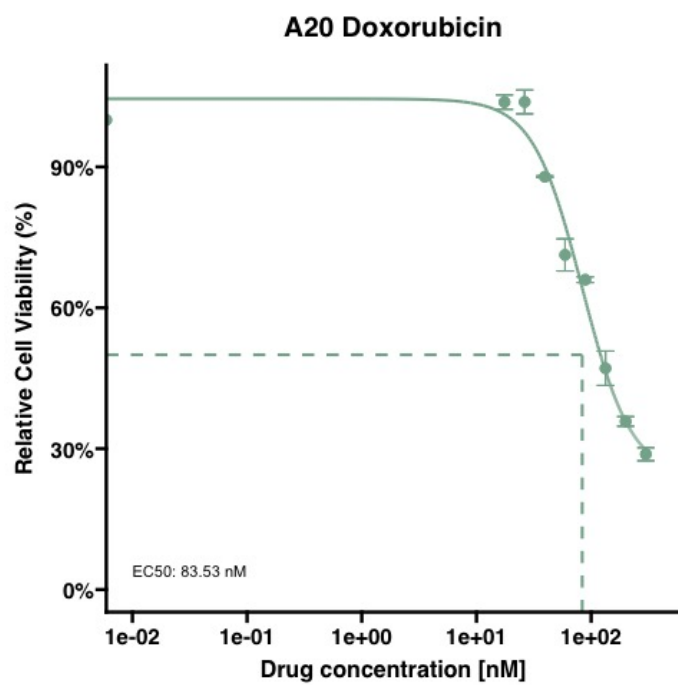


Figure 1: A20 Doxorubicin Results

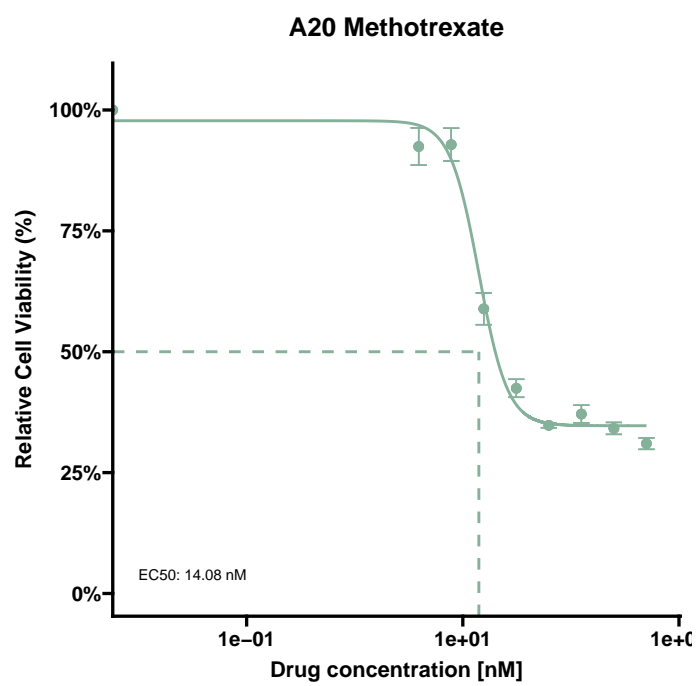


Figure 2: A20 Methotrexate Results

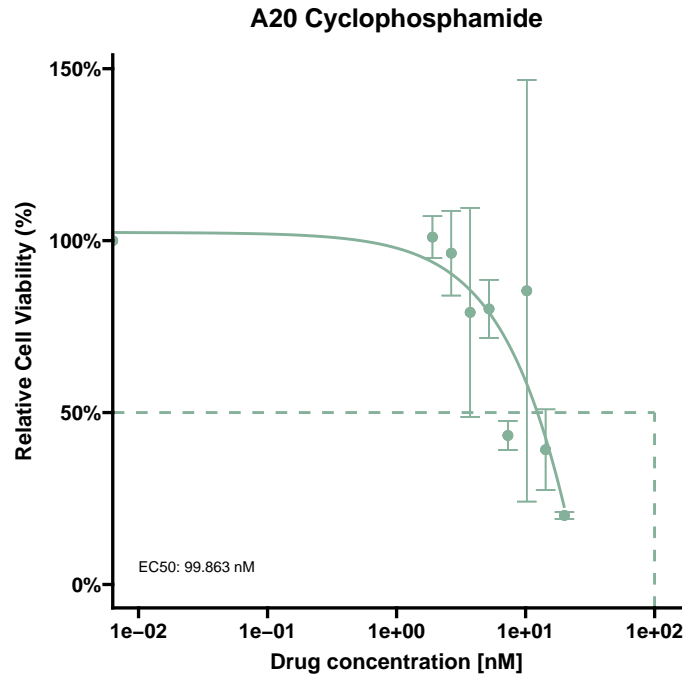


Figure 3: A20 Cyclophosphamide Results

Plate ID: EC50\_231013\_Plate2

**Thursday 19-10-2023**

**Cell Culture**

**A20**

- Split 1/5
- AR-2

**N4**

- Split 1/5
- AR-2

**Friday 20-10-2023**

**Cell Culture**

**A20**

- Split 1/5
- AR-2

**N4**

- Split 1/5
- AR-2

## Sunday 22-10-2023

### Cell Culture

#### RAMOS Barcode

- LN Tank 6, Rack 5, Box 6
- Thawed Ramos-BC 1, 5, 6

## Monday 23-10-2023

### Cell Culture

#### N4

- ~40-50% confluent
- Added 10mL AR-2

#### A20

- ~60-70% confluent
- Added 10mL AR-2

#### Ramos-BC 1

- ~50% confluent
- Added 10mL AR-2

#### Ramos-BC 5

- ~50% confluent
- Added 10mL AR-2

#### Ramos-BC 6

- ~50% confluent
- Added 10mL AR-2

### Seed Feeders (FDC)

- Count:  $4.19 \times 10^5$
- Dilute to  $4.8 \times 10^4$
- Resuspended in 5mL AR-2, added 3.73 mL
- Seeded 1x 6-well plate
  - 2 mL/well
- Seeded 1x 96-well plate (15 wells)
  - 90  $\mu$ L/well

## Tuesday 23-10-2023

### Cell Culture

#### N4

- No split, maintained

## A20

- No split, maintained

## Ramos-BC 1

- ~60% confluent
- Split 1/3
- AR-2

## Ramos-BC 5

- ~60% confluent
- Split 1/3
- AR-2

## Ramos-BC 6

- ~60% confluent
- Split 1/3
- AR-2

## FDC

- Looked attached, healthy
- Media fine

## Transduce A20/N4 with WILDseq Library

1. Made 2mL cell suspension of  $4 \times 10^6$  cells/mL
  - A20:  $1.91 \times 10^6$  count
  - N4:  $1.3 \times 10^5$  count (took 2 mL of CS from flask)
2. Add 0.5mL of cell suspension to 4 wells of 12 well plate
3. Add virus mixture (listed below) to respective wells:

Well ID	Virus	Media	HEPES
A1	25 $\mu$ L	900 $\mu$ L	50 $\mu$ L
A2	10 $\mu$ L	930 $\mu$ L	50 $\mu$ L
A3	1 $\mu$ L	948 $\mu$ L	50 $\mu$ L
A4	0 $\mu$ L	950 $\mu$ L	50 $\mu$ L

4. Centrifuged cells at 600 x g for 1.5 hours at 32C
5. Resuspended cells and transferred to 6 well plate
  - N4 were put on feeder plate which was seeded 23-10-2023

## Wednesday 25-10-23

### Attempt 1: A20/N4 Barcoding - Media change

- Changed media
- AR-2

## Made New Media: AR-3

1. Made new media: AR-3

Solution	ID code	Volume	% Total volume
Advanced RMPI	—	500 mL	78%
Glutamax	—	6.5 mL	20%
Pen-Strep	—	6.5 mL	1%
FBS	—	128 mL	1%

## Thursday 26-10-23

### Attempt 1: A20/N4 Barcoding - Sort

- DNW
- No cells were positive
- Possible causes:
  - Virus volume too low
  - Spin speed too low
  - packaging envelope incorrect for mouse cells (GALV)

## Friday 27-10-23

### Cell Culture

#### A20

- Maintained
- No Split

#### N4

- Split 1/4
- AR-3

#### Ramos-BC 1

- ~60% confluent
- Split 1/4
- AR-3

#### Ramos-BC 5

- ~60% confluent
- Split 1/4
- AR-3

#### Ramos-BC 6

- ~60% confluent
- Split 1/4
- AR-3

## EC50\_231027

### Overview

- Treated A20 with Doxorubicin, Vincristine, Cyclophosphamide, Methotrexate
- Concentrations were the same as EC50\_231013
- Repeat for biological replicates
- Made an error for plate 1

### EC50 Protocol:

1. Create 12 ml cell suspension of  $2.5 \times 10^6$  cells
  - A20 Cell count:  $2.5 \times 10^6$
  - 8.33mL AR-3 added to 1.66mL of flask CS
2. Add 90  $\mu$ L of CS to respective wells of a 96-well plate
3. Dilute drug to starting concentration indicated in DMSO
4. Serial dilute into DMSO 8 times
5. Add 10uL of diluted drug to respective eppendorf containing 990uL media
6. Add 10 $\mu$ L of drugged media to respective well

### EC50 Plate ID: EC50\_231027\_Plate1

- Doxorubicin top/Vincristine bottom

Doxorubicin Stock: 10mM DMSO


Drug	Dilution	Concentration	Units	Starting_Concentration
Doxorubicin	1	300	nM	300000
Doxorubicin	2	200	nM	200000
Doxorubicin	3	133	nM	133000
Doxorubicin	4	88.9	nM	88900
Doxorubicin	5	59.3	nM	59300
Doxorubicin	6	39.5	nM	39500
Doxorubicin	7	26.3	nM	26300
Doxorubicin	8	17.6	nM	17600
Doxorubicin	9	11.7	nM	11700

Vincristine stock: 30 $\mu$ M DMSO

Drug	Dilution	Concentration	Units	Starting_Concentration
Vincristine	1	10	nM	10000
Vincristine	2	7.14	nM	7140
Vincristine	3	5.1	nM	5100
Vincristine	4	3.64	nM	3640
Vincristine	5	2.6	nM	2600
Vincristine	6	1.86	nM	1860
Vincristine	7	1.33	nM	1330
Vincristine	8	0.949	nM	949
Vincristine	9	0.678	nM	678

# 96 Well Plate Map

## EC50\_231013\_plate1

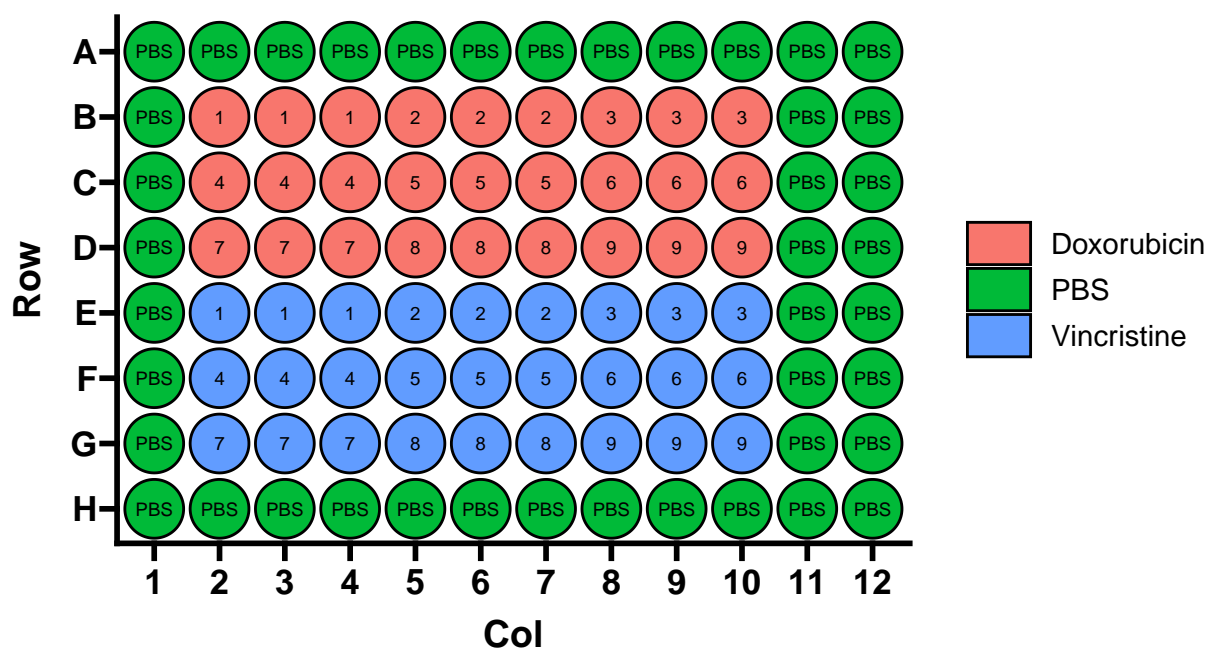

NOTE: Made an error by adding Methotrexate and Cyclophosphamide treatments to PLate 1 (as well as listed treatments) - Disregard results

## EC50 Plate ID: EC50\_231027\_Plate2

- Methotrexate top/Cyclophosphamide bottom

Methotrexate Stock: 10mM DMSO

Drug	Dilution	Concentration	Units	Starting_Concentration
Methotrexate	1	500	nM	500000
Methotrexate	2	250	nM	250000
Methotrexate	3	125	nM	125000
Methotrexate	4	62.5	nM	62500
Methotrexate	5	31.2	nM	31200
Methotrexate	6	15.6	nM	15600
Methotrexate	7	7.81	nM	7810
Methotrexate	8	3.91	nM	3910
Methotrexate	9	1.95	nM	1950


Cyclophosphamide stock: 20 mM DMSO



Drug	Dilution	Concentration	Units	Starting_Concentration
Cyclophosphamide	1	20	uM	20000
Cyclophosphamide	2	14.3	uM	14300
Cyclophosphamide	3	10.2	uM	10200
Cyclophosphamide	4	7.29	uM	7290
Cyclophosphamide	5	5.21	uM	5210
Cyclophosphamide	6	3.72	uM	3720
Cyclophosphamide	7	2.66	uM	2660
Cyclophosphamide	8	1.9	uM	1900
Cyclophosphamide	9	1.36	uM	1360

## 96 Well Plate Map

### EC50\_231013\_plate2

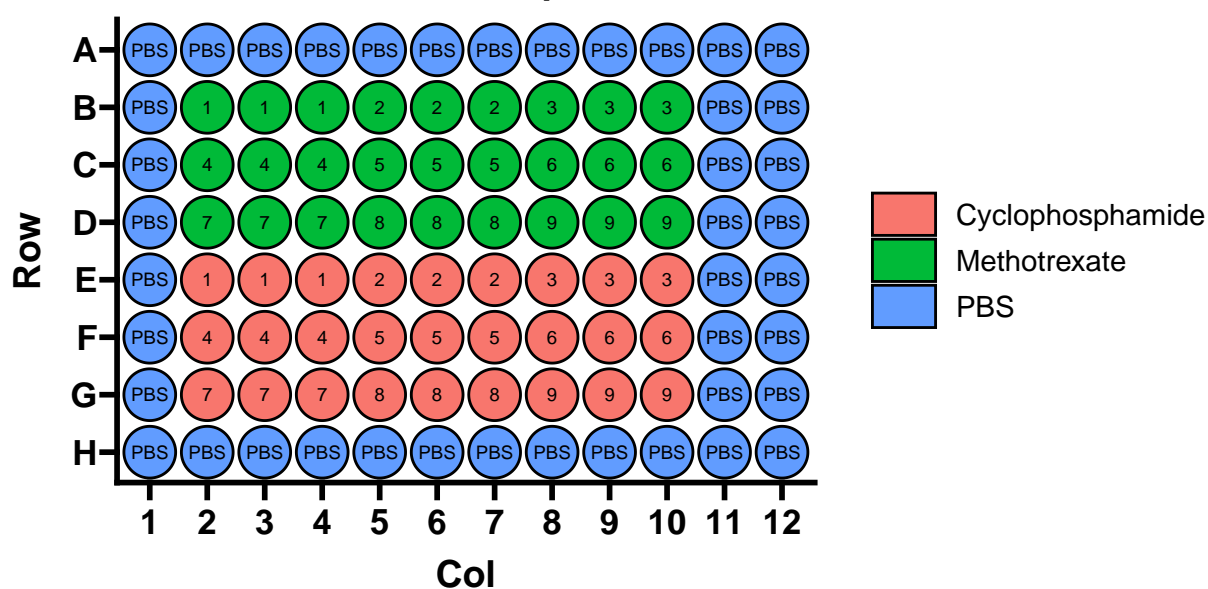

Monday 27-10-23

Cell Culture

A20

- Split 1/2
- AR-3

N4

- Split 1/2
- AR-3

Ramos-BC 1

- Split 1/2
- AR-3

## Ramos-BC 5

- Split 1/2
- AR-3

## Ramos-BC 6


- Split 1/2
- AR-3

## Seed Feeders

- Repeating barcoding of A20/N4
- Count:  $3.64 \times 10^5$
- Dilute to  $4.8 \times 10^4$
- Resuspended in 5mL AR-2, added 2.58 mL
- Seeded 1x 6-well plate
  - 2 mL/well
- Seeded 1x 96-well plate (15 wells)
  - 90  $\mu$ L/well

## Collect EC50\_231027

Plate ID: EC50\_231027\_Plate

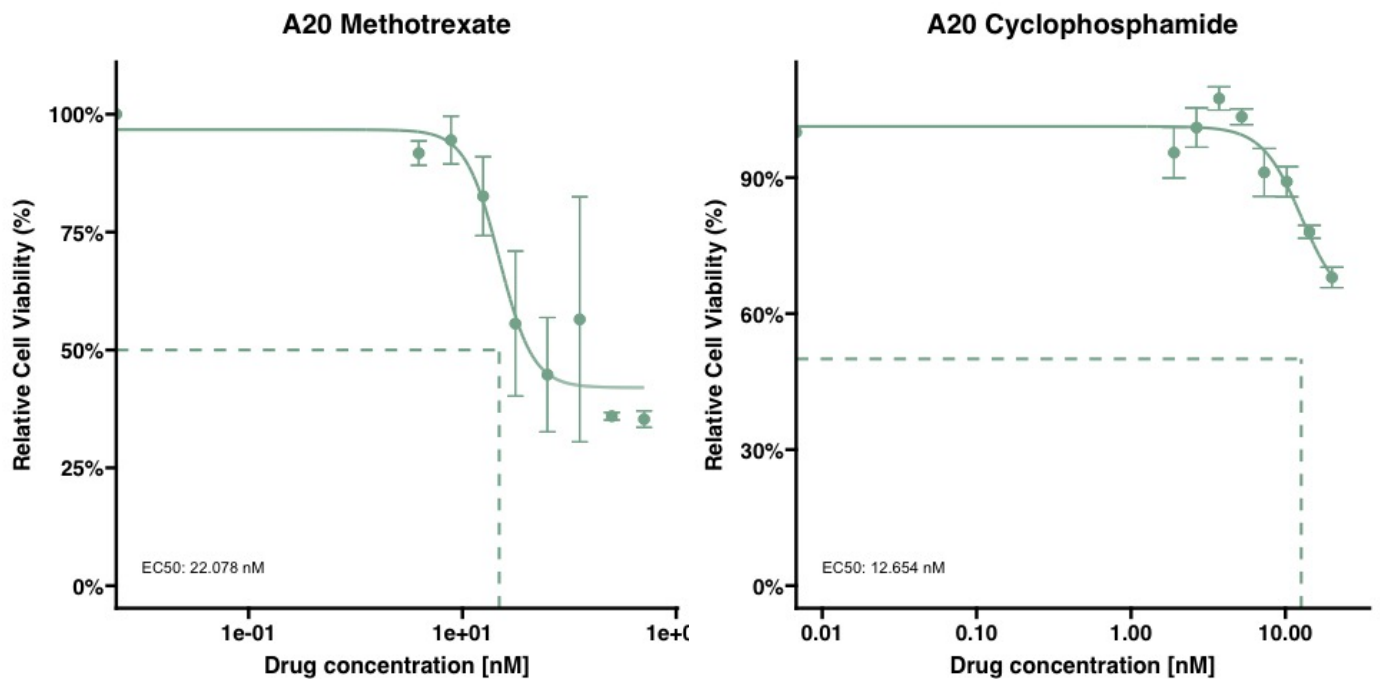

Tuesday 31-10-2023

## Attempt 2: Transduce A20/N4 with WILDseq Library

1. Made 2mL cell suspension of  $4 \times 10^6$  cells/mL
  - A20:  $2.18 \times 10^6$  count
  - N4:  $3.06 \times 10^5$  count (took 2 mL of CS from flask)

2. Add 0.5mL of cell suspension to 4 wells of 12 well plate
3. Add virus mixture (listed below) to respective wells:

Well ID	Virus	Media	HEPES
A1	250 $\mu$ L	225 $\mu$ L	25 $\mu$ L
A2	150 $\mu$ L	150 $\mu$ L	25 $\mu$ L
A3	50 $\mu$ L	50 $\mu$ L	25 $\mu$ L
A4	0 $\mu$ L	0 $\mu$ L	25 $\mu$ L

4. Centrifuged cells at 1000 x g for 1.5 hours at 0C

- Mistakenly set centrifuge to 0C

5. Resuspended cells and transferred to 6 well plate

- N4 were put on feeder plate which was seeded 23-10-2023

## Drug Pressure Experiment 1: RAMOS Day 1

### Overview

- Designing combination treatment to reflect common BL treatments in clinic
- Plan to analyse differences in barcode representation between baseline, DMSO, and treatment-resistant
- Drug concentrations:

Drug	Treatment Day	Concentration	Stock
Doxorubicin	Day 2	80nM	100 $\mu$ M working stock
Vincristine	Day 1	10nM	30 $\mu$ M
Methotrexate	Day 1	20nM	100 $\mu$ M working stock
Cyclophosphamide	Day 2/3/4	20 $\mu$ M	20 $\mu$ M working stock


- Plan:
  1. Collect sample at baseline
  2. Complete 4 day treatment course of RAMOS-BC cells (with corresponding DMSO controls)
  3. Allow treatment(T) cells to recover and collect sample
  4. Complete 4 day treatment course again and collect sample
- Possible future experiments:
  - EC50's of T vs WT cells (are T RAMOS cells drug resistant?)

## Day 1

### Seeded RAMOS BC

- Seeded 3x 6 well plates with RAMOS-BC 1, 5, 6 respectively

## 96 Well Plate Map

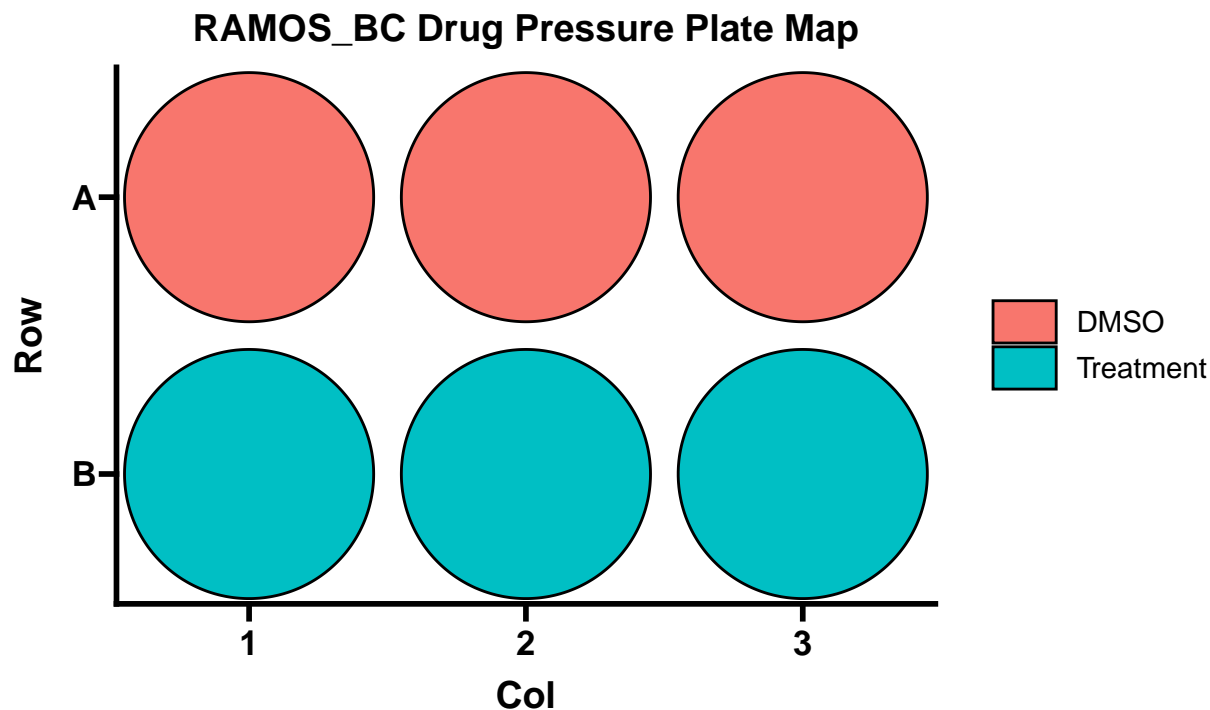

- Seeded  $0.25 \times 10^6$  cells/well in 3mL AR-3
- Treatment wells

Drug	Treatment Day	Concentration	Stock	Stock Volume	Media Volume
DMSO	Day 1	-	-	14.4 $\mu$ L	27 mL AR-3
Vincristine	Day 1	10nM	30 $\mu$ M	9 $\mu$ L	27 mL AR-3
Methotrexate	Day 1	20nM	100 $\mu$ M working stock	5.4 $\mu$ L of WS	27 mL AR-3

### Collected RAMOS Baseline

- Took 1mL CS from respective flask and transferred to eppendorf
- Spun down 2500 x g at 4C
- Removed Media
- Snap froze in LN
- Put in -80 CM Box 1

**NOVEMBER 2023**

**Wednesday 01-11-2023**

### Drug Pressure Experiment 1: RAMOS Day 2

- Changed media
- Spun down at 600 x g for 5 min
- Resuspended in Day 2 Treatment:

Drug	Treatment Day	Concentration	Stock	Stock Volume	Media Volume
DMSO	Day 1	-	-	48.60 $\mu$ L	27 mL AR-3
Doxorubicin	Day 1	80nM	100 $\mu$ M working stock	21.60 $\mu$ L of working stock	27 mL AR-3
Cyclophosphamide	Day 1	20 $\mu$ M	20mM	27 $\mu$ L	27 mL AR-3

## Attempt 2: Barcoding A20/N4

- Sanity check
- Ran samples on E6 flow
- A20 negative
- Suggests envelope protein is not compatible with mouse cells
- N4 dead
- Re-make WILDseq virus with VSVG packaging envelope

## Thursday 02-11-2023

### Cell Culture

#### A20

- Split 1/2
- AR-3

#### N4

- Split 1/2
- AR-3

#### Ramos-BC 1

- Split 1/4
- AR-3

#### Ramos-BC 5

- Split 1/4
- AR-3

#### Ramos-BC 6

- Split 1/4
- AR-3

### Seeded Feeders

- 1 vial
- Split into 2 T75 flasks

## Drug Pressure Experiment 1: RAMOS Day 3

- Changed media
- Spun down at 600 x g for 5 min
- Resuspended in Day 3 Treatment:

Drug	Treatment Day	Concentration	Stock	Stock Volume	Media Volume
DMSO	Day 1	-	-	27 $\mu$ L	27 mL AR-3
Cyclophosphamide	Day 1	20 $\mu$ M	20mM	27 $\mu$ L	27 mL AR-3

## Friday 03-11-2023

### Cell Culture

#### A20

- Split 1/2
- AR-3

#### N4

- Split into two new feeder flasks
  - Seeded 02-11-2023
  - N4.1 / N4.2
  - 1 flask to transduce, 1 flask to freeze down
- AR-3

#### Ramos-BC 1

- Add 10mL AR-3

#### Ramos-BC 5

- Add 10mL AR-3

#### Ramos-BC 6

- Add 10mL AR-3

### Drug Pressure Experiment 1: RAMOS Day 4

- Treatment cells too sick to continue
- Changed media
- Spun down at 600 x g for 5 min
- Resuspended both DMSO and Treatment cells in fresh AR-3
- Allow to grow

## Monday 06-11-2023

### Cell Culture

#### A20

- Split 1/2
- AR-3

#### N4.1

- Healthy/growing
- Leave till tomorrow

## N4.2

- Healthy/growing
- Leave till tomorrow

## Ramos-BC 1

- Froze down
- 3x vials
- FBS +10% DMSO

## Ramos-BC 5

- Froze down
- 3x vials
- FBS +10% DMSO

## Ramos-BC 6

- Froze down
- 3x vials
- FBS +10% DMSO

## Seed Feeders

- For WILDseq transduction of N4 (tuesday?)
- 1x 4 wells of 6-well plate, 15wells of 96-well plate
- Dilute in 8ml
- 2mL/well of 6well
- 90 $\mu$ L/well of 96well plate

## Drug Pressure Experiment 1: RAMOS Day 7

- Split DMSO cells 1/6
- Add 1mL AR-3 to treatment cells
  - treatment cells still sick, allow to grow for a few more days
  - Perhaps transfer to smaller well?

**Tuesday 07/11/2023**

## Made New Media: AR-4

1. Made new media: AR-4

Solution	ID code	Volume	% Total volume
Advanced RMPI	—	500 mL	78%
Glutamax	—	6.5 mL	20%
Pen-Strep	—	6.5 mL	1%
FBS	—	128 mL	1%

## Cell Culture

### A20

- Split 1/6
- AR-4

#### N4.1

- Healthy/growing
- Leave till tomorrow

#### N4.2

- Froze down
- 2x vials
- FBS +10% DMSO

### Drug Pressure Experiment 1: RAMOS Day 5

- Split D 1/6

#### Attempt 3: Transduce N4 with WILDseq Library

1. Made 2mL cell suspension of  $4 \times 10^6$  cells/mL
  - N4: took 10ml of 20ml CS, spun down and resuspended in 2mL AR-4
2. Add 0.5mL of cell suspension to 4 wells of 12 well plate
3. Add virus mixture (listed below) to respective wells:

Well ID	Virus	Media	HEPES
A1	250 $\mu$ L	225 $\mu$ L	25 $\mu$ L
A2	150 $\mu$ L	150 $\mu$ L	25 $\mu$ L
A3	50 $\mu$ L	50 $\mu$ L	25 $\mu$ L
A4	0 $\mu$ L	0 $\mu$ L	25 $\mu$ L

4. Centrifuged cells at 1000 x g for 1.5 hours at 32C
5. Resuspended cells and transferred to 6 well plate
  - N4 were put on feeder plate which was seeded 06-11-2023

### Wednesday 08/11/2023

#### Attempt 3: Barcoding N4

- Changed Media
- Cells look healthy ## Drug Pressure Experiment 1: RAMOS Day 6
- Checked treatment cells with trypan blue
- Everything was dead
- Ended experiment: will reconsider approach
  - Treat Ramos with individual agents as well as combination
  - Treat for 3 days and recover or 12 days and recover(?)

### Thursday 09/11/2023

#### Attempt 3: Barcoding N4

- Booked sort for tomorrow
- Changed media
- Sanity check on E6 Fortessa (very low MOI) ## Cell Culture



## A20

- Split 5mL CS into 15mL of AR-4 ### N4
- Split 5mL CS into 30mL of AR-4

## Friday 10/11/2023

### Attempt 3: Barcoding N4

- Sorted N4
- 1000 cells into wells A1-A3 of 96-well plate seeded with feeders
  - Seeded 06/11/2023
- 500 cells into wells B1-B3
  - Seeded 06/11/2023

## Cell Culture

### RAMOS-BC1

- Thawed vial of RAMOS-BC1 into 15mL of AR-4 ### N4
- Threw away ### A20
- Threw away

## Monday 13/11/2023



### Start RAMOS-Drug Pressure v2


1. Count RAMOS-BC1 CS
  - CS:  $1 \times 10^6$  cells/mL
2. Add 122.8 $\mu$ L CS per well
3. Make drug dilutions according to following table:
  - Make working stock of Doxorubicin and Methotrexate
  - Dilute 10mM stock to 100 $\mu$ M
    - 10 $\mu$ L stock in 990  $\mu$ L DMSO

DMSO: 16.67  $\mu$ L DMSO in 10mL Combo: Added doxorubicin, vincristine, methotrexate, and cyclophosphamide to 10mL at same volumes indicated below

Drug	Concentration	Stock	Dilution
Doxorubicin	80nM	100 $\mu$ M working stock (10mM Stock)	5 $\mu$ L WS in 10mL
Vincristine	10nM	30 $\mu$ M	1.67 $\mu$ L stock in 10mL
Methotrexate	20nM	100 $\mu$ M working stock (10mM Stock)	5 $\mu$ L WS in 10mL
Cyclophosphamide	20 $\mu$ M	20 $\mu$ M working stock	5 $\mu$ L stock in 10mL

4. Add 3 mL of drugged media to respective wells

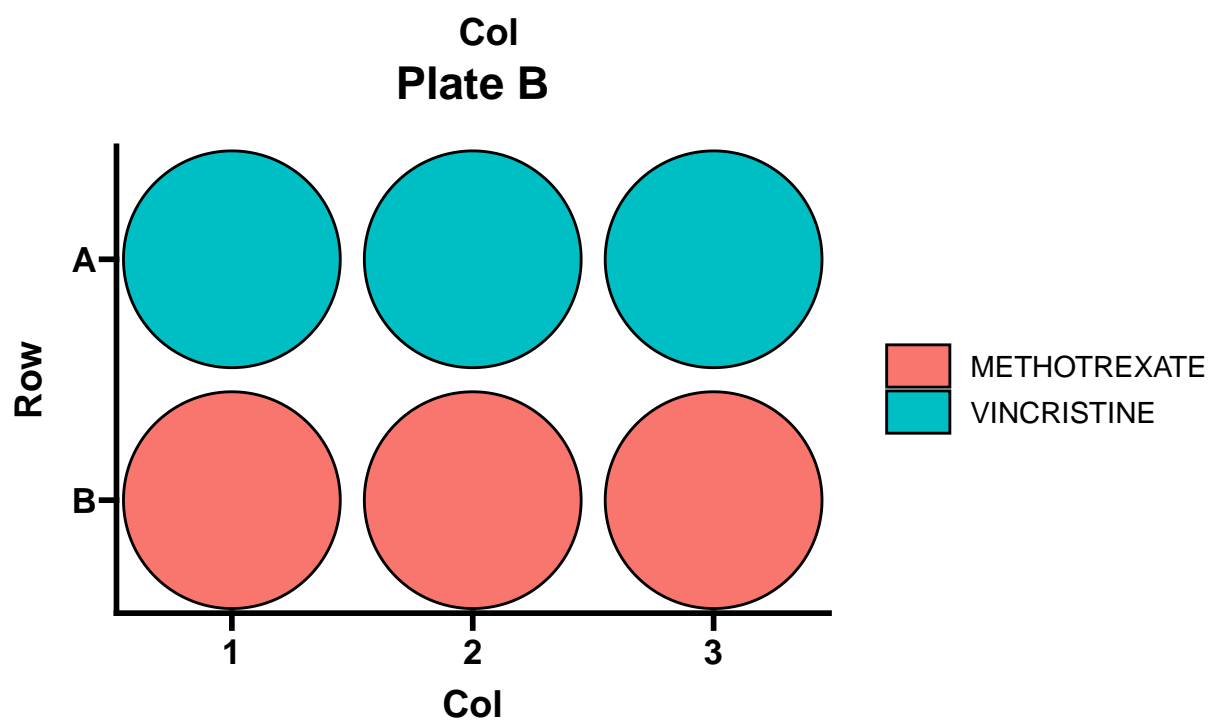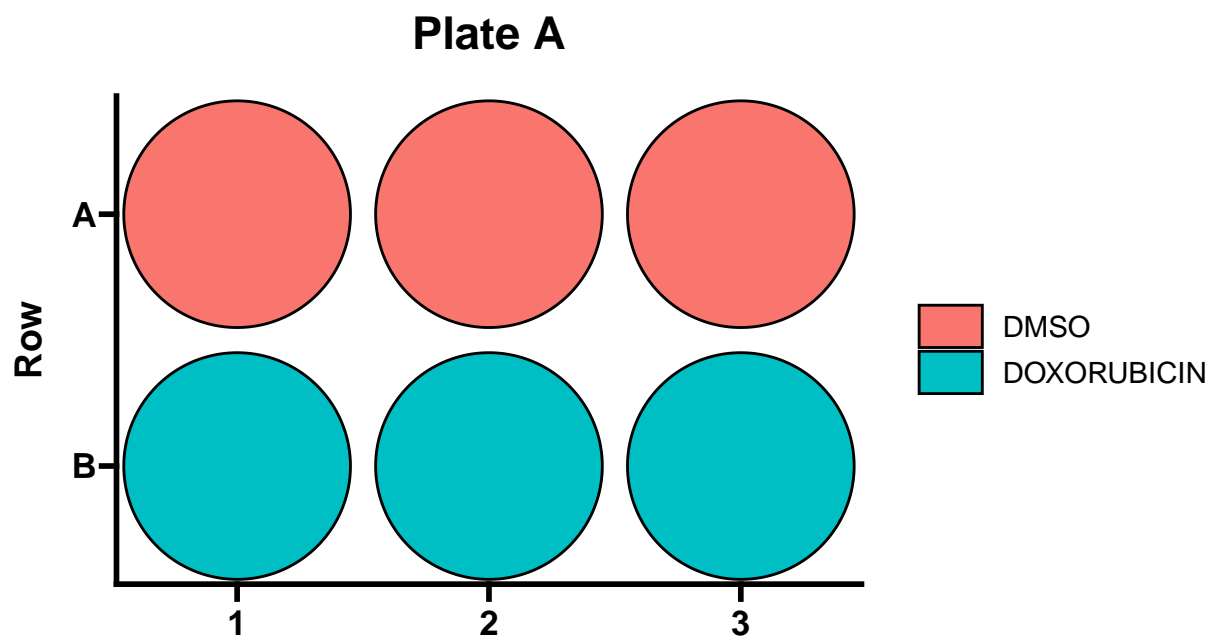

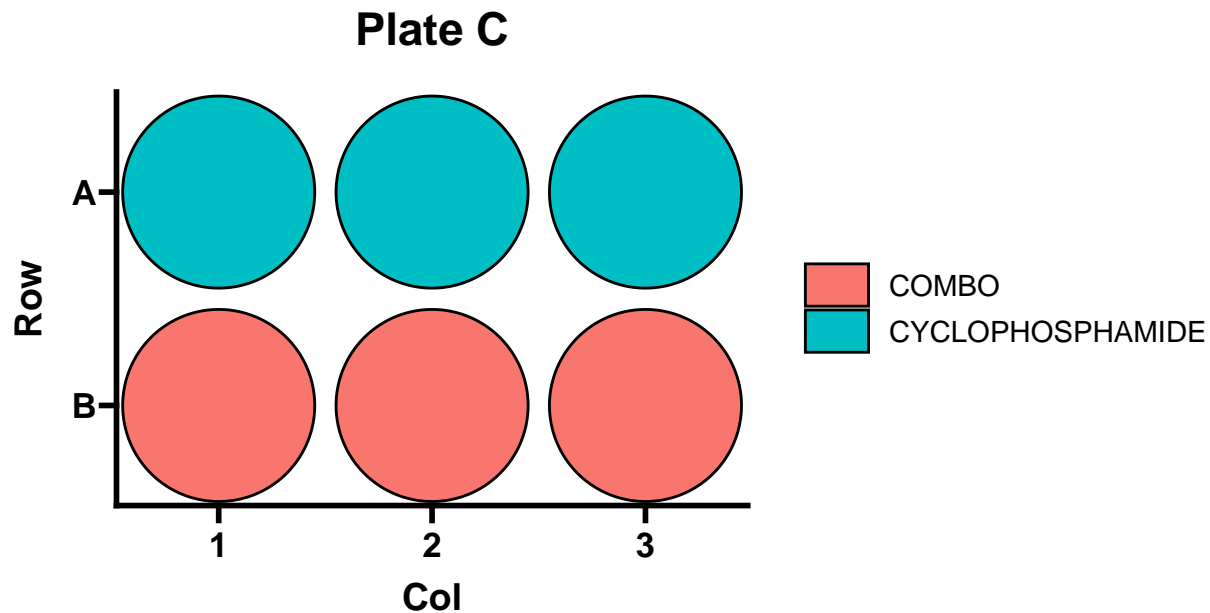

## Cell Culture

### N4-BC Attempt 3

- Add media to wells with cells
- Lots of death, some clusters of live cells ### RAMOS-BC1
- Split 4mL CS into 20mL AR-4

### PIL COURSE

**Thursday 13/11/2023**

### RAMOS-DPv2

Change media: 1. Transfer into 15ml Falcon 2. Spun down @ 500xg for 6min 3. Resuspended in 1mL 4. Transferred into new 12-well plate

**Friday 17/11/2023**

## Cell Culture

**INFECTION** - RAMOS-BC1 and RAMOS-DPv2 infected - Threw out cultures and AR-4 - N4 seemingly unaffected

**Monday 20/11/2023**

## WILDseq - Mouse Experiment 1

### Overview

- Chris Steel injected mice with barcoded N2 cells
- He applied various treatment conditions and tracked tumour growth over time
- Collected tumours at humane endpoint for barcode analysis and scRNA-seq ### RNA extraction WILDseq-ME1
- Performed RNA extraction of tissue/cell samples according to the following kit:

Sample ID	Treatment/Sample	ng/uL	i7 index	i5 index
1	Cyclophosphamide	773.7	N701	S502
2	Cyclophosphamide	545.3	N702	S502
3	Cyclophosphamide	903.6	N703	S502
4	Cyclophosphamide	1056.1	N704	S502
5	Cyclophosphamide	959.1	N705	S502
6	Combination	730.1	N706	S502
7	Combination	602.3	N707	S502
8	Combination	449.6	N710	S502
9	Combination	1002.1	N701	S503
10	Combination	1929.2	N702	S503
11	Methotrexate	861.5	N703	S503
12	Methotrexate	1110.4	N704	S503
13	Methotrexate	1171.0	N705	S503
14	Methotrexate	1347.8	N706	S503
15	Methotrexate	891.4	N707	S503
16	Vehicle	374.5	N710	S503
17	Vehicle	911.2	N701	S505
18	Vehicle	829.8	N702	S505
19	Vehicle	600.1	N703	S505
20	Vehicle	750.3	N704	S505
21	Baseline	401.5	N705	S505
22	Baseline	443.3	N706	S505
23	Baseline	373.4	N707	S505
24	Baseline	444.5	N710	S505
25	Baseline	267.6	N701	S506
26	BLLW 14K	378.7	N702	S506
27	BLLW 2K	362.7	N703	S506
28	BLLW 1K	563.8	N704	S506
29	Methotrexate (outlier)	348.1	N705	S506
30	Combo	380.8	N706	S506